



GE HealthCare

SIGNA™ Architect

(Platform Configuration)

Preinstallation Manual

5670003
Revision 15

Contents

- 1 Introduction 9**
 - 1.1 Preinstall Manual Introduction 9
 - 1.1.1 Document Purpose 9
 - 1.1.2 Intended User 9
 - 1.1.3 Who Should Read This Manual 10
 - 1.1.4 Related Publications 12
 - 1.1.5 Document Overview 12
 - 1.2 Symbols Key 13
- 2 General System Level 14**
 - 2.1 System Level Requirements for Installing into Existing MR Suite 14
 - 2.2 System components 16
 - 2.2.1 Magnet Room 16
 - 2.2.2 Equipment Room 16
 - 2.2.3 Control Room 16
 - 2.2.4 Accessories 17
 - 2.2.5 System Overview 17
 - 2.3 MR Suite Minimum Room Size Requirements 18
 - 2.4 MR System Seismic Requirements 20
 - 2.5 Structure-borne Vibration Control Specifications 21
 - 2.6 MR Suite Magnetic Field Specifications 24
 - 2.6.1 Magnetic Fringe Field 24
 - 2.6.2 Interference from Changing Magnetic Fields 31
 - 2.6.3 Electrical Current 35
 - 2.6.4 Non-MR System Equipment Sensitivity to Magnetic Fields 36
 - 2.7 Multiple MR System Requirements 37
 - 2.7.1 Multiple Magnets 37
 - 2.7.2 Shared Equipment Rooms 37
 - 2.8 MR Suite Temperature and Humidity 39
 - 2.8.1 Temperature and Humidity Requirements 39
 - 2.8.2 Equipment Heat Output Specifications 40
 - 2.9 Facility Coolant Requirements 41
 - 2.9.1 Integrated Cooling Cabinet (ICC) Coolant Requirements 41
 - 2.9.2 Requirements for Emergency Backup Facility Coolant (Optional) 45
 - 2.10 MR Suite Electrical Requirements 47
 - 2.10.1 General Electrical Requirements 47
 - 2.10.2 GE HealthCare supplied Main Disconnect Panel (MDP) Specifications for M7100ZA and M7100ZB 50
 - 2.10.3 Customer-supplied Main Disconnect Panel (MDP) Requirements (exempt countries only*) 52
 - 2.10.4 Emergency Power Backup Specifications (Optional) 55
 - 2.11 MR System Shipping and Receiving 56

2.11.1 Receiving Requirements56

2.11.2 Facility Delivery Route Requirements56

2.11.3 MR System Component Shipping Specifications57

2.11.4 Temperature and Humidity Storage Requirements61

3 Magnet Room 62

3.1 Magnet Room Introduction.....62

3.2 Magnet Room Structural Requirements.....64

 3.2.1 Overview.....64

 3.2.2 Environmental Steel Limits64

 3.2.3 Vibration Requirements65

3.3 Magnetic Shielded Room Requirements.....67

3.4 Penetration Panel Wall Opening Requirements67

3.5 Finished Room Requirements70

 3.5.1 Ferrous Materials in the Magnet Room70

 3.5.2 Walls.....70

 3.5.3 Magnet Preinstallation Markings.....70

 3.5.4 Penetration Panel Closet74

 3.5.5 Penetration Panel Closet Specifications75

 3.5.6 Doors, Magnet Access Openings, and Patient Viewing Windows.....76

 3.5.7 Finished Ceiling.....76

 3.5.8 Magnet Room Floors77

 3.5.9 Storage Cabinets82

3.6 Magnet Room Equipment Specifications.....83

 3.6.1 Magnet (MAG) Assembly Specifications83

 3.6.2 Patient Table (PT) Specifications.....84

 3.6.3 Optional Surgical Suite Transmobile Table86

 3.6.4 Magnet Rundown Unit (MRU) Specifications and Requirements86

 3.6.5 Oxygen Monitor Option Sensor Specifications87

3.7 Magnet Room Lighting Requirements.....88

4 Equipment Room 89

4.1 Equipment Room Overview89

4.2 Main Disconnect Panel (MDP) Requirements and Specifications92

 4.2.1 Requirements92

 4.2.2 Specifications92

4.3 Power, Gradient, RF Cabinet (PGR) Specifications94

4.4 Penetration Cabinet (PEN) Specifications.....96

4.5 Integrated Cooling Cabinet (ICC) Specifications.....98

4.6 Secondary Penetration Wall (SPW) Specifications102

4.7 Magnet Monitor (MON) Requirements and Specifications103

 4.7.1 Requirements103

 4.7.2 Specifications103

- 4.8 Magnetic Resonance Elastography (MRE) Specifications (Optional Equipment).....105
 - 4.8.1 Requirements.....105
 - 4.8.2 Specifications.....105
- 4.9 Multi-Nuclear Spectroscopy (MNS) Specifications (Optional Equipment).....106
- 5 Control Room 108**
 - 5.1 Operator Workspace Equipment Specifications.....108
 - 5.1.1 Operator Workspace Assembly108
 - 5.1.2 Operator Workspace (OW) (Optional Equipment)109
 - 5.1.3 Global Operator Cabinet (GOC)110
 - 5.1.4 Host Display.....111
 - 5.1.5 Host Keyboard.....111
 - 5.1.6 Pneumatic Patient Alert.....112
 - 5.2 Oxygen Monitor (OXY) Specifications (Optional Equipment)113
 - 5.3 Physiological Acquisition Transceiver (PAT) Charging Station Specification (Optional Equipment)114
- 6 Digital Service and Connectivity Requirements 115**
 - 6.1 InSite RSvP (Remote Service Platform) Requirements.....115
 - 6.1.1 InSite RSvP Connectivity Requirements115
- 7 MR System Interconnects 116**
 - 7.1 MR System Interconnects Specifications116
 - 7.1.1 Component Designator Definitions116
 - 7.1.2 Available Cable Lengths.....117
 - 7.1.3 Multi-Nuclear Spectroscopy (MNS) Option120
 - 7.1.4 BrainWave Option121
 - 7.1.5 Magnetic Resonance Elastography (MRE) Option121
 - 7.1.6 Storage Requirements for Excess Gradient Cable121
 - 7.2 MR System Interconnects Routing Requirements125
 - 7.2.1 General Requirements125
 - 7.2.2 Magnet Room Requirements126
 - 7.2.3 Recommended Cable Groupings127
 - 7.2.4 Cable Tray Requirements and Examples.....127
 - 7.2.5 Penetration Panel Closet Requirements130
 - 7.2.6 Equipment Room Requirements.....131
 - 7.3 Facility-Supplied System Interconnects Specifications132
- 8 Appendix..... 135**
 - 8.1 Glossary135
 - 8.2 MR Site Vibration Test Guidelines137
 - 8.2.1 Test Measurements137
 - 8.2.2 Equipment (Spectral Analyzer) Set-Up137
 - 8.2.3 Data Collection.....138
 - 8.2.3.1 Ambient Baseline Condition138
 - 8.2.3.2 Normal Condition138
 - 8.2.4 Presentation/Interpretation of Results138

8.3 Sample Calculation AC Power Equipment Minimum Distance	142
8.4 Selecting Anchor Size.....	144
8.5 Sample control schematic for customer-supplied MDP.....	145
Revision History.....	147

Language Policy

DOC0371395 - Global Language Procedure

ПРЕДУПРЕЖ ДЕНИЕ (BG)	Това ръководство е налично само на китайски (ZH-CN), английски, френски, немски, японски, корейски, полски, португалски (PT-BR), руски, испански и виетнамски. Ако доставчикът на услуги на даден клиент изисква език, който е различен от тези езици, отговорност на клиента е да предостави преводачески услуги.
警告 (ZH-CN)	本手册仅提供中文 (ZH-CN)、英文、法语、德语、意大利语、日语、韩语、波兰语、葡萄牙语 (PT-BR)、俄语、西班牙语和越南语版本。如果客户的服务提供商需要其他语言，则客户有责任提供翻译服务。
警告 (ZH-HK)	本手冊僅提供中文 (ZH-CN)、英文、法文、德文、意大利文、日文、韓文、波蘭文、葡萄牙文 (PT-BR)、俄文、西班牙文及越南文版本。如客戶的服務供應商需要這些語言以外的版本，則相關客戶有責任提供有關的翻譯服務。
警告 (ZH-TW)	此手冊僅提供中文 (ZH-CN)、英文、法文、德文、義大利文、日文、韓文、波蘭文、葡萄牙文 (PT-BR)、俄文、西班牙文和越南文版本。假如客戶的服務提供者所需語言版本不在所列語言之中，客戶需自行負責提供翻譯服務。
UPOZORENJE (HR)	Ovaj je priručnik dostupan samo na kineskom (ZH-CN), engleskom, francuskom, njemačkom, talijanskom, japanskom, korejskom, poljskom, portugalskom (PT-BR), ruskom, španjolskom i vijetnamskom jeziku. Ako klijentov serviser zahtijeva jezik koji nije jedan od tih jezika, odgovornost je klijenta pružiti uslugu prevođenja.
VÝSTRAHA (CS)	Tato příručka je k dispozici pouze v čínštině (ZH-CN), angličtině, francouzštině, němčině, italštině, japonštině, korejštině, polštině, portugalštině (PT-BR), ruštině, španělštině a vietnamštině. Pokud poskytovatel služeb zákazníka vyžaduje jiný jazyk než tyto jazyky, je odpovědností zákazníka poskytovat překladatelské služby.
ADVARSEL (DA)	Denne vejledning findes kun på kinesisk (ZH-CN), engelsk, fransk, tysk, italiensk, japansk, koreansk, polsk, portugisisk (PT-BR), russisk, spansk og vietnamesisk. Hvis en kundes tjenesteudbyder kræver et andet sprog end disse sprog, er det kundens ansvar at levere oversættelsestjenester.
WAARSCHUWING (NL)	Deze handleiding is alleen beschikbaar in het Chinees (ZH-CN), Engels, Frans, Duits, Italiaans, Japans, Koreaans, Pools, Portugees (PT-BR), Russisch, Spaans en Vietnamees. Als de serviceprovider van een klant een andere taal dan deze talen vereist, is het de verantwoordelijkheid van de klant om vertaalservices te leveren.
WARNING (EN)	This manual is available in Chinese (ZH-CN), English, French, German, Italian, Japanese, Korean, Polish, Portuguese (PT-BR), Russian, Spanish, and Vietnamese only. If a customer's service provider requires a language other than these languages, it is the customer's responsibility to provide translation services.
HOIATUS (ET)	See juhend on saadaval ainult hiina (ZH-CN), inglise, prantsuse, saksa, itaalia, jaapani, korea, poola, portugali (PT-BR), vene, hispaania ja vietnami keeles. Kui kliendi teenusepakkujal on vaja juhendit mõnes muus keeles, on tõlketeenuste osutamine kliendi kohustus.
VAROITUS (FI)	Tämä opas on saatavilla vain kiinaksi (ZH-CN), englanniksi, ranskaksi, saksaksi, italiaksi, japaniksi, koreaksi, puolaksi, portugaliksi (PT-BR), venäjäksi, espanjaksi ja vietnamiksi. Jos asiakkaan palveluntarjoaja edellyttää muuta kuin näitä kieliä, käänöspalveluiden tarjoaminen on asiakkaan vastuulla.
ATTENTION (FR)	Ce manuel est disponible uniquement en allemand, anglais, chinois (ZH-CN), coréen, espagnol, français, italien, japonais, polonais, portugais (PT-BR), russe et vietnamien. Si le prestataire de services d'un client nécessite que le manuel soit rédigé dans une autre langue que celles mentionnées ci-dessus, il incombe au client de le faire traduire.
WARNUNG (DE)	Dieses Handbuch ist nur auf Chinesisch (ZH-CN), Englisch, Französisch, Deutsch, Italienisch, Japanisch, Koreanisch, Polnisch, Portugiesisch (PT-BR), Russisch, Spanisch und Vietnamesisch verfügbar. Wenn ein Dienstleister des Kunden dieses in einer anderen Sprache als die genannten benötigt, liegt es in der Verantwortung des Kunden, Übersetzungsdienstleistungen zu erbringen.

ΠΡΟΕΙΔΟΠΟΙΗΣΗ (EL)	Αυτό το εγχειρίδιο είναι διαθέσιμο μόνο σε Κινεζικά (ZH-CN), Αγγλικά, Γαλλικά, Γερμανικά, Ιταλικά, Ιαπωνικά, Κορεατικά, Πολωνικά, Πορτογαλικά (PT-BR), Ρωσικά, Ισπανικά και Βιετναμέζικα. Εάν ο πάροχος υπηρεσιών ενός πελάτη απαιτεί γλώσσα που δεν συμπεριλαμβάνεται σε αυτές τις γλώσσες, αποτελεί ευθύνη του πελάτη να παρέχει υπηρεσίες μετάφρασης.
FIGYELMEZTETÉS (HU)	Ez a kézikönyv az alábbi nyelveken érhető el: angol, francia, japán, kínai (ZH-CN), koreai, lengyel, német, olasz, orosz, portugál (PT-BR), spanyol és vietnámi. Ha az ügyfél szolgáltatója ezektől eltérő nyelvű kézikönyvet szeretne, akkor az ügyfél feladata, hogy gondoskodik a megfelelő fordításról.
AÐVÖRUN (IS)	Þessi handbók er aðeins fánleg á kínversku (ZH-CN), ensku, frönsku, þýsku, ítölsku, japönsku, kóresku, pólsku, portúgölsku (PT-BR), rússnesku, spænsku og víetnömsku. Ef þjónustuaðili viðskiptavinar þarfnast annars tungumáls en þessara tungumála er það á ábyrgð viðskiptavinarins að veita þýðingarþjónustu.
AVVERTENZA (IT)	Questo manuale è disponibile solo in lingua cinese (ZH-CN), inglese, francese, tedesco, italiano, giapponese, coreano, polacco, portoghese (PT-BR), russo, spagnolo e vietnamita. Qualora un fornitore di servizi del cliente richieda una lingua diversa dall'inglese, sarà responsabilità del cliente fornire il servizio di traduzione corrispondente.
警告 (JA)	このマニュアルは、中国語 (ZH-CN)、英語、フランス語、ドイツ語、イタリア語、日本語、韓国語、ポーランド語、ポルトガル語 (PT-BR)、ロシア語、スペイン語、およびベトナム語のみで提供されています。お客様のサービスプロバイダがこれらの言語以外の言語を必要とする場合は、お客様の責任において翻訳サービスを提供してください。
경고 (KO)	이 설명서는 중국어(중국어-중국), 영어, 프랑스어, 독일어, 이탈리아어, 일본어, 한국어, 폴란드어, 포르투갈어(포르투갈어-브라질), 러시아어, 스페인어, 베트남어로만 제공됩니다. 고객의 서비스 제공자가 이 언어를 제외한 다른 언어를 요구하는 경우, 번역 서비스를 제공하는 것은 고객의 책임입니다.
BRĪDINĀJUMS (LV)	Šī rokasgrāmata ir pieejama tikai ķīniešu (ZH-CN), angļu, franču, vācu, itāliešu, japāņu, korejiešu, poļu, portugāļu (PT-BR), krievu, spāņu un vjetnamiešu valodā. Ja klientu apkalpošanas speciālistam ir nepieciešama cita valoda, kas atšķiras no šeit norādītajām, klienta pienākums ir nodrošināt tulkotājas pakalpojumu.
ĮSPĖJIMAS (LT)	Šis vadovas pateikiamas tik kinų (ZH-CN), anglų, prancūzų, vokiečių, italų, japonų, korėjiečių, lenkų, portugalų (PT-BR), rusų, ispanų ir vietnamiečių kalbomis. Jei klientų paslaugų teikėjui reikalinga kita nei šios kalba, už vertimo paslaugų suteikimą atsako klientas.
ADVARSEL (NO)	Denne håndboken er bare tilgjengelig på kinesisk (ZH-CN), engelsk, fransk, tysk, italiensk, japansk, koreansk, polsk, portugisisk (PT-BR), russisk, spansk og vietnamesisk. Hvis en kundes tjenesteleverandør krever et annet språk, er det kundens ansvar å levere en oversettelsestjeneste.
OSTRZEŻENIE (PL)	Niniejsza instrukcja jest dostępna wyłącznie w języku chińskim (ZH-CN), angielskim, francuskim, niemieckim, włoskim, japońskim, koreańskim, polskim, portugalskim (PT-BR), rosyjskim, hiszpańskim i wietnamskim. Jeśli usługodawca klienta wymaga języka, który nie został wymieniony powyżej, obowiązkiem klienta jest zapewnienie usług tłumaczeniowych.
ATENÇÃO (PT-BR)	Este manual está disponível somente em chinês (ZH-CN), inglês, francês, alemão, italiano, japonês, coreano, polonês, português (PT-BR), russo, espanhol e vietnamita. Se o prestador de serviços de um cliente necessitar de um idioma diferente dos mencionados, o fornecimento dos serviços de tradução é de responsabilidade do cliente.
ATENÇÃO (PT-PT)	Este manual está disponível apenas em alemão, chinês (ZH-CN), coreano, espanhol, francês, inglês, italiano, japonês, polaco, português (PT-BR), russo e vietnamita. Se o fornecedor de serviços de um cliente necessitar de um idioma diferente dos listados aqui, é da responsabilidade do cliente assegurar os serviços de tradução.
ATENȚIE (RO)	Acest manual este disponibil numai în limbile chineză (ZH-CN), engleză, franceză, germană, italiană, japoneză, coreeană, poloneză, portugheză (PT-BR), rusă, spaniolă și vietnameză. Dacă furnizorul de servicii al unui client solicită o limbă diferită față de aceste limbi, este responsabilitatea clientului să furnizeze servicii de traducere.

ОСТОРОЖНО! (RU)	Настоящее руководство доступно только на китайском (ZH-CN), английском, французском, немецком, итальянском, японском, корейском, польском, португальском (PT-BR), русском, испанском и вьетнамском языках. Если поставщику услуг заказчика требуется руководство на каком-либо другом языке, перевод руководства на необходимый язык осуществляется стороной заказчика.
UPOZORENJE (SR)	Ovaj priručnik dostupan je samo na kineskom (ZH-CN), engleskom, francuskom, nemačkom, italijanskom, japanskom, korejskom, poljskom, portugalskom (PT-BR), ruskom, španskom i vijetnamskom jeziku. Ako korisnik kao pružalac usluge zahteva neki drugi jezik od navedenih, njegova je dužnost da obezbedi prevod.
UPOZORNENIE (SK)	Táto príručka je dostupná len v nasledovných jazykoch: čínština (ZH-CN), angličtina, francúzština, nemčina, taliančina, japončina, kórejščina, poľština, portugalčina (PT-BR), ruština, španielčina a vietnamčina. Ak poskytovateľ služieb daného zákazníka požaduje iný ako tieto jazyky, za poskytnutie prekladateľských služieb zodpovedá zákazník.
ATENCIÓN (ES)	Este manual está disponible solo en chino (ZH-CN), inglés, francés, alemán, italiano, japonés, coreano, polaco, portugués (PT-BR), ruso, español y vietnamita. Si el proveedor de servicios de un cliente requiere un idioma distinto de estos idiomas, es responsabilidad del cliente proporcionar los servicios de traducción.
VARNING (SV)	Den här manualen finns endast tillgänglig på kinesiska (ZH-CN), engelska, franska, tyska, italienska, japanska, koreanska, polska, portugisiska (PT-BR), ryska, spanska och vietnamesiska. Om en kunds tjänsteleverantör behöver ett annat språk än dessa är det kundens ansvar att ordna med översättningstjänster.
OPOZORILO (SL)	Ta priručnik je na voljo v kitajščini (ZH-CN), angleščini, francoščini, nemščini, italijanščini, japonščini, korejščini, poljščini, portugalščini (PT-BR), ruščini, španščini in vietnamščini. Če kupčev ponudnik storitev potrebuje drug jezik, mora za prevod poskrbeti kupec.
DİKKAT (TR)	Bu kılavuz yalnızca Çince (ZH-CN), İngilizce, Fransızca, Almanca, İtalyanca, Japonca, Korece, Lehçe, Portekizce (PT-BR), Rusça, İspanyolca ve Vietnamca dillerinde mevcuttur. Müşteri servis sağlayıcısı bu dillerden başka bir dil talep ederse çeviri hizmeti sağlamak müşterinin sorumluluğundadır.
ЗАСТЕРЕЖЕННЯ (UK)	Цей посібник доступний лише китайською (ZH-CN), англійською, французькою, німецькою, італійською, японською, корейською, польською, португальською (PT-BR), російською, іспанською та в'єтнамською мовами. Якщо постачальник послуг замовника використовує мову, яку не вказано у цьому переліку, послуги з перекладу має забезпечити замовник.

1 Introduction

1.1 Preinstall Manual Introduction



(Applies to all subsections within this section)

1.1.1 Document Purpose

This preinstallation manual provides the necessary information to prepare a site for system installation. Specifically, this manual provides information:

1. To define system requirements and interactions.
2. For the effective arrangement and interconnection of system components.
3. The customer is responsible for:
 - 3.1. Compliance with all local and national codes and regulations
 - 3.2. Siting requirements for customer-specific site procedures (medical, MR, safety, and so on)
 - 3.3. Any special architectural requirements (for example, seismic codes)

The implementation of all requirements and adherence to all specifications in this manual is the responsibility of the customer or its architect and engineers. Refer any questions to the GE HealthCare Project Manager of Installation (PMI).

1.1.2 Intended User

The primary users of this manual are the customer, the customer's architectural planner, and/or the customer's contractors.

1.1.3 Who Should Read This Manual

The following personnel must be aware of the content listed in the following sections:

Table 1-1 Personnel Index















Section	Personnel						
							
	Architect	General Contractor	Customer	Electrician	Plumber	HVAC	RF Vendor
1.1 Preinstall Manual Introduction	X	X	X	X	X	X	X
2.1 System Level Requirements for Installing into Existing MR Suite	X	X					
2.2 System components	X	X					
2.3 MR Suite Minimum Room Size Requirements	X						
2.4 MR System Seismic Requirements	X	X					
2.5 Structure-borne Vibration Control Specifications	X						
2.6 MR Suite Magnetic Field Specifications	X		X				
2.7 Multiple MR System Requirements	X						
2.8 MR Suite Temperature and Humidity	X	X			X	X	
2.9 Facility Coolant Requirements	X	X			X	X	
2.10 MR Suite Electrical Requirements	X			X			
2.11 MR System Shipping and Receiving	X						
3.1 Magnet Room Introduction	X	X					
3.2 Magnet Room Structural Requirements	X						
3.2.1 Overview	X						
3.2.2 Environmental Steel Limits	X						
3.2.3 Vibration Requirements	X	X				X	
3.3 Magnetic Shielded Room Requirements	X						
3.4 Penetration Panel Wall Opening Requirements	X						


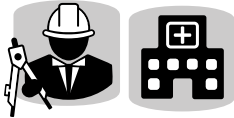

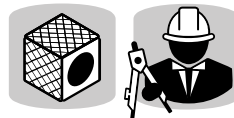
Table 1-1 Personnel Index (Table continued)

Section	Personnel						
							
	Architect	General Contractor	Customer	Electrician	Plumber	HVAC	RF Vendor
3.5.1 Ferrous Materials in the Magnet Room	X	X	X				
3.5.2 Walls	X						
3.5.3 Magnet Preinstallation Markings	X						
3.5.4 Penetration Panel Closet	X						
3.5.5 Penetration Panel Closet Specifications	X						
3.5.6 Doors, Magnet Access Openings, and Patient Viewing Windows	X						
3.5.7 Finished Ceiling	X						
3.5.8 Magnet Room Floors	X	X					
3.5.9 Storage Cabinets	X	X					
3.6 Magnet Room Equipment Specifications	X	X					
3.7 Magnet Room Lighting Requirements	X			X			
4 Equipment Room chapter	X	X					
5 Control Room chapter	X	X					
6 Digital Service and Connectivity chapter	X			X			
7.1 MR System Interconnects Specifications	X	X		X			
7.2 MR System Interconnects Routing Requirements	X	X		X			
7.3 Facility-Supplied System Interconnects Specifications	X	X		X	X	X	
8.1 Glossary	X	X	X	X	X	X	X
8.2 MR Site Vibration Test Guidelines	X						
8.3 Sample Calculation AC Power Equipment Minimum Distance	X						
8.4 Selecting Anchor Size	X	X					
8.5 Sample control schematic for customer-supplied MDP	X			X			

1.1.4 Related Publications

The preinstallation requirements in the following publications are applicable to all systems. This document and all documents referenced herein shall be provided to the Responsible Organization or Operator as a supplement to the product instructions for use and/or technical description.

Table 1-2 Additional Preinstallation Requirements

Publication Number	Title	Personnel who must be aware of the content
5850262-1EN	<i>Acoustic Room Details</i>	
5850261-1EN	<i>International Electrotechnical Commission (IEC) Electro-magnetic Compatibility (EMC)</i>	
5850263-1EN	<i>Magnet Room Venting</i>	
5850260-1EN	<i>RF Shielded Room</i>	







1.1.5 Document Overview

This manual describes requirements and specifications for the following:

1. General system requirements that apply to the entire MR suite
2. Shipping and delivery
3. Magnet Room
4. Equipment Room
5. Control Room
6. Interconnects within and between the rooms listed above

1.2 Symbols Key

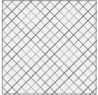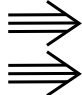Table 1-3 Symbols Key

Symbol/Unit	Definition
	Center of gravity
	Magnet isocenter
	Service area
	Airflow
	Space for airflow and cables
	Valve

2 General System Level

2.1 System Level Requirements for Installing into Existing MR Suite



When planning for the installation of this system in an existing GE HealthCare MR suite or a non-GE HealthCare MR suite, all requirements in this manual must be met because these rooms are considered new installations.

When upgrading from a Signa OpenSpeed system, extensive building updates are required. The facility must remove any Magnet Room isolation system.

1. If the existing MR suite contains a GE HealthCare system, the vibration environmental assessment must be done using the High Speed (magnetic field) Stability tool.



NOTE

The customer may have to hire a vibration consultant based on the results of the analysis.

2. Some GE HealthCare MR suites with overhead cabling may need to be modified to meet current cable tray requirements, such as minimum width and position. See [7.2 MR System Interconnects Routing Requirements on page 125](#).
3. Remove, cover, or fill in abandoned ducts or troughs from the Equipment and Magnet Rooms. Access/computer room flooring in the Equipment Room can either be removed or assessed and reinforced to support heavier cabinets.
4. Structural vibration levels may be higher at some frequencies than other MR Systems, which may increase acoustic levels. Refer to [2.5 Structure-borne Vibration Control Specifications on page 21](#).
5. The VibroAcoustic damping kit must be surface mounted (if the floor is recessed, it must be filled in and level). For upgrades that reuse the existing magnet, recessed floors can be filled with magnet spacers (25 mm (1 in.) aluminum) that were included when the magnet shipped.
6. RF vendor responsibilities:
 - 6.1. The old dock anchor cannot be reused. It must be removed and the hole filled in. The new anchor is reset after the magnet is installed. For upgrades that reuse the existing magnet, contact the PMI for further details about the potential reuse of the old dock anchor.
 - 6.2. The RF shield ceiling must support the cable routing mechanism and cables. Reinforce RF shield ceiling (see [7.2 MR System Interconnects Routing Requirements on page 125](#)).
 - 6.3. Two penetration panel openings are required and must meet the requirements in: *RF Shielded Room Requirements*, 5850260-1EN.
 - 6.4. RF shield attenuation must comply with: *RF Shielded Room Requirements*, 5850260-1EN

7. Cryogen vent may need to be relocated to align with the Magnet Cryogen Vent opening. The cryogen vent must meet all cryogen venting requirements (see *Magnet Room Venting Requirements*, 5850263-1EN).

2.2 System components



(Applies to all subsections within this section)

This system consists of the following components:

2.2.1 Magnet Room

1. 3.0T Magnet and Magnet Enclosure (MAG) and Vibroacoustic Damping Kit
2. Rear Pedestal (PED)
3. Patient Transport Table (PT)
4. Optional: Surgical suite table
5. Magnet Rundown Unit (MRU)



NOTE

An optional remote MRU may be located outside the Magnet Room.

2.2.2 Equipment Room

1. Main Disconnect Panel (MDP) (may be customer-supplied in exempt regions)
2. Power, Gradient, RF Cabinet (PGR)
3. Integrated Cooling Cabinet (ICC) (ICC refers to either ICC or pICC)
4. Penetration Cabinet (PEN)
5. Secondary Penetration Wall (SPW)
6. Magnet Monitor (MON)
7. Optional:
 - Magnetic Resonance Elastography (MRE)
 - Multi-Nuclear Spectroscopy (MNS)
 - MR Guided Focus Ultrasound (FUS)

2.2.3 Control Room


1. Operator Workspace equipment (OW)
2. Pneumatic Patient Alert System (PA1)
3. Optional:
 - Oxygen Monitor (OXY)
 - Physiological Acquisition Transceiver (PAT) Charging Station

2.2.4 Accessories

1. Patient accessories, including RF coils, phantoms, cushions, sponges, straps, and wedges
2. Gating accessories, including patient cardiac leads, peripheral gating probe, and respiratory bellows

2.2.5 System Overview

Figure 2-1 System overview (example)

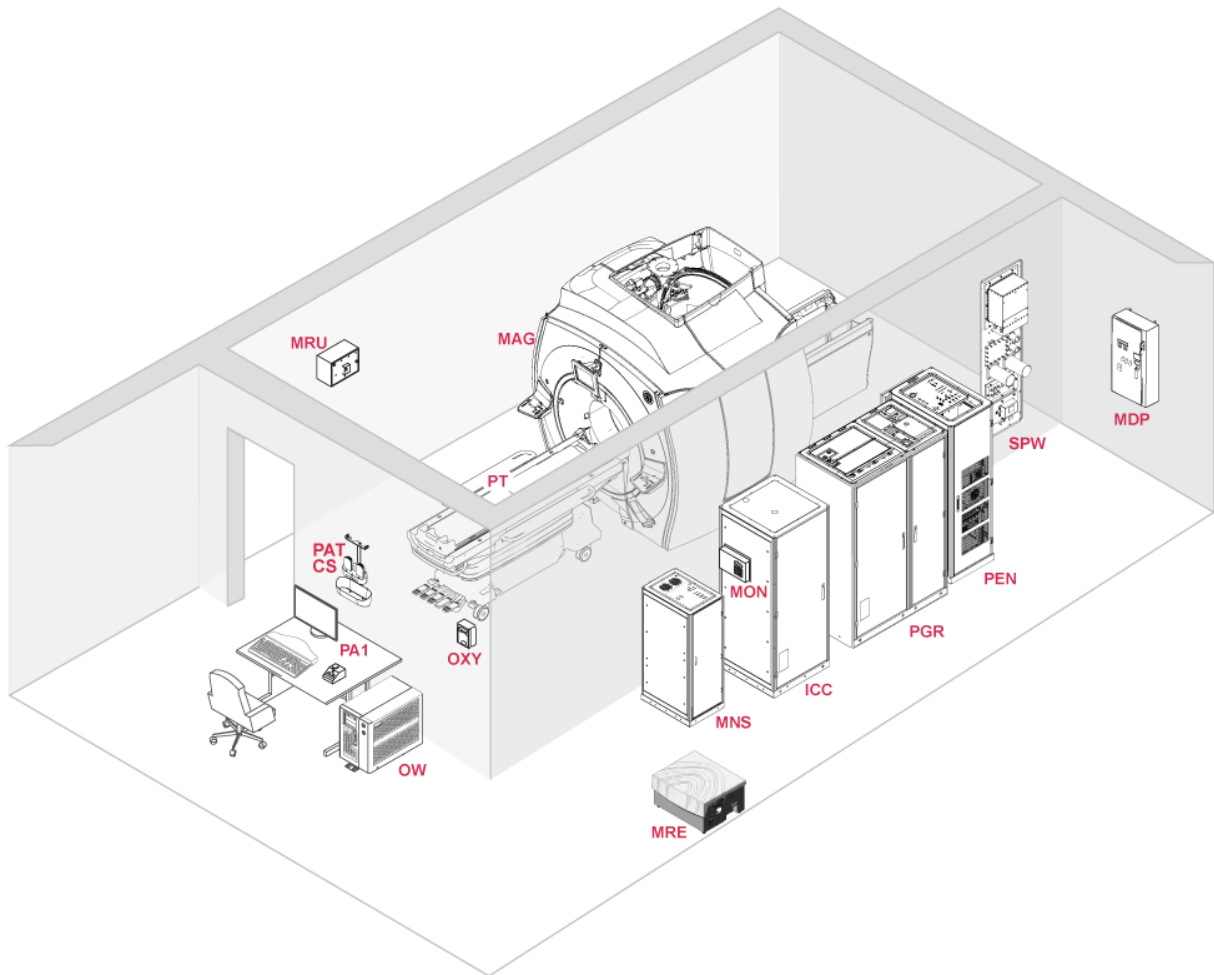

NOTE

MNS, MRE, PAT Charging Station and OXY shown above are optional components of the system.

2.3 MR Suite Minimum Room Size Requirements



Room dimensions shown in the table below are the minimum finished room space requirements to safely install and service the MR System. Minimum dimensions are for service only. Room size may increase due to the items listed below, which are not included in the minimum area dimensions:

1. Building code requirements (for example, exit routes, door placement, seismic mounting requirements, local and national electrical codes, and so on).
2. Equipment and Magnet Room evacuation routes to comply with facility emergency procedures.
3. System requirements, including cable run locations, cryogen venting, patient observation requirements, and penetration panel placements.
4. Penetration panel closet and all associated areas.
5. GE HealthCare optional equipment, such as MRE, MNS, accessories, and so on.
6. Non-GE HealthCare equipment options (such as additional AC or water cooling equipment in the Equipment Room).
7. Clinical workflow considerations.
8. Accessory storage. Refer to *Customer Site Storage Requirements*, 5182674 (available in the Customer Documentation Portal), or contact the GE HealthCare Project Manager of Installation (PMI) for any additional accessory storage requirements.
9. Magnetic field containment, for example, the 5 gauss line to the room. If fringe field containment is needed, see [2.6 MR Suite Magnetic Field Specifications on page 24](#).
10. The minimum service area shown must be kept clear of permanent or installed cabinetry, the MRU, the penetration closet, millwork, shelving, coil storage fixtures, furniture, and so on.


Table 2-1 Room Dimensions to Satisfy Minimum Service Area Requirements

Configuration	Equipment Room ¹			Magnet Room ²			Control Room		Total System Area m ² (ft ²)
	W x D mm (in.)	Area m ² (ft ²)	Ceiling Height mm (in.)	W x D mm (in.)	Area m ² (ft ²)	Finished Ceiling Height mm (in.)	W x D mm (in.)	Area m ² (ft ²)	
Minimum Room Size (2050 mm Scan Range)	2774 x 2728 (109.2 x 107.4)	7.6 (81.4)	2896 (114)	3645 x 6331 (143.5 x 253)	20.8 (224.1)	2500 (98.5)	1524 x 2134 (60 x 84)	3.2 (35)	31.6 (340.5)
With MR Guided Focused Ultrasound	Contact the Project Manager of Installation (PMI) for FUS site planning details.								
With MNS (m ² (ft ²))		Add: 1.4 (15)							

Table 2-1 Room Dimensions to Satisfy Minimum Service Area Requirements (Table continued)

Configuration	Equipment Room ¹			Magnet Room ²			Control Room		Total System Area m ² (ft ²)
	W x D mm (in.)	Area m ² (ft ²)	Ceiling Height mm (in.)	W x D mm (in.)	Area m ² (ft ²)	Finished Ceiling Height mm (in.)	W x D mm (in.)	Area m ² (ft ²)	
¹ See Figure 4-1 Typical Minimum Equipment Room with Service Clearances (with ICC) on page 90 or Figure 4-2 Typical Minimum Equipment Room with Service Clearances (with pICC) on page 91 for specific dimensions ² See Figure 2-2 Minimum Magnet Service Area (Top View) on page 19 for specific dimensions									


Figure 2-2 Minimum Magnet Service Area (Top View)

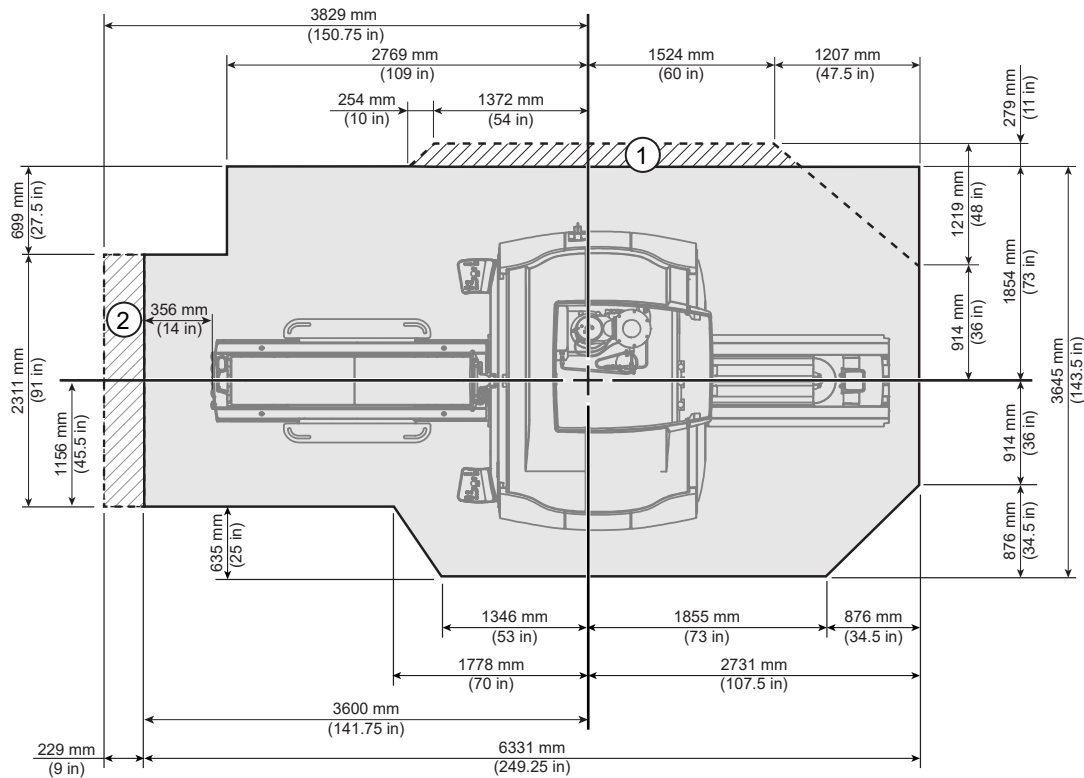

Item	Description	Item	Description
1	Alternate design tradeoff area	2	Recommended customer workflow area

The following requirements apply to the Magnet Room finished ceiling:

1. Ideal Magnet Room suspended ceiling height is 2667 mm (105 in.). Minimum Magnet Room suspended ceiling height is 2500 mm (98.5 in.). See [Figure 2-3 Area for Minimum Magnet Ceiling Height \(Top View\) on page 20](#).
2. The ceiling service area should be kept clear of overhead items, including soffits, HVAC, plumbing components, and brackets. Permanent or installed objects in this area may prevent or delay magnet service or operation.

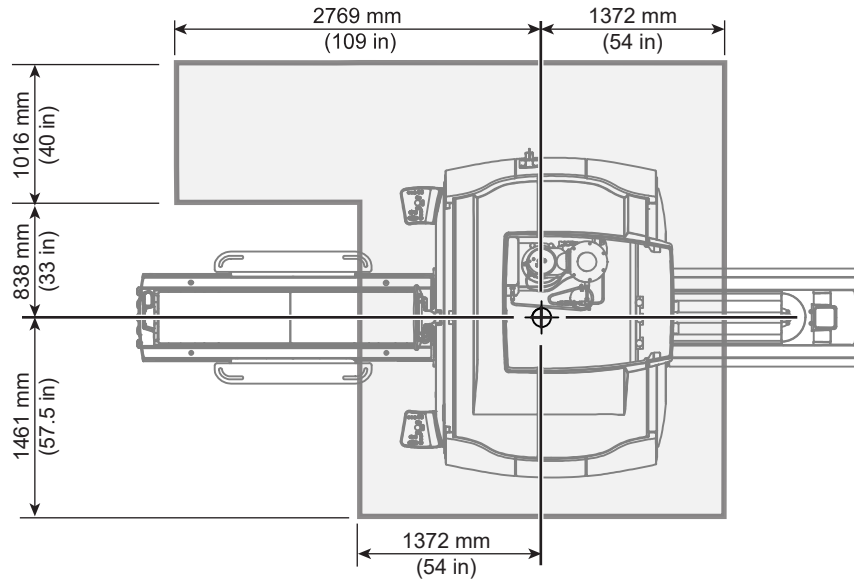Figure 2-3 Area for Minimum Magnet Ceiling Height (Top View)

NOTE
 The shaded area within solid lines indicates floor to ceiling minimum height of 2500 mm (98.5 in).

NOTE
 If the ceiling height is between **2500 mm (98.5 in.)** and **2667 mm (105 in.)**, the flexible main lead extension for low ceiling height (2.5M Low Ceiling Kit-Passive, M7000GM) is required for ramping the magnet. Contact the GE PMI and GE Service Field Engineer for further evaluation.

2.4 MR System Seismic Requirements



Contact the Project Manager of Installation with any questions.

1. The customer is responsible for seismic anchoring of GE HealthCare components.
2. Center of gravity, weight, physical dimensions, and attachment points are provided for seismic calculations. Refer to the specifications or illustrations for each component (see [Magnet Room Equipment Specifications on page 83](#), [Equipment Room on page 89](#), and [Control Room on page 108](#)).

2.5 Structure-borne Vibration Control Specifications



Structure-borne acoustic issues tend to occur at MR installations above the ground floor of the facility. Two options to mitigate structure-borne acoustic transmission are:


1. GE HealthCare provides a VibroAcoustic Damping kit (which must be surface mounted). Contact the GE HealthCare Project Manager of Installation for information.
2. The customer may design and implement a custom solution in addition to the VibroAcoustic Damping kit. See [Figure 2-4 Vibration Transmitted through VibroAcoustic Mat for UA Magnet on page 21](#) and [Figure 2-5 Vibration Transmitted through VibroAcoustic Mat for AR Magnet on page 22](#) for the plot of spectral vibration transmitted through the VibroAcoustic mat into the floor. If required, the customer should consult an acoustic engineer for a solution to further attenuate this transmitted vibration).



NOTE


The amount of vibration attenuation provided by the VibroAcoustic Damping kit will be site dependent.

Figure 2-4 Vibration Transmitted through VibroAcoustic Mat for UA Magnet

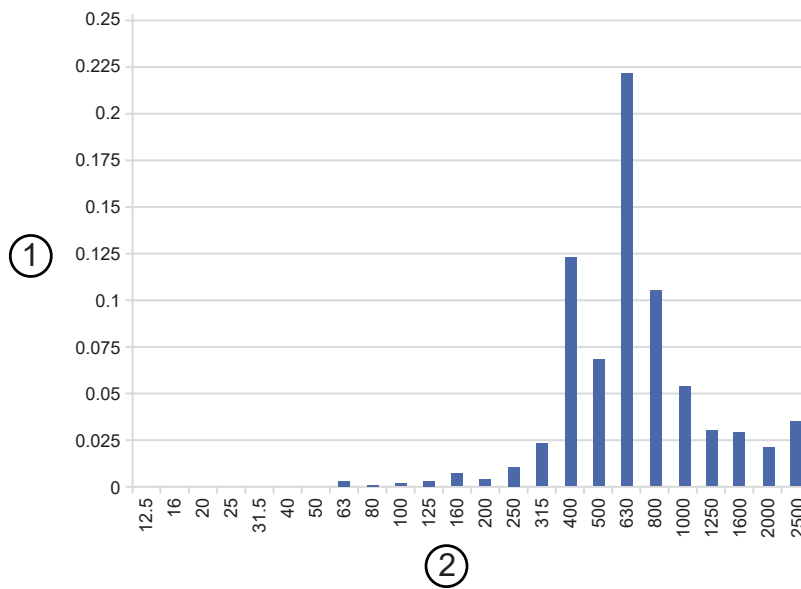


Item	Description
1	Acceleration [m/s ²]
2	1/3 Octave Frequency [Hz]

Figure 2-5 Vibration Transmitted through VibroAcoustic Mat for AR Magnet

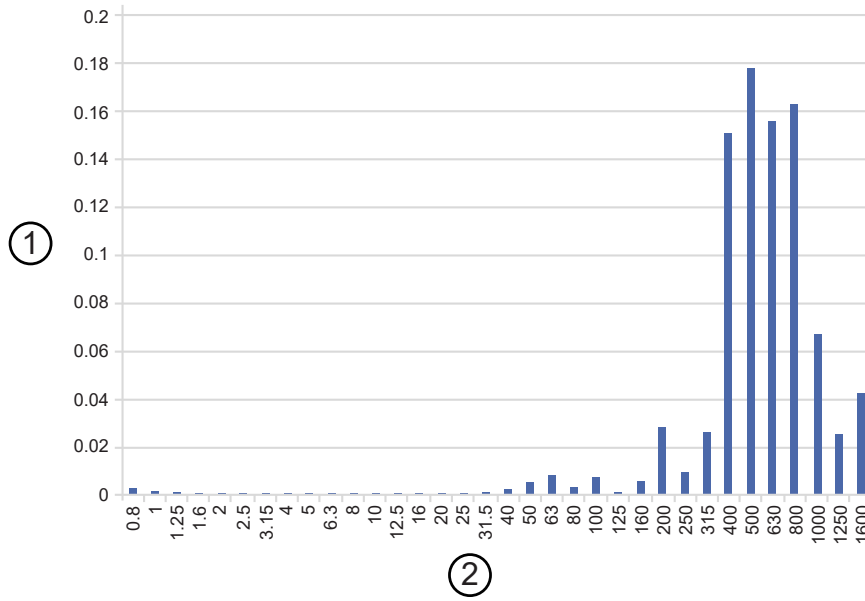

Item	Description
1	Acceleration [m/s ²]
2	1/3 Octave Frequency [Hz]

Figure 2-6 Low Frequency Magnet Floor Vibration (Vibration Amplitude at Each Foot) for AR and UA magnet

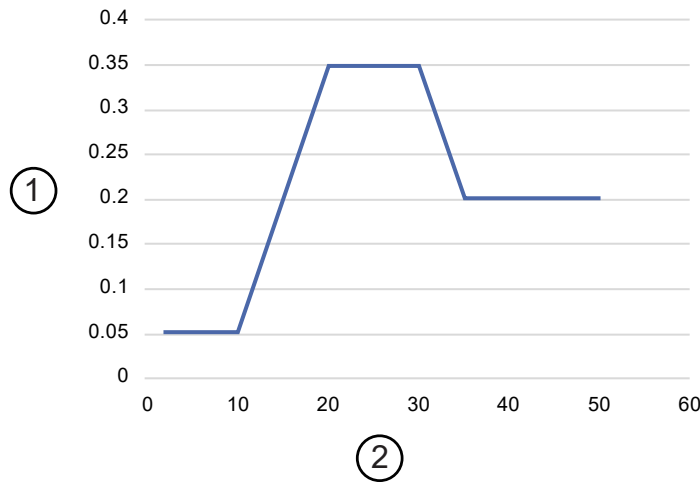

Item	Description
1	Vibration Amplitude (m/s ²)
2	Frequency (Hz)

Frequency (Hz)	Amplitude (m/s ²)
2	0.05
10	0.05
20	0.35
30	0.35
35	0.20
50	0.20


Low Frequency Magnet Floor Vibration Notes:



NOTE

1. Illustrations above define the potential vibration level that may pass into the customer site. [Figure 2-4 Vibration Transmitted through VibroAcoustic Mat for UA Magnet on page 21](#) and [Figure 2-5 Vibration Transmitted through VibroAcoustic Mat for AR Magnet on page 22](#) is the high frequency audible vibration. [Figure 2-6 Low Frequency Magnet Floor Vibration \(Vibration Amplitude at Each Foot\) for AR and UA magnet on page 23](#) is low frequency vibration that may dynamically displace the floor.
2. Vibration transfer may be the result of customer specific building construction as low levels of vibration transmit into the building through airborne and structure-borne paths. Customer MR clinicians recognize the vibration defined in the illustrations above is typically short bursts of vibration repeated multiple times as the scan progresses.
3. The customer should consider the impact of this vibration for the evaluation and design solution.

2.6 MR Suite Magnetic Field Specifications

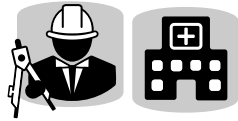

(Applies to all subsections within this section)

2.6.1 Magnetic Fringe Field

The following illustrations show the static magnet isogauss plot lines for the magnet. This information must be used to evaluate potential site interaction of GE HealthCare equipment with other non-GE HealthCare equipment, interaction with ferrous materials on the site, and to locate personnel and equipment within the site.


For UA magnet: The 0.5 mT (5G) line can expand to 6.0 m (19.7 ft.) axially and 4.5 m (14.8 ft.) radially for up to 8 seconds in the rare event of a quench.

For AR magnet: The 0.5mT (5G) line can expand to 6.5 m (21.33 ft.) axially and 3.5 m (11.48 ft.) radially for up to 2 seconds in the rare event of a quench.

The isogauss plots show an idealized magnetic field relative to magnet isocenter. The actual field strength can be affected by any of the following:


- Magnetic shielding
- Earth's magnetic field
- Other magnetic fields
- Stationary or moving metal

Figure 2-7 Magnetic Fringe Field Side View (UA series magnet)

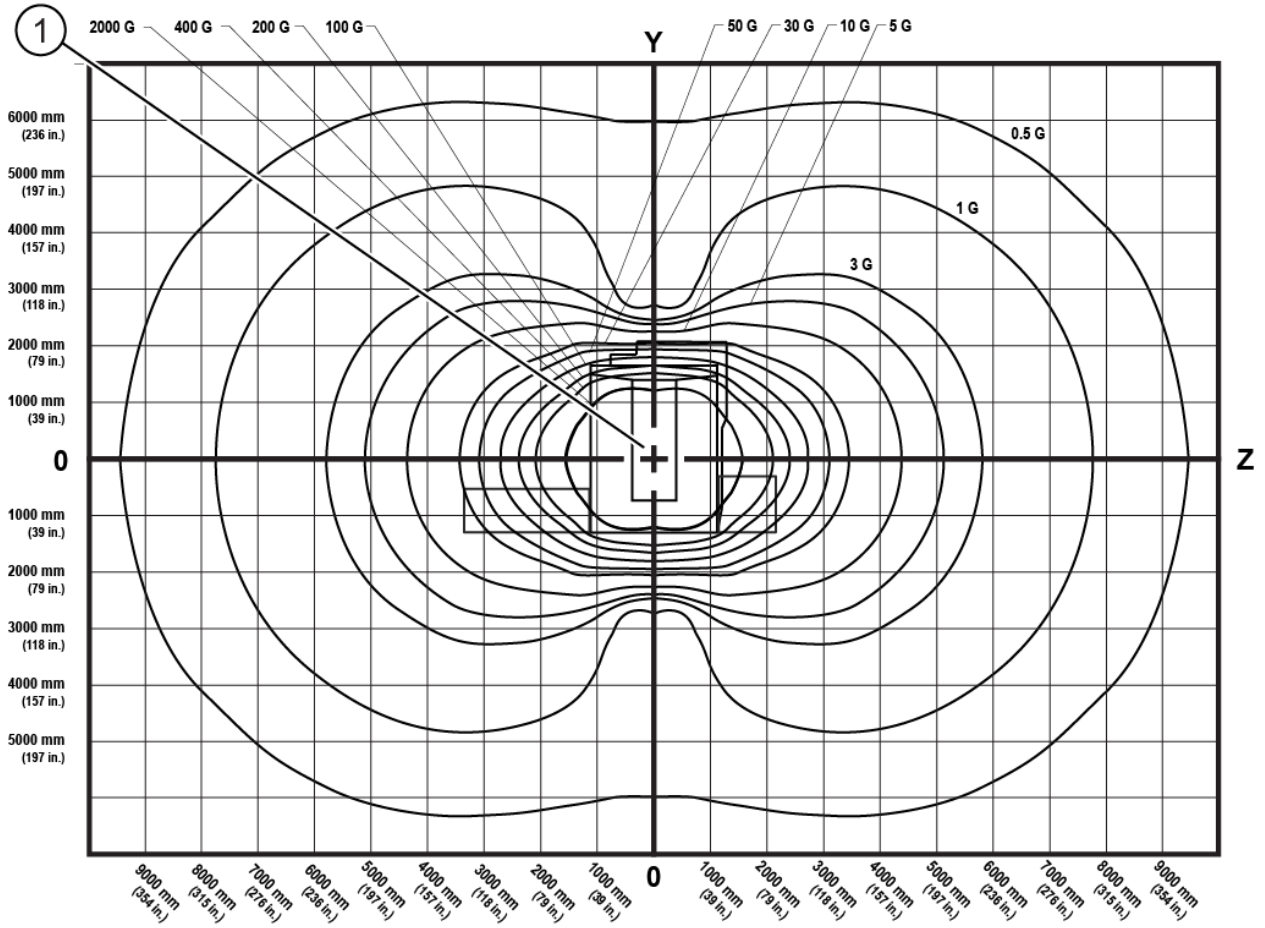


Item	Description
1	Magnet isocenter

Figure 2-8 Magnetic Fringe Field Top View (UA series magnet)

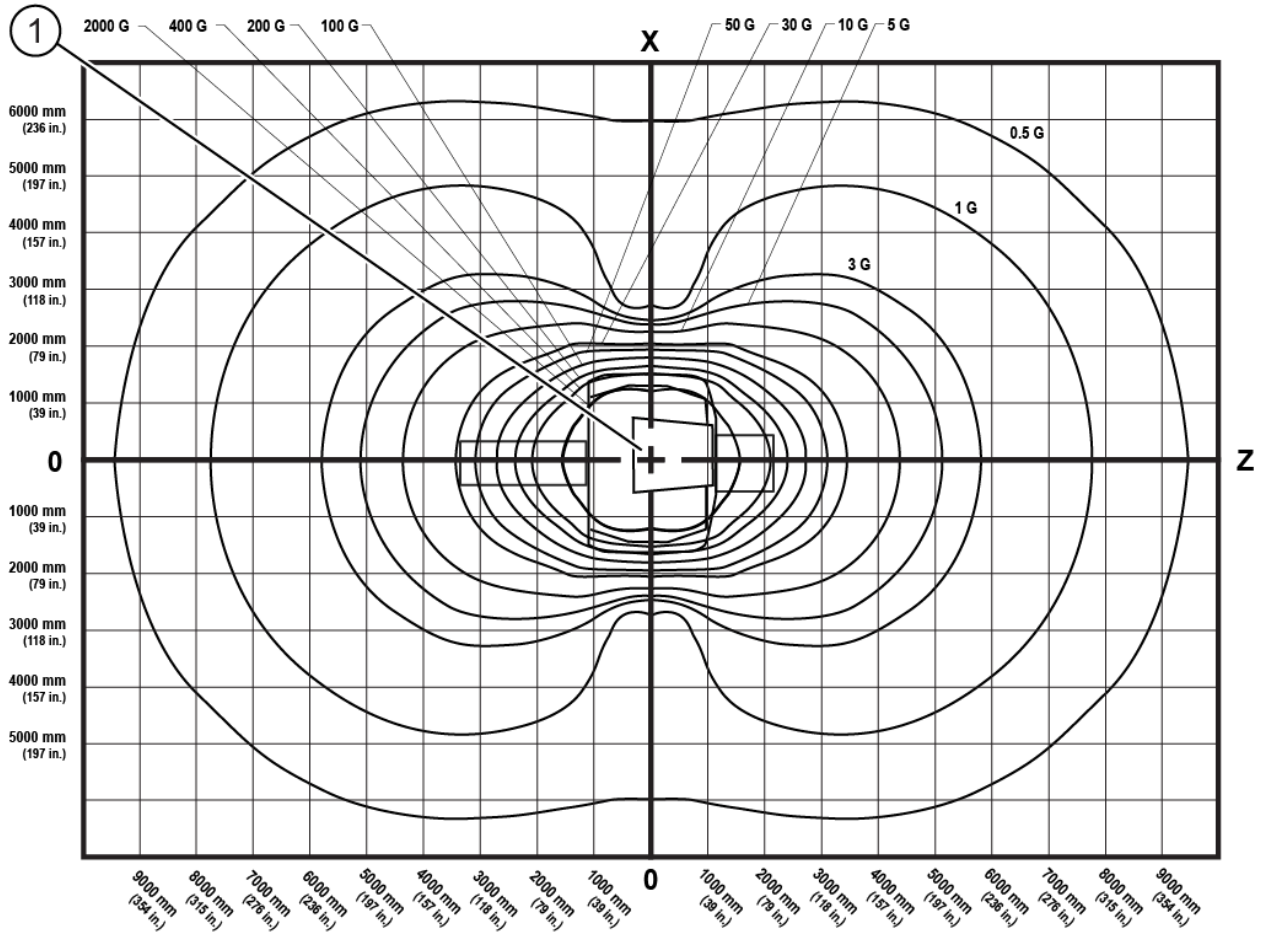


Item	Description
1	Magnet isocenter

Figure 2-9 Magnetic Fringe Field Front View (UA series magnet)

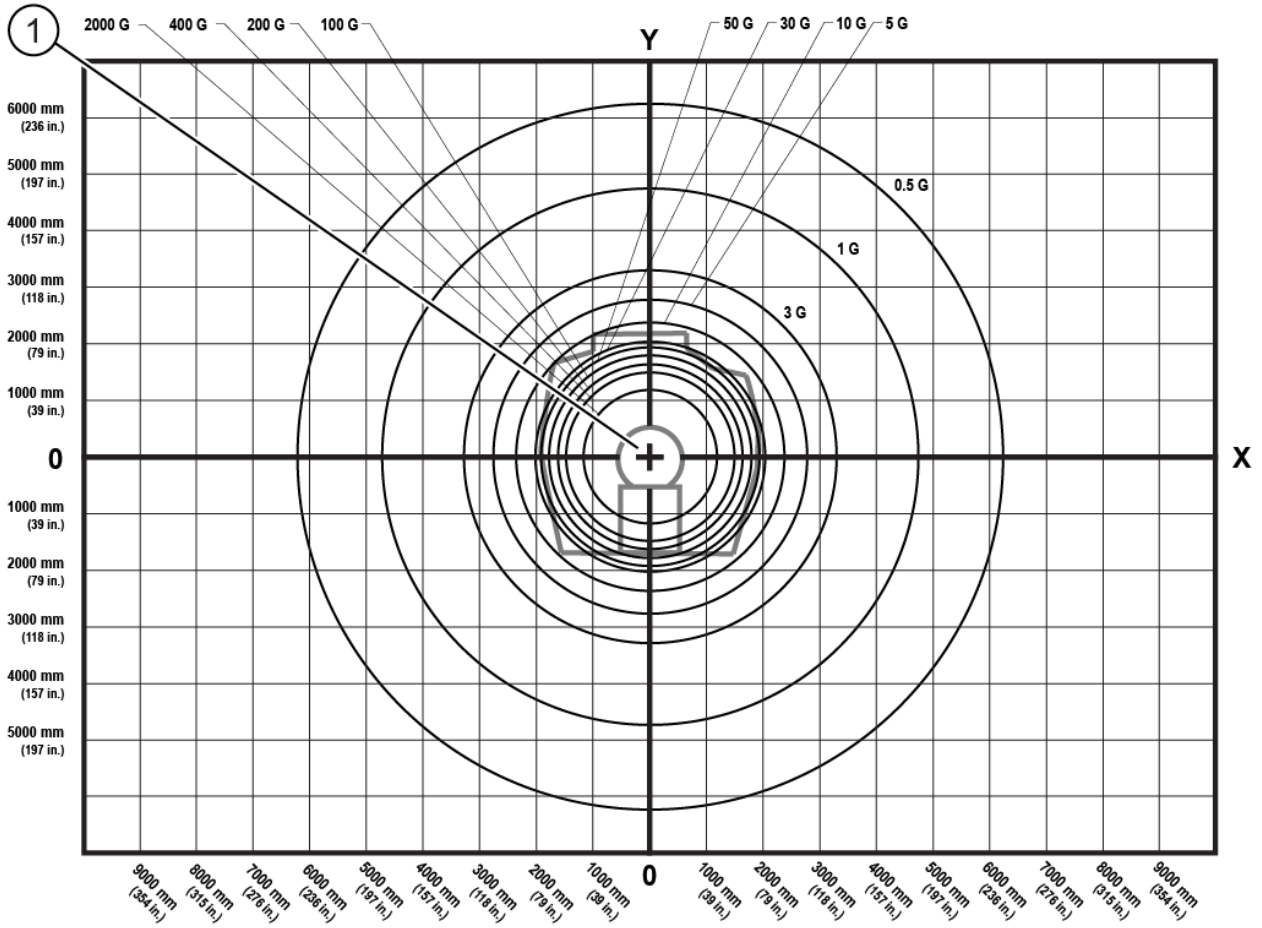


Item	Description
1	Magnet isocenter

Figure 2-10 Magnetic Fringe Field Side View (AR series magnet)

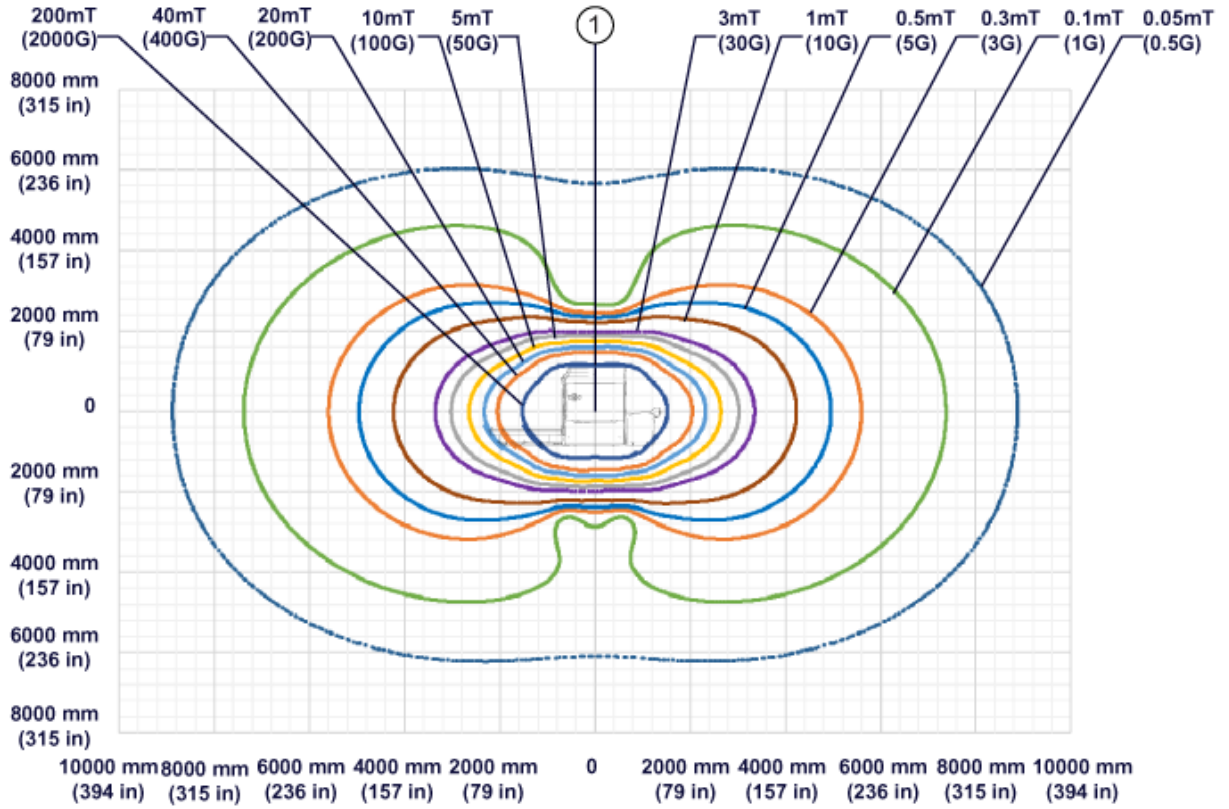


Item	Description
1	Magnet isocenter

Figure 2-11 Magnetic Fringe Field Top View (AR series magnet)

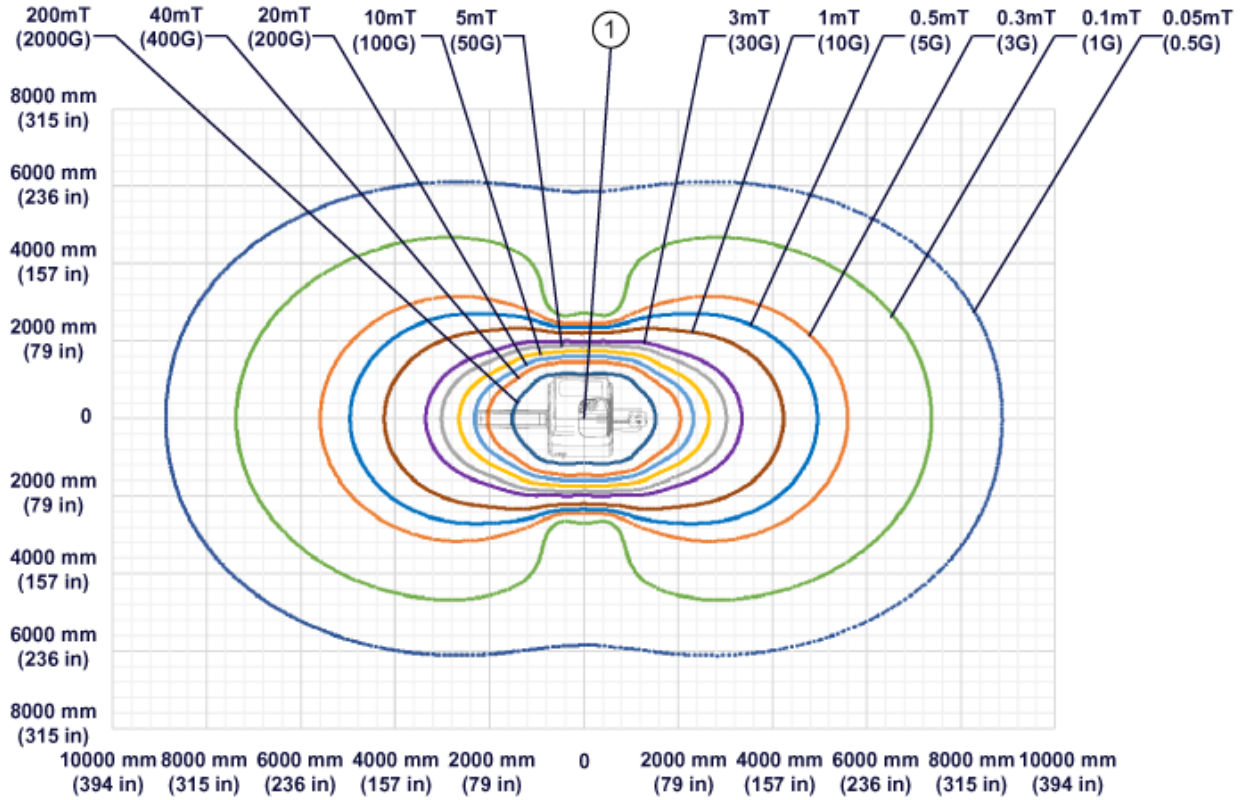

Item	Description
1	Magnet isocenter

Figure 2-12 Magnetic Fringe Field Front View (AR series magnet)

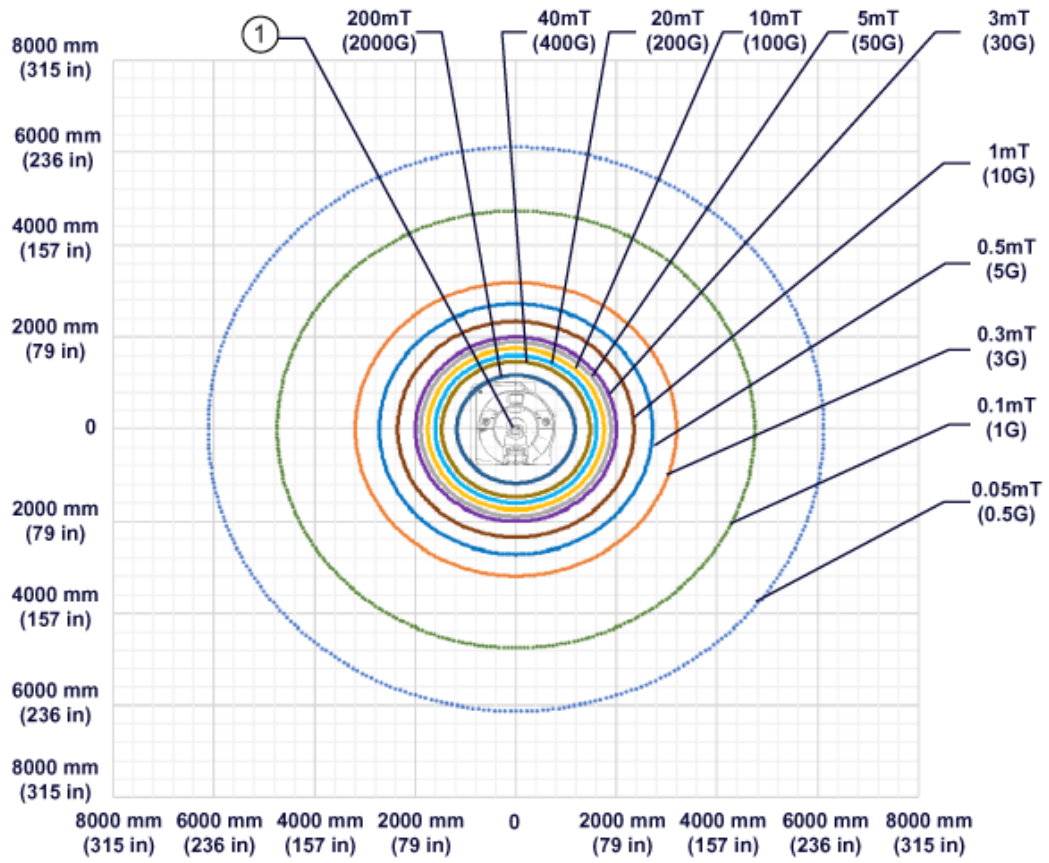


Item	Description
1	Magnet isocenter

2.6.2 Interference from Changing Magnetic Fields

Metal objects moving within the magnet sensitivity lines can produce a field disturbance during clinical imaging. If the metal object is moving it will produce a fluctuating dipole type of field which causes image artifacts. As an example, a car driven inside the moving metal line will act as a dipole and produce a time varying field which changes the magnet's main field during the scanning. The same vehicle may park within the moving metal line and remain parked during clinical scanning without impact to the main field.

Figure 2-13 Magnet Moving Metal Sensitivity Line Plot (Side View, UA series magnet)

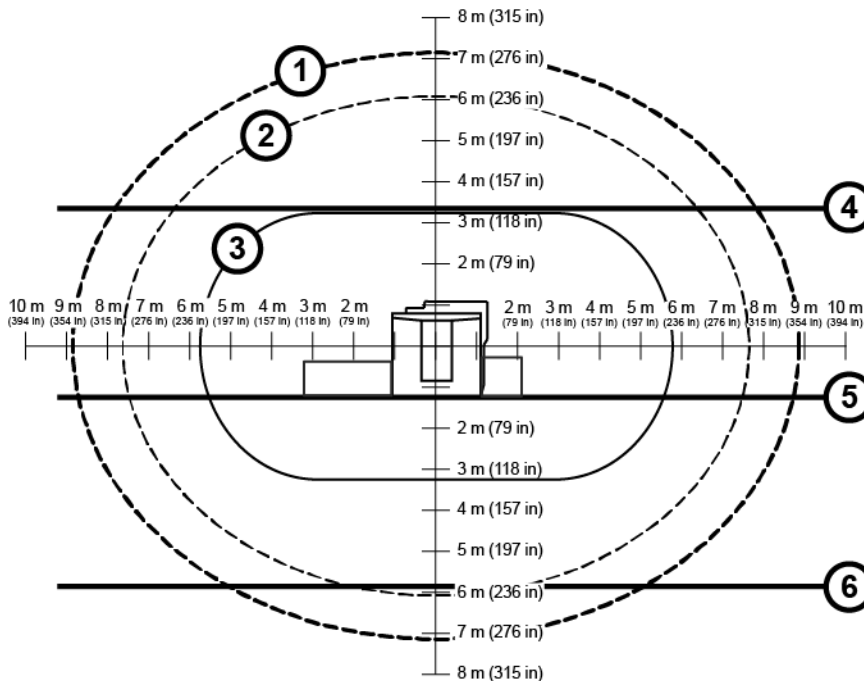


Item	Description	Item	Description
1	Trucks, Buses	4	Floor Above
2	Cars, Pickups, Vans, Ambulances	5	Magnet Room Floor
3	3 Gauss line	6	Floor Below



NOTE

The magnet isocenter, which is 1070 mm (42.1 in.) above the floor, is the origin of both the x-axis and the y-axis.

Figure 2-14 Magnet Moving Metal Sensitivity Line Plot (Side View, AR series magnet)

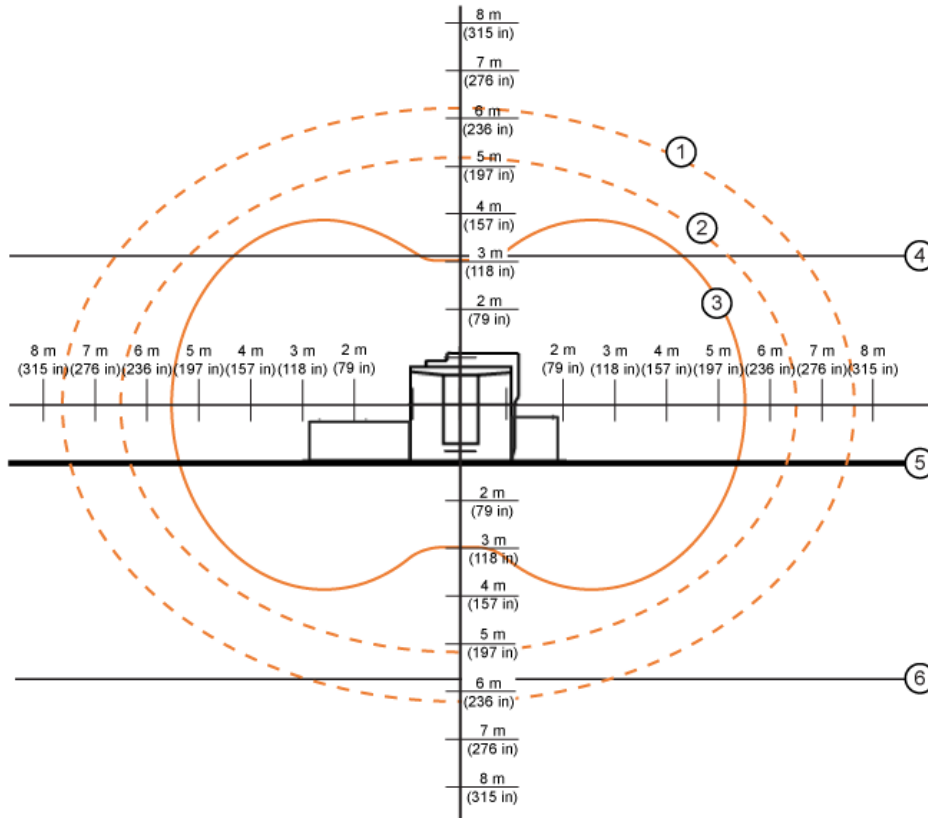


Item	Description	Item	Description
1	Trucks, Buses	4	Floor Above
2	Cars, Pickups, Vans, Ambulances	5	Magnet Room Floor
3	3 Gauss line	6	Floor Below



NOTE


The magnet isocenter, which is 1070 mm (42.1 in.) above the floor, is the origin of both the x-axis and the y-axis.

Figure 2-15 Magnet Moving Metal Sensitivity Line Plot (Top View, UA series magnet)

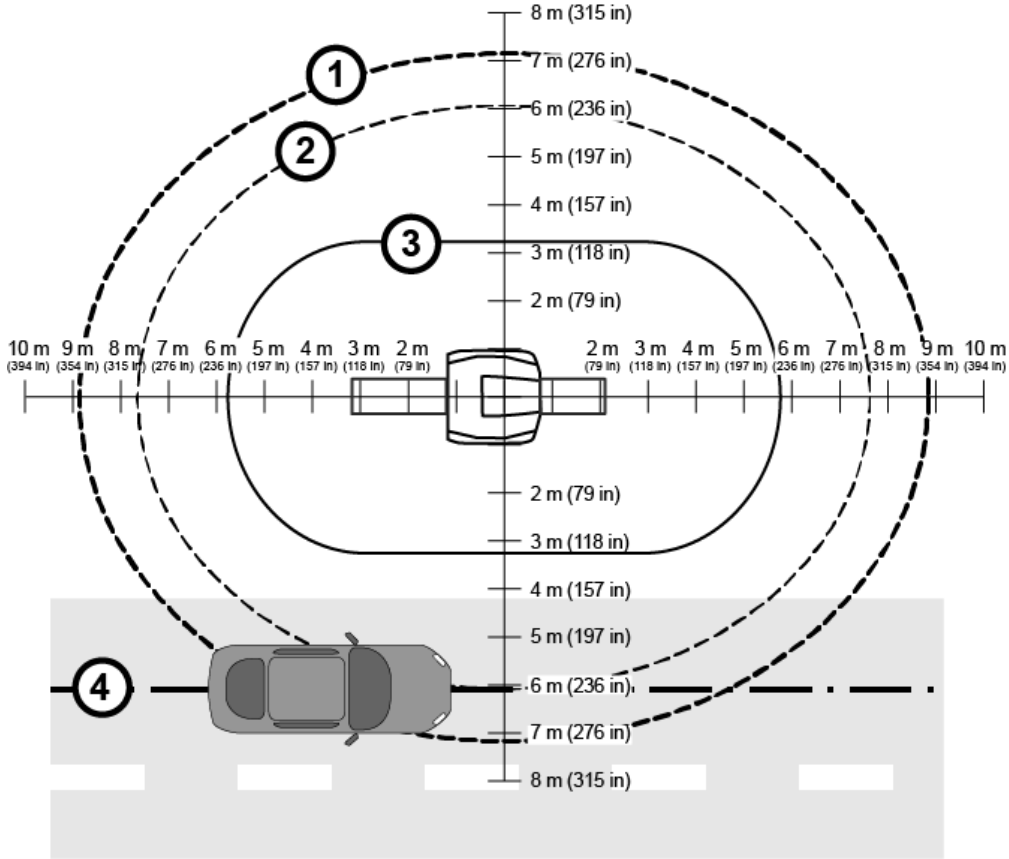

Item	Description	Item	Description
1	Trucks, Buses	3	3 Gauss line
2	Cars, Pickups, Vans, Ambulances	4	Center of Driving Lane

Figure 2-16 Magnet Moving Metal Sensitivity Line Plot (Top View, AR series magnet)

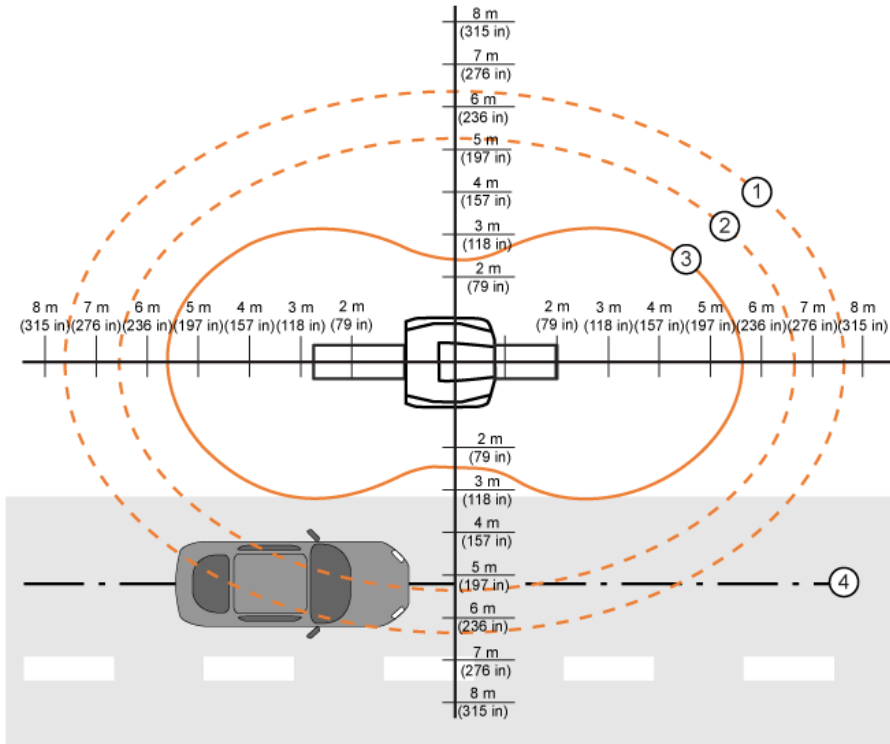

Item	Description	Item	Description
1	Trucks, Buses	3	3 Gauss line
2	Cars, Pickups, Vans, Ambulances	4	Center of Driving Lane

Table 2-2 Magnet Moving Metal Requirements

Metal Objects Category	Definition Of Distance Location	Magnet Minimum Distance Radial X Axial ¹ m (ft.)
Objects 45.36 - 181.44 kg (100 - 400 lb.)	Distance from isocenter radial x axial	0.3 mT (3 G) line
Cars, Minivans, Vans, Pickup Trucks, Ambulances	Distance from isocenter measured to center of driving or parking lane radial x axial	6.05 x 7.65 (20 x 25) (UA series magnet) 5.29 x 6.58 (17.4 x 21.6) (AR series magnet)
Bus, Trucks (Utility, Dump, Semi)	Distance from isocenter measured to center of driving or parking lane radial x axial	7.10 x 8.90 (23.3 x 29.2) (UA series magnet) 6.17 x 7.69 (20.3 x 25.2) (AR series magnet)
Objects > 181.44 kg (400 lb.), Elevators, Trains, Subways	Place a directional probe (for example, flux gate sensor) at isocenter of proposed magnet location aligned along the Z-axis. Measure peak-to-peak magnetic field change (DC).	Z-axis EMI disturbances must not exceed 560 nanotesla (5.6 milligauss) peak-to-peak. (UA series magnet) Z-axis EMI disturbances must not exceed 520 nanotesla (5.2 milligauss) peak-to-peak. (AR series magnet)


Notes:

1. Radial distances are magnet X and Y axis. Axial distances are magnet Z axis.

2.6.3 Electrical Current

1. Electrical current in high voltage power lines, transformers, motors, or generators near the magnet may affect magnetic field homogeneity.
2. Magnetic field interference at 50 or 60 Hz must not exceed 1.622 μT (16.22 mG) or 1.710 μT (17.10 mG) RMS respectively at the magnet location (see [Figure 2-17 Magnet Allowable Milligauss vs. Line Frequency for AC Equipment on page 35](#)).
3. The following equation can be used as a general guide in determining allowable current in feeder lines at a given distance from the magnet isocenter:
 - 3.1. For 3.0T Magnet: $I = (8.55X^2)/S$, where:
 - 3.2. I = Maximum allowable RMS single phase current (in amps) or maximum allowable RMS line current (in amps) in three phase feeder lines
 - 3.3. S = Separation (in meters) between single phase conductors or greatest separation between three phase conductors
 - 3.4. X = Minimum distance (in meters) from the feeder lines to isocenter of the magnet

Figure 2-17 Magnet Allowable Milligauss vs. Line Frequency for AC Equipment

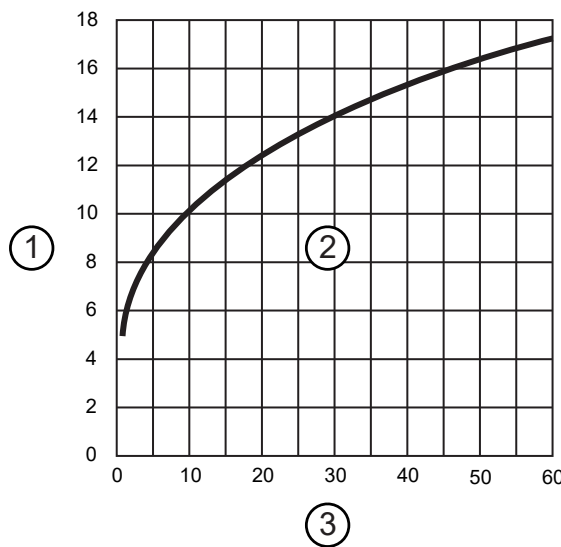

Item	Description	Item	Description
1	Milligauss	3	Excitation Frequency (Hz)
2	No Impact on Imaging		

Refer to [Sample Calculation AC Power Equipment Minimum Distance on page 142](#) for additional examples.

2.6.4 Non-MR System Equipment Sensitivity to Magnetic Fields

Site plans must include consideration for magnetic field interaction with all customer equipment.

This section lists equipment known to be sensitive to high magnetic fields.

Use the table for reference only. The Tesla (Gauss) limits in the table are approximate for that type of equipment. Refer to OEM manuals for the equipment at your site to determine the actual Tesla (Gauss) limits.

Table 2-3 Magnetic Proximity Limits (For Reference Only)

mT (Gauss) Limit	Equipment		
0.05mT (0.5 G)	Nuclear camera		
0.1mT (1 G)	Positron Emission Tomography scanner	Video display (tube)	
	Linear Accelerator	CT scanner	
	Cyclotrons	Ultrasound	
	Accurate measuring scale	Lithotripter	
	Analog image intensifiers	Electron microscope	
	Bone Densitometers		
0.3mT (3 G)	Power transformers	Main electrical distribution transformers	
0.5mT (5 G)	Cardiac pacemakers	Biostimulation devices	
	Neurostimulators		
1mT (10 G)	Magnetic computer media	Telephone switching stations	
	Hard copy imagers	Water cooling equipment	
	Line printers	HVAC equipment	
	Video Cassette Recorder (VCR)	Major mechanical equipment room	
	Film processor	Credit cards, watches, and clocks	
	X-ray tubes		
	Large steel equipment, including:		
	Emergency generators	Air conditioning equipment	
	Commercial laundry equipment	Fuel storage tanks	
	Food preparation area	Motors greater than 5 horsepower	
5mT (50 G)	Metal detector for screening	Telephones	
	LCD panels		
No Limit	Digital Detectors		

2.7 Multiple MR System Requirements




(Applies to all subsections within this section)

2.7.1 Multiple Magnets

When installing multiple magnets, the 0.3 mT (3 G) lines must not intersect or the magnets will be interactive. Contact the GE HealthCare Project Manager of Installation (PMI) for any questions regarding magnetic field interaction.

Magnet Rooms cannot share walls.

Figure 2-18 Two Magnet Installation (No Interaction)

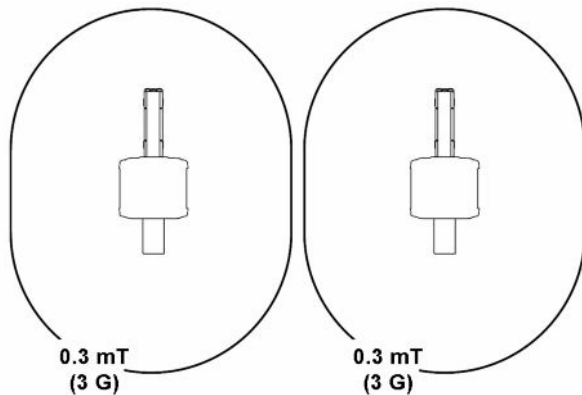

2.7.2 Shared Equipment Rooms

When installing multiple MR Systems in a shared Equipment Room, of the same field strength, the following conditions must be met:

1. Refer to [Figure 2-19 Cabinet Separations in Shared Equipment Room \(Top View\) on page 38](#) for cabinet location examples.

The RF cabinet of the existing installed system must be separated from the RF bay of the PGR by at least 2000 mm (79 in.).

This requirement does not apply if the existing 3.0T system is one of the following: SIGNA Pioneer, SIGNA Hero, SIGNA Performer, SIGNA Pilot, Discovery MR750, Discovery MR750w, SIGNA Architect, SIGNA Premier, or SIGNA PET/MR.

2. Penetration panels must be separated by at least 3000 mm (118 in).
3. Cables from different MR Systems must not be routed together.
4. Two systems cannot share common power or ground feeds.



NOTE

* The 2000 mm (79 in.) requirement does not apply if the existing system is one of those listed in the requirements in item 1 above.

Figure 2-19 Cabinet Separations in Shared Equipment Room (Top View)


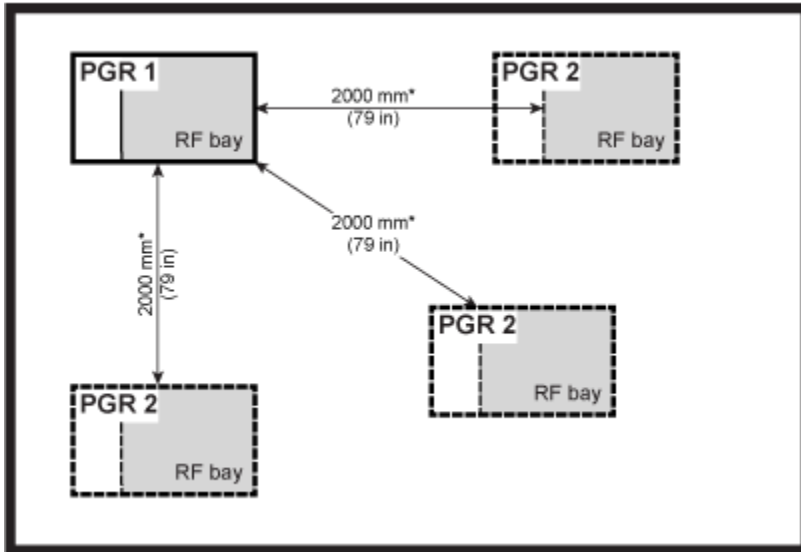


Figure 2-20 Panel Separations in Shared Equipment Room (Top View)



2.8 MR Suite Temperature and Humidity

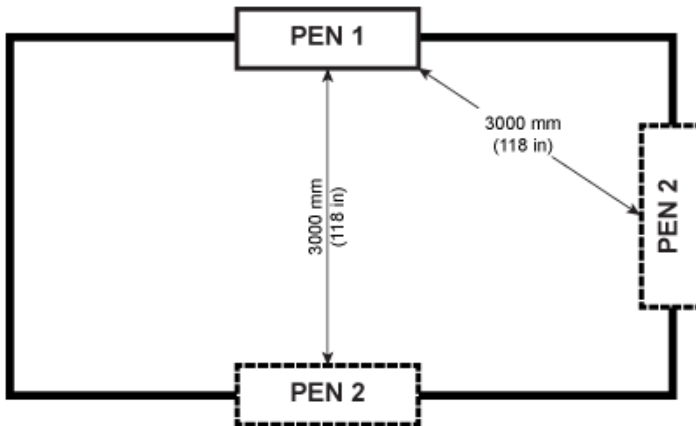

(Applies to all subsections within this section)

This section provides temperature and humidity requirements for the MR suite.



NOTE

Make sure the HVAC system has the correct capacity for the room size, equipment heat output, and environmental conditions to maintain correct temperature and humidity for the protection of the patient.

Specific construction requirements for each room can be found in the following chapters:

- Magnet Room
- Equipment Room
- Control Room

2.8.1 Temperature and Humidity Requirements

1. The customer is responsible for HVAC system design, purchase, and installation.
2. The temperature and humidity requirements must not be exceeded at any point during the day (both working or non-working hours).
3. A separate thermostat must be provided for the Magnet Room.

Table 2-4 Room Temperature and Humidity Requirements

Room	Temperature		Humidity	
	Range °C (°F)	Change °C/Hr (°F/Hr) ¹	Range %RH	Change %RH/Hr ²
Equipment Room (at Inlet to Equipment)	15-32 (59-89.6) ³	3 (5)	30-70	5
Magnet Room	15-21 (59-69.8)	3 (5)	30-60	5
Operator Room	15-32 (59-89.6)	3 (5)	30-70	5
Notes:				
1. Operating temperature gradient limits shall be between -3°C/Hr (-5°F/Hr) and 3°C/Hr (5°F/Hr), when averaged over 1 hour.				
2. Operating humidity gradient limits shall be between -5% RH/hour and 5% RH/hour, when averaged over 1 hour.				
3. Maximum ambient temperature is derated by 1°C per 300 m above 2000 m (not to exceed 2600 m).				

2.8.2 Equipment Heat Output Specifications

This section details the heat output for specific components. These heat outputs define the minimum, maximum and an assumed average condition over a 12-hour period. Actual heat output and room temperature may vary due to environmental factors, room insulation, clinical usage, and any non-GE HealthCare equipment used in the MR suite. Also, due to large variations in heat loads, the HVAC system may require unloaders, hot gas bypass, and reheat to maintain humidity levels.

Table 2-5 System Heat Output for Air Cooling

Component	Magnet Room W (BTU/hr)			Equipment Room W (BTU/hr)			Control Room W (BTU/hr)		
	Maximum	Average	Idle	Maximum	Average	Idle	Maximum	Average	Idle
Magnet (MAG) and Patient Transport Table (PT)	2400 (8189)	1200 (4095)	561 (1915)						
Penetration Panel Cabinet (PEN)	300 (1024)	150 (512)	0	3135 (10697)	1568 (5349)	1568 (5349)			
Secondary Penetration Wall (SPW)	0			0					
Main Disconnect Panel (MDP)				264 (901)	132 (450)	132 (450)			
Power Gradient RF Cabinet (PGR)				6137 (20940)	3068 (10470)	2500 (8530)			
Cryocooler Compressor (CRY)				500 (1706)	500 (1706)	500 (1706)			
Integrated Cooling Cabinet (ICC)				1000 (3410)	600 (2046)	250 (853)			
Magnet Monitor (MON)				240 (819)	240 (819)	240 (819)			
Operator Workspace Equipment (OW)							1450 (4947)		

Table 2-6 System Options Heat Output for Air Cooling

Component	Magnet Room W (BTU/hr)			Equipment Room W (BTU/hr)			Control Room W (BTU/hr)		
	Maximum	Average	Idle	Maximum	Average	Idle	Maximum	Average	Idle
Multi-Nuclear Spectroscopy Cabinet (MNS)				7205 (24584)					
MR Elastography (MRE)				141 (480)					

2.9 Facility Coolant Requirements



(Applies to all subsections within this section)



IMPORTANT

Equipment Failure. A continuous supply of facility liquid coolant to the Integrated Cooling Cabinet (ICC) is required at all times for correct system operation. Failure to provide liquid coolant with the requirements listed in this section may cause equipment failure.

2.9.1 Integrated Cooling Cabinet (ICC) Coolant Requirements

1. The facility must provide an uninterrupted supply of liquid coolant to the Integrated Cooling Cabinet (ICC) at magnet delivery. Coolant circuit must be operational at magnet delivery.
2. The facility must provide pipe/hose, filter, and connectors to the ICC.
3. The customer must provide and install an in-line flow meter on either the supply or return facility coolant hose. The flow meter must be capable of visually displaying volumetric flow between 76 and 189 L/min (20 and 50 GPM) and configured for the properties of the cooling fluid in use within the ICC coolant circuit.


Table 2-7 Facility Liquid Coolant Requirements

Parameter	Requirements
Availability	Continuous
Antifreeze or treated process water	No more than 50% propylene glycol-water (PGW) or ethylene glycol-water (EGW)
Minimum Flow	114 L/min (30 GPM)
Maximum Flow	132 L/min (35 GPM)
Maximum Pressure Drop in ICC at Minimum Flow	2.2 bar (32.0 psi) with 50% propylene glycol-water; 1060 kg/m ³ density 1.4 bar (20.3 psi) with pure water; 994 kg/m ³ density
Maximum Pressure Drop in ICC at Maximum Flow	2.9 bar (42 psi) with 50% propylene glycol-water; 1060 kg/m ³ density 1.9 bar (27.6 psi) with pure water; 994 kg/m ³ density
Temperature rise at Minimum Flow	10.4°C (18.7°F) with 50% propylene glycol-water; 3346 J/(kg K) specific heat; 1060 kg/m ³ density; 70 kW heat
Temperature rise at Maximum Flow	9.0°C (16.1°F) with 50% propylene glycol-water; 3346 J/(kg K) specific heat; 1060 kg/m ³ density; 70 kW heat
Maximum Inlet Pressure to ICC	6 bar (87 psi)
Chiller Size	Minimum 70 kW
Condensation Protection	Facility Plumbing to the ICC must be properly routed and insulated to prevent equipment damage or safety hazards
Minimum Continuous Heat Load	7.5 kW
Inlet Temperature	5 to 15°C (41 to 59°F) measured at the inlet to the ICC

Table 2-7 Facility Liquid Coolant Requirements (Table continued)


Parameter	Requirements
Customer supplied feeder hose (from main water supply to ICC)	38.1 mm (1.5 in.) minimum hose inside diameter
Hose connections to the ICC	38.1 mm (1.5 in.) male NPT

Figure 2-21 Allowable Facility Water Temperature and Flow

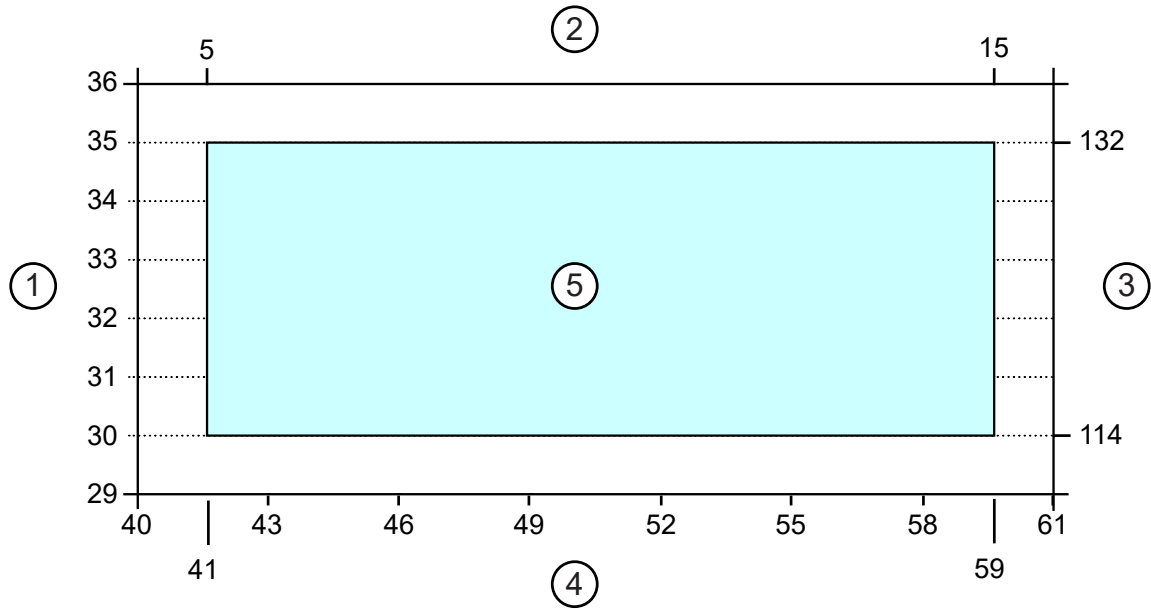


Item	Description	Item	Description
1	Flow [gpm]	4	Temperature [°F]
2	Temperature [°C]	5	0-50% Propylene or Ethylene Glycol
3	Flow [L/min]		

Figure 2-22 Pressure Drop through ICC (50% PGW, 1060 kg/m³)

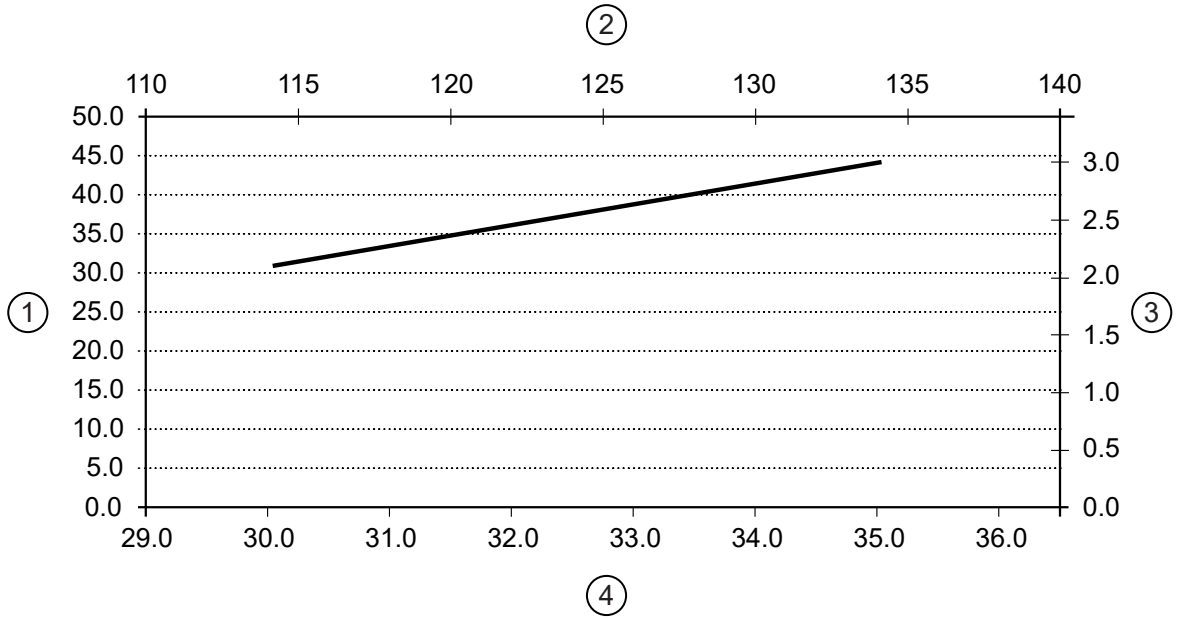

Item	Description	Item	Description
1	Pressure [psig]	3	Pressure [bar]
2	Flowrate [L/min]	4	Flowrate [gpm]

Figure 2-23 Temperature Rise through ICC (50% PGW, 3346 J/kg-K, 1060 kg/m³)

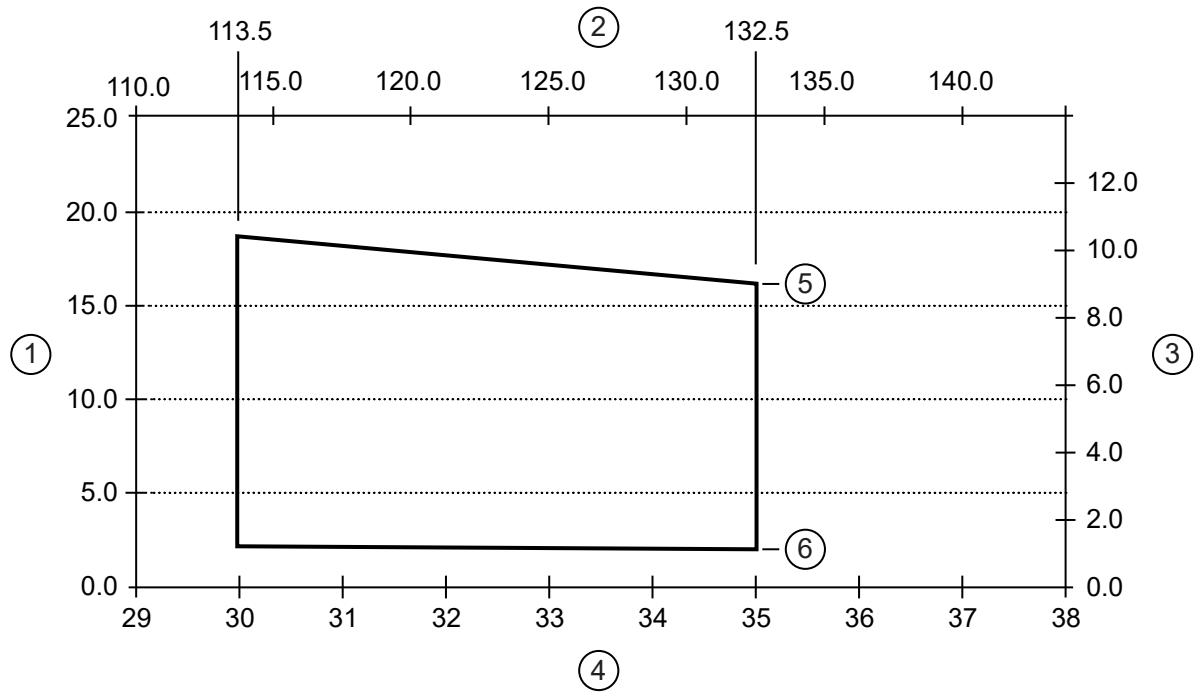

Item	Description	Item	Description
1	Temperature Rise [°F]	4	Flowrate [gpm]
2	Flowrate [L/min]	5	Max Heat Load: 70kW
3	Temperature Rise [°C]	6	Min Heat Load: 7.5kW


Table 2-8 Facility Water Quality Requirements

Parameter	Requirement
pH Value	6.5 to 8.2 at 25 °C (77 °F)
Electrical Conductivity	< 0.8 mmho/cm
Chloride Ion	< 200 ppm
Sulfate Ion	< 200 ppm
M-Alkalinity	< 100 ppm
Total Hardness	< 200 ppm
Calcium Hardness	< 150 ppm
Ionic Silica	< 50 ppm
Iron	< 1.0 ppm
Copper	< 0.3 ppm
Sulfide Ion	None, not detectable
Ammonium Ion	< 1.0 ppm
Residual Chlorine	< 0.3 ppm
Free Carbon Dioxide	< 4.0 ppm

Table 2-8 Facility Water Quality Requirements (Table continued)

Parameter	Requirement
Stability Index	6.0 to 7.0
Suspended Matter	< 10 ppm
Particle Size	< 100 micron (with field changeable filter)

Figure 2-24 MR System Water Cooling Block Diagram

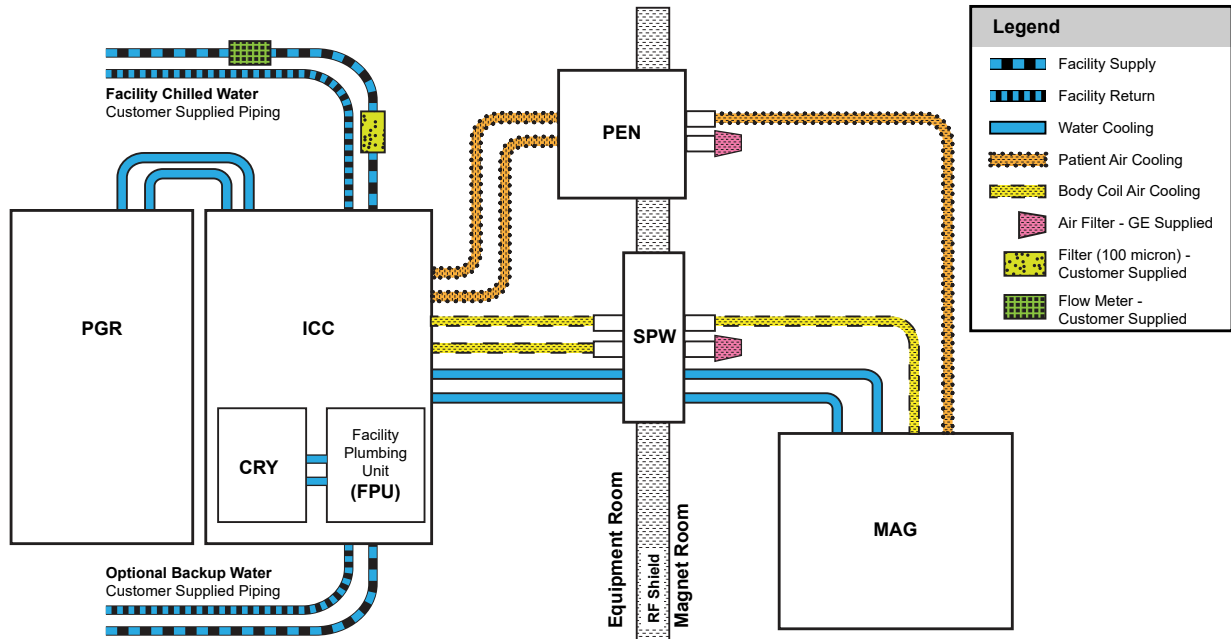

2.9.2 Requirements for Emergency Backup Facility Coolant (Optional)

The customer must balance the cost of cryogenes and local controls with the cost of emergency backup facility coolant. There are two options for emergency backup, either total ICC backup or Cryocooler Compressor backup.

The facility is responsible for either option.

The facility may provide an optional emergency backup coolant supply in one of the following configurations:


1. Total ICC backup:
 - 1.1. The facility is responsible for the connection of all hoses of a backup system.
 - 1.2. Coolant must meet all other ICC coolant requirements listed in [2.9.1 Integrated Cooling Cabinet \(ICC\) Coolant Requirements on page 41](#).

2. Cryocooler Compressor backup only:

Backup coolant may be routed to the ICC backup water port for the Cryocooler compressor at the location indicated in [Figure 2-24 MR System Water Cooling Block Diagram on page 45](#) with the following requirements:

- 2.1. The emergency coolant supply must be drained to the facility through customer supplied 12.5 mm (0.5 in.) hose and must not back-feed to the ICC.
- 2.2. Coolant must meet all other ICC coolant requirements listed in [2.9.1 Integrated Cooling Cabinet \(ICC\) Coolant Requirements on page 41](#).
- 2.3. The facility is responsible for coolant, 12.5 mm (0.5 in.) pipe/hose, filters, and connectors to supply the coolant to the ICC emergency facility coolant supply line (12.5 mm (0.5 in.) hose with a female threaded connector). The emergency coolant supply line is available from the Facility Plumbing Unit on the ICC. The hose comes out from the bottom rear of the ICC or the top of the pICC.
- 2.4. The charts below show the coolant flow rate and temperature requirements for the Cryocooler Compressor:

Figure 2-25 Cryocooler Water Cooling Requirements

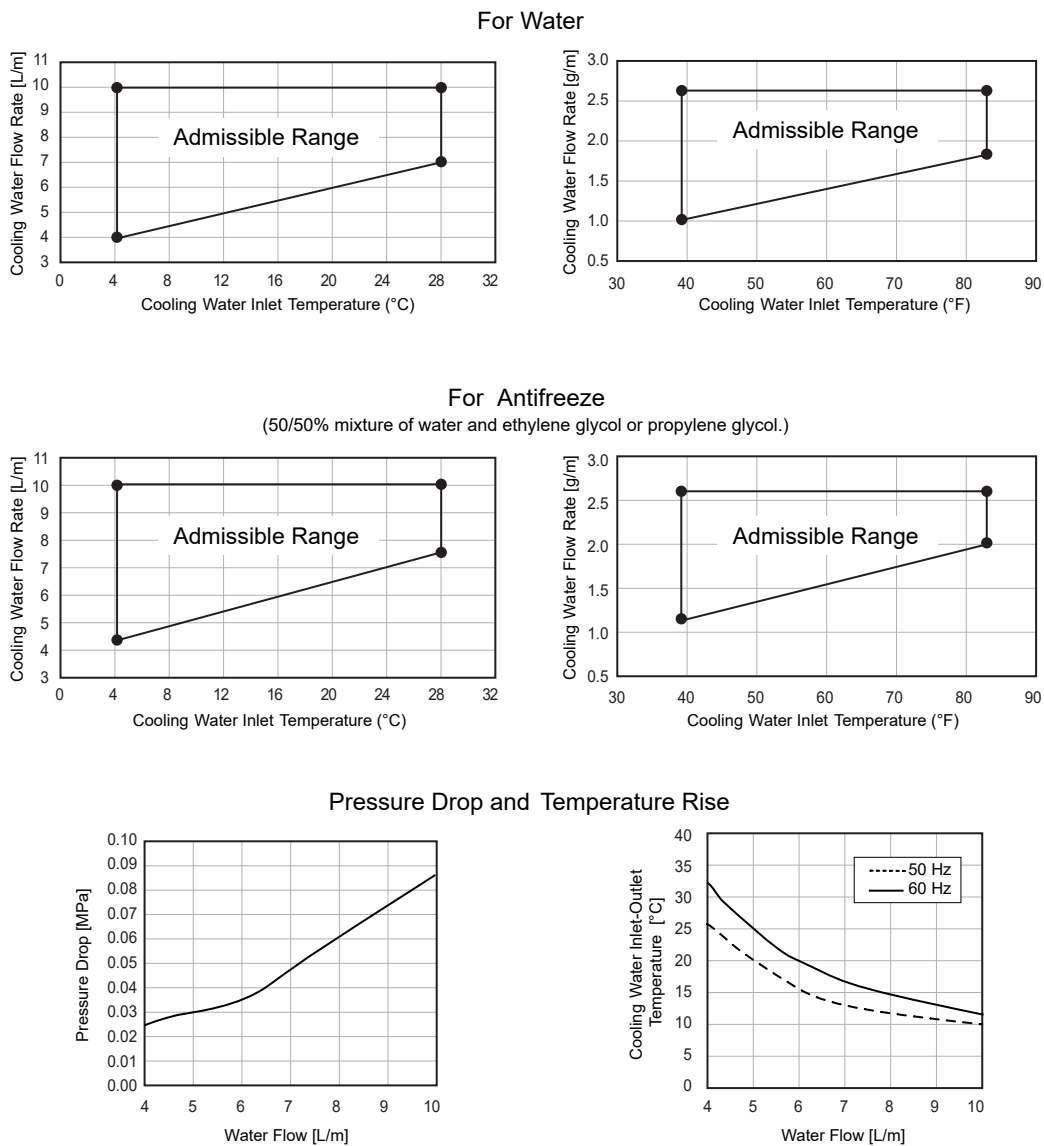

2.10 MR Suite Electrical Requirements



(Applies to all subsections within this section)

2.10.1 General Electrical Requirements

1. Customer is required to install a Main Disconnect Panel (MDP):
 - 1.1. For GE HealthCare supplied MDP, M7100ZA and M7100ZB Design setup, see [2.10.2 GE HealthCare supplied Main Disconnect Panel \(MDP\) Specifications for M7100ZA and M7100ZB on page 50](#).
 - 1.2. For Customer-supplied MDP:
 - 1.2.1. Customer-supplied MDP may not be permissible in all regions. Contact your GE HealthCare Project Manager of Installation (PMI) to verify local requirements.
 - 1.2.2. MDP Design Requirements: [2.10.3 Customer-supplied Main Disconnect Panel \(MDP\) Requirements \(exempt countries only*\) on page 52](#).
 - 1.2.3. MDP Design Setup [Figure 2-27 Customer-supplied MDP MR System Main Disconnect Panel \(MDP\) Setup on page 54](#).
2. At least one remote Emergency Off push-button shall be installed in a location that is visible and accessible to the device operator (Control Room or Magnet Room). The push-button shall be normally closed and require operator action to release after activation (for example, twist and pull). GE HealthCare recommends installing two remote Emergency Off push buttons, installed in the Control Room and Magnet Room.
3. The facility must provide system power to the MDP.
4. All associated transformers and cables must be correctly sized for system power requirements.
5. The facility must connect the GE-supplied cabling from the MDP to the Cryocooler Compressor (F-50SH) in the Integrated Cooling Cabinet.
6. Runs E0009, E3030, M3030 and E4005 are supplied by GE HealthCare. All other wiring shown in [Figure 2-26 GE HealthCare supplied Main Disconnect Panel \(MDP\) Setup for M7100ZA and M7100ZB on page 51](#) and [Figure 2-27 Customer-supplied MDP MR System Main Disconnect Panel \(MDP\) Setup on page 54](#) must be customer-supplied and installed. A customer-supplied substitute for E0009 can be used if the supplied run is shorter than required.
7. All feeder circuits require dedicated ground wires.

Table 2-9 Facility Power Requirements

Component	Parameter	Requirements	
At Main Disconnect Panel (MDP)	Voltage / Frequency	480 VAC	60 ±3 Hz
		415 VAC	50 ±3 Hz, 60 ±3 Hz
		400 VAC	50 ±3 Hz, 60 ±3 Hz
		380 VAC	50 ±3 Hz, 60 ±3 Hz

Table 2-9 Facility Power Requirements

Component	Parameter	Requirements	
	Daily Voltage Variation	Customer to provide +10% / -10% from nominal at MDP input under all line and load conditions. This includes variation of power source and transmission losses up to the MDP.	
	Phase	Input power to the MDP may use one of the following configurations: <ul style="list-style-type: none"> A 3 phase solidly grounded WYE with Ground (3 Wire + Ground) A neutral conductor is not required for MR System operation. If a neutral conductor is present, it can be terminated on the neutral bus provided in the GE HealthCare-supplied MDP. A 3 phase floating DELTA with Ground (3 Wire + Ground). Do not connect a corner grounded DELTA source. Note: Some UPS options may require a neutral (refer to manufacturer documentation for requirements).	
	Phase Balance	Difference between the highest phase line-to-line voltage and the lowest phase line-to-line voltage must not exceed 2%	
	Power Quality	Recommended THD-V of less than 2.5%	
	Facility Zero Voltage Reference Ground	<ul style="list-style-type: none"> The facility ground for the MR System must originate at the system power source (that is, transformer or first access point of power into the facility) and be continuous to the MR System Main Disconnect Panel (MDP) in the room. Main facility ground conductor to Main Disconnect Panel (MDP) must be appropriately sized insulated copper wire. The main facility ground to the Main Disconnect Panel (MDP) must meet local codes. 	
	Power Availability	Continuous facility power is required at all times for operation of the Cryocooler (CRY) to minimize cryogen consumption.	
Service receptacle in Magnet Room	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	Receptacle required for small power tools. Local voltage and portable transformers for voltage values.
	Phase	1	
	Maximum Current	20A (North America) 16A (International)	
Pneumatic Patient Alert	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	The Control Box must be mounted within reach of the operator and within 1.5 m (5 ft.) of an electrical outlet.
	Phase	1	
	Maximum Current	20A (North America) 16A (International)	
Magnet Rundown Unit (MRU)	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	Connection type: Hardwired or permanently wired directly to facility power, no plugs or connectors allowed. 25 mm (1 in.) PVC Schedule 40 Conduit recommended Availability: Continuous Circuit Breaker: Dedicated AC disconnect required for both live and neutral connections
	Phase	1	

Table 2-9 Facility Power Requirements (Table continued)

Component	Parameter	Requirements	
	Maximum Current	1A	
Magnet Monitor (MON)	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	Power at the outlet must be continuously available.
	Phase	1	
	Maximum Current	3A	
Optional MRE Re-soundant Acoustic Driver	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	
	Phase	1	
	Maximum Current	20A (North America) 16A (International)	
Oxygen Monitor (OXY) Option	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	Connection type: Hardwired in unit
	Phase	1	
	Maximum Current	0.9A	
Optional Physiological Acquisition Transceiver (PAT) Charging Station	Voltage / Frequency	100-120 VAC 60 Hz (North America) 200-240 VAC 50/60 Hz (International)	The charging station must be mounted within 3000 mm (118 in.) of an electrical outlet.
	Phase	1	
	Maximum Current	1.5A	

Table 2-10 System Power Demand

Equipment	Power Draw (kVA)
PDU Continuous Power	79
PDU 50 millisecond peak power	120
Cryocooler Compressor Continuous Power	9
Total system Continuous Power	88
Total system 50 millisecond peak power	129

2.10.2 GE HealthCare supplied Main Disconnect Panel (MDP) Specifications for M7100ZA and M7100ZB

The customer is responsible for determining the suitability of the GE HealthCare supplied MDP with respect to governing electrical codes.

The GE HealthCare MDP consists of the following:

- A 3-pole circuit breaker rated for the current of the PDU circuit.
 - 1.1. 150A for M7100ZA
 - 1.2. 160A for M7100ZB
- A 3-pole 200A main circuit breaker rated for the total current of all the sub-breaker circuits.
- A 3-pole 25A circuit breaker rated for the current of the Cryocooler Compressor (CRY).
- All circuit breakers have a short circuit current interrupting rating of 25000 Amps at 480V.
- Auto restart on the Cryocooler circuit following loss and restoration of power.
- Two remote Emergency Off Buttons to be installed external to the MDP. Emergency Off removes power from all outputs when activated. MDP supports maximum E-off cable length of 100 meters when remote EPO push buttons are installed in the field.
- Manual restart push button to restore power after a power outage.
- Terminal blocks that can accept wire sizes for M7100ZA are listed in [Table 2-12 Range of Standard connector accepted for M7100ZA on page 52](#).
- Terminal blocks that can accept wire sizes for M7100ZB are listed in [Table 2-13 Range of Standard connector accepted for M7100ZB on page 52](#).
- Provision for terminating facility incoming neutral wire on the neutral terminal block.
- Multiple ground terminal blocks as required by panel design.
- The global MDP (M7100ZA, countries outside of Europe) is listed and labeled by Nationally Recognized Testing Lab (NRTL) and bears UL and CE mark. The global MDP is certified as per “UL 508A and IEC 61439-2” standards.
- The CE MDP (M7100ZB) bears manufacturer’s CE mark and is certified as per “IEC 61439-2” standards.
- Power on indicators.
- Two isolated, normally open contact pairs that open when e-OFF is pressed or facility power is interrupted for use with optional accessories.
- Capability for single point lock-out/tag-out for the entire system (Mains Disconnect / Input Breaker) and a means to lock-out/tag-out each output breaker independently. All LOTO points support a standard sized hasp for lock-out.

Figure 2-26 GE HealthCare supplied Main Disconnect Panel (MDP) Setup for M7100ZA and M7100ZB

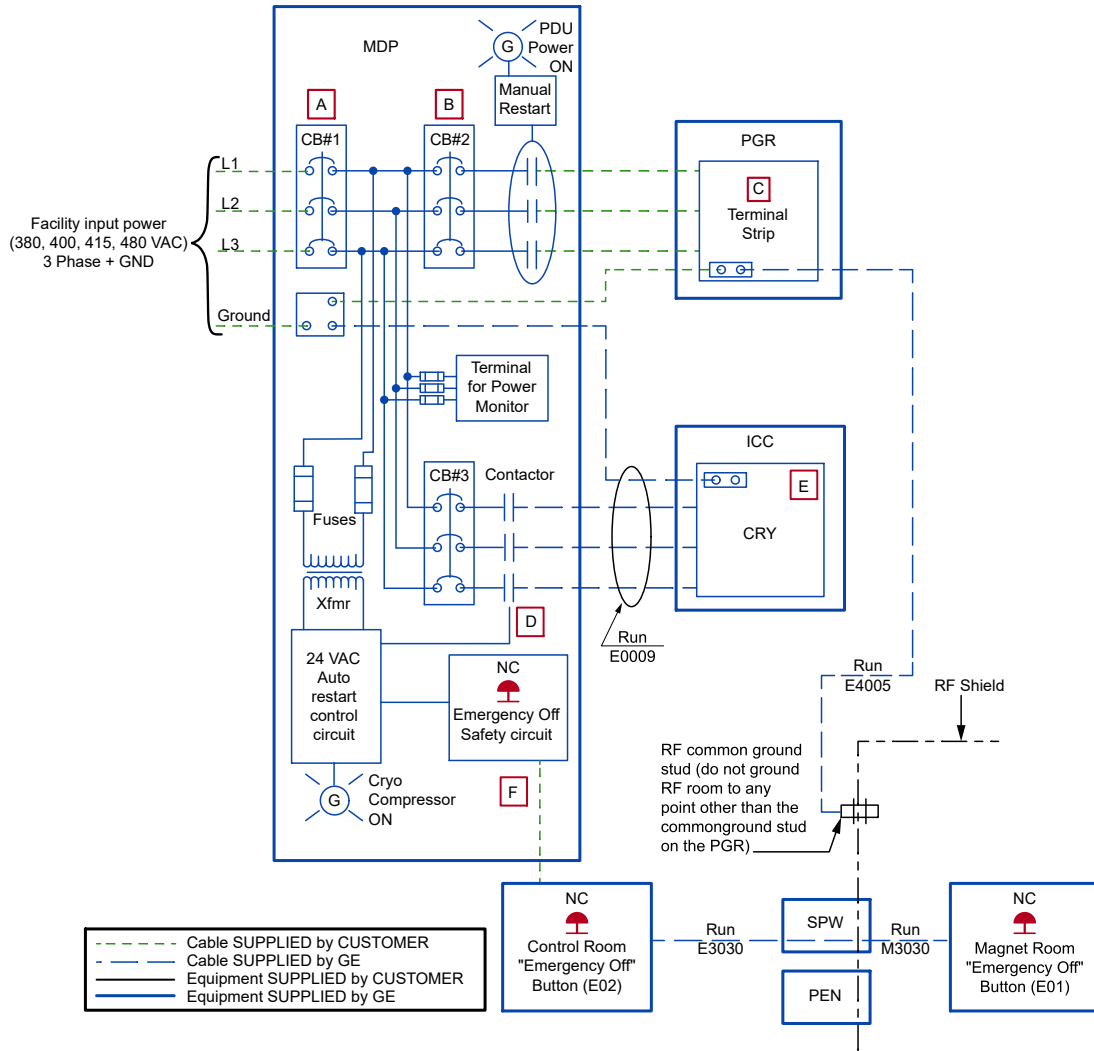


Table 2-11 Specifications for GE HealthCare supplied MDPs

Catalog No.	M7100ZA	M7100ZB
Description	Global MDP 380V - 480V (countries outside of Europe)	CE MDP 380V - 415V; with RCD
Mains Input Breaker, CB1	200A	200A
PDU Breaker, CB2	150A	160A
Compressor Breaker, CB3	25A	25A

The labeled connections shown in [Figure 2-26 GE HealthCare supplied Main Disconnect Panel \(MDP\) Setup for M7100ZA and M7100ZB on page 51](#) can accept the range of standard stranded conductors shown in [Table 2-12 Range of Standard connector accepted for M7100ZA on page 52](#) and [Table 2-13 Range of Standard connector accepted for M7100ZB on page 52](#). All wire types, color, and sizing are to be selected in accordance with governing electrical codes.

Table 2-12 Range of Standard connector accepted for M7100ZA

Item	Phase		Ground	
	sq. mm	AWG/kcmil	sq. mm	AWG/kcmil
A	6 to 120	10 to 250	16 to 120	6 to 250
B	10 to 120	8 to 250	35 to 120	2 to 250
C	2.5 to 70	14 to 2/0	2.5 to 70	14 to 2/0
D	2.5 to 4	14 to 10	2.5 to 16	14 to 10
E	2.5 to 6	14 to 10	2.5 to 6	14 to 10
F	0.5 to 10	22 to 10	-	-

Table 2-13 Range of Standard connector accepted for M7100ZB

Item	Phase		Ground	
	sq. mm	AWG/kcmil	sq. mm	AWG/kcmil
A	6 to 185	10 to 350	16 to 185	6 to 350
B	10 to 120	8 to 250	16 to 185	6 to 350
C	2.5 to 70	14 to 2/0	2.5 to 70	14 to 2/0
D	2.5 to 10	14 to 8	2 to 21	16 to 4
E	2.5 to 6	14 to 10	2.5 to 6	14 to 10
F	0.5 to 10	22 to 10	-	-

2.10.3 Customer-supplied Main Disconnect Panel (MDP) Requirements (exempt countries only*)

**NOTE**

* The requirements listed below apply only to exempt countries. Please contact your GE HealthCare Project Manager of Installation (PMI) for the list of exempt countries.

WARNING**PERSONNEL INJURY OR EQUIPMENT DAMAGE**

Customer supplied MDP must have correctly sized wires and rated components to meet the MR System Power Requirements.

**NOTE**

Refer to [8.5 Sample control schematic for customer-supplied MDP on page 145](#).

1. MDP shall provide Auto-Restart to the Cryocooler Compressor.
2. Manual Restart Capability
 - 2.1. The MDP shall disconnect the PDU circuits upon power loss.
 - 2.2. The MDP shall require a manual restart on the PDU circuits when power is reapplied after an outage.
3. Emergency Off Circuit

- 3.1. The MDP shall have an emergency off control circuit that disables power to the entire MR System.
- 3.2. The emergency off circuit shall be actuated by a push button on the panel, and shall also be capable of being actuated by remotely located push button(s).
- 3.3. Manual reset of the emergency off circuit shall be required to restore power to the entire system.
4. The MDP shall include three breakers, as specified in [Figure 2-27 Customer-supplied MDP MR System Main Disconnect Panel \(MDP\) Setup on page 54](#) and [Table 2-14 Customer Supplied MDP - Required Breaker Size on page 54](#). Output breakers will feed the terminal block on the top of the PGR cabinet and Cryocooler Compressor (located in the Integrated Cooling Cabinet (ICC)).
5. Lock-out/Tag-out:
 - 5.1. The MDP shall provide single point lock-out/tag-out for the entire system and a means to lock-out/tag-out each output breaker independently.
 - 5.2. The lock-out/tag-out feature shall accommodate a standard sized lock hasp.
 - 5.3. The lock-out/tag-out features shall be accessible from the outside of the panel, without the need to open the panel door(s).
6. The MDP shall have a Power ON indicator (Green light) on the panel.
7. The MDP shall meet national/local regulations.
8. The MDP shall provide terminations for all grounds entering, leaving and residing within the panel.
9. The MDP shall provide terminations of appropriate size for all power wiring entering and leaving the panel. All wire types, color, and sizing are to be selected in accordance with governing electrical codes. Refer to [Figure 2-27 Customer-supplied MDP MR System Main Disconnect Panel \(MDP\) Setup on page 54](#)

Figure 2-27 Customer-supplied MDP MR System Main Disconnect Panel (MDP) Setup

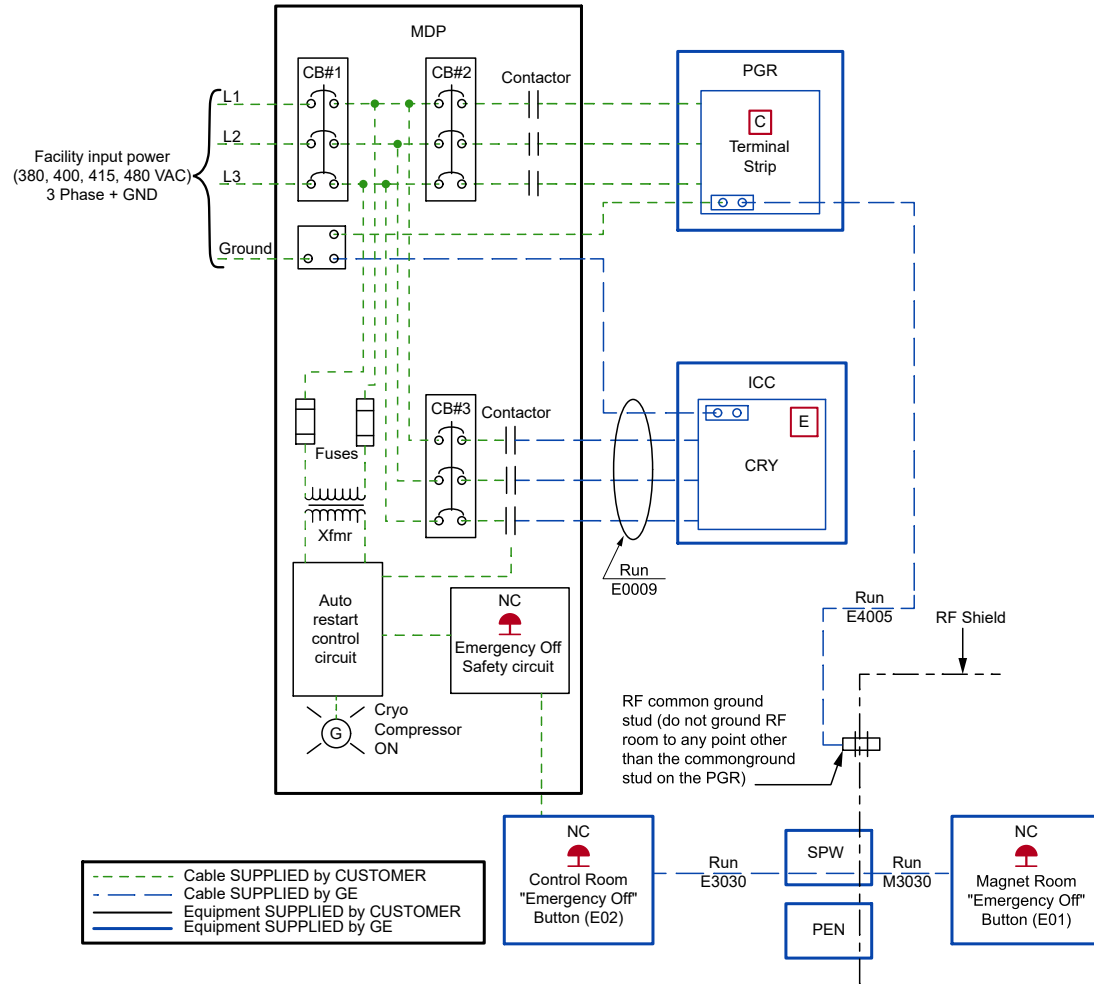


Table 2-14 Customer Supplied MDP - Required Breaker Size

Circuit Breaker	Breaker Size
CB1	200A minimum
CB2	150-160A
CB3	25-30A

Table 2-15 Customer Supplied MDP - Range of Standard Stranded Conductors Accepted

Item	Phase		Ground	
	sq mm	AWG/kcmil	sq mm	AWG/kcmil
A	NA			
B	NA			
C	2.5 to 70	14 to 2/0	2.5 to 70	14 to 2/0
D	NA			
E	2.5-6	14-10	2.5-6	14-10
F	NA			

2.10.4 Emergency Power Backup Specifications (Optional)

The following facility backup power is recommended for continuous operation of the cryocooler compressor and Magnet Monitor:




NOTE

If the compressor must operate on emergency backup power, it still requires chilled water defined in the [2.9.2 Requirements for Emergency Backup Facility Coolant \(Optional\)](#) on page 45.

- Dedicated, single power supply to the compressor
- Magnet Monitor emergency power (110V / 220V, 3A). Refer to [Magnet Monitor \(MON\) Requirements and Specifications](#) on page 103.
- Emergency Off Circuit (E-Off) for the emergency backup to the compressor. LOTO is required for the power source between the generator and compressor.
- A transfer switch to remove the primary power source from the compressor when in emergency backup power mode.

Table 2-16 Specifications for Emergency Power to Cryocooler Compressor

Parameter	Requirements
Power Line Voltage	AC 380, 400, 415V ($\pm 10\%$) / 50 Hz, 3 phase (3 Wire + Ground) AC 480V ($\pm 10\%$) / 60 Hz, 3 phase (3 Wire + Ground) Commercial Power Source <div style="text-align: center;"> <p>WARNING</p>  <p>Do not use an inverter for the main power source.</p> </div>
Operating Current	Max. 13A (Both 50 and 60 Hz)
Starting current	75 / 80 A (50 / 60 Hz)
Maximum Fuse or Circuit Breaker Size	30A
Power Requirement	Minimum 9 kVA Note: The manufacturer recommends a connection capable of 12 kVA.
Power Consumption	Max. 8.3 kW / Steady State 7.5 kW at 60 Hz Max. 7.2 kW / Steady State 6.5 kW at 50 Hz

2.11 MR System Shipping and Receiving



(Applies to all subsections within this section)



IMPORTANT

All shipping dimensions and weights are approximate and may vary based on ship-to location, required rigging, or other requirements. Some shipping or access routes may have requirements in addition to those listed in this section. Contact the GE HealthCare Project Manager of Installation (PMI) to verify magnet shipping, rigging, and access.

2.11.1 Receiving Requirements

1. The customer must provide an area for unloading system components from the truck and delivering to the MR suite



NOTE

Contact GE HealthCare project manager for magnet handling document to be used by rigging companies.

2. The customer is responsible for ensuring:
 - 2.1. All floors along the route will support the weight of the magnet (GE HealthCare recommends a structural analysis)
 - 2.2. Doors or other openings are sufficiently wide to allow passage
 - 2.3. Sufficient room is provided for any required rigging tools

2.11.2 Facility Delivery Route Requirements

The following table lists the delivery dimensions of system components. Upon delivery, verify the component dimensions and weight. The delivery route must be planned to accommodate the dimensions listed.

Table 2-17 Delivery Route Requirements

Component	Width		Depth		Height		Weight		Notes
	mm	in.	mm	in.	mm	in.	kg	lb.	
Magnet	See Magnet shipping dimensions in Table 2-18 MR System Component Shipping Specifications on page 57 .								
Cryogen	Dimensions vary depending on dewar type used. Verify with cryogen supplier.								
PGR Cabinet	2000	78.75	890	35.03	2108	83	1510	3329	Cabinets are moved with dollies attached to each side (adding 520 mm (20 in.) to the width and 82 kg (180 lb.) to the weight). Cabinets must be raised to
ICC	1100 With dolly: 1710	43.3 With dolly: 67.3	842	33.1	1970	77.6	560 With dolly: 642	1235 With dolly: 1415	

Table 2-17 Delivery Route Requirements (Table continued)

Component	Width		Depth		Height		Weight		Notes
	mm	in.	mm	in.	mm	in.	kg	lb.	
pICC	830 With dolly: 1440	32.7 With dolly: 56.7	811	31.9	1983	78.1	424 With dolly: 506	935 With dolly: 1115	remove the pallet but may be lowered almost to floor level while moving. The ICC or pICC can be moved on its casters.
PEN Cabinet	1120	44.1	1001	39.4	1930	76	313	691	

2.11.3 MR System Component Shipping Specifications

MR System component shipping dimensions and weight are listed below:

Table 2-18 MR System Component Shipping Specifications

Component	W x D x H		Weight		Notes
	mm	in.	kg	lb.	
UA Magnet (as shipped with lifting bars)	2393 x 2388 x 2317	94.2 x 94 x 91.2	7621	16800	Domestic - Tarp
UA Magnet (Crated for International Shipping)	2438 x 2896 x 2438	96 x 114 x 96	8528	18800	International - crate/pallet
AR Magnet (as shipped with lifting bars)	2311 x 2388 x 2250	91 x 94 x 88.6	5882	12968	Domestic - Tarp
AR Magnet (Crated for International Shipping)	2438 x 2896 x 2438	96 x 114 x 96	6633	14622	International - crate/pallet
Magnet Accessory Equipment	1219 x 1219 x 711	48 x 48 x 28	182	400	Crate
VibroAcoustic Damping Kit for UA magnet	914 x 1651 x 305	36 x 65 x 12	261	575	Box on pallet
Cryocooler Compressor	660 x 711 x 1067	26 x 28 x 42	125	275	Pallet with box cover
Rear Pedestal Assembly with Rear Split Bridge Assembly, Low Profile Carriage Cover	864 x 1473 x 1219	34 x 58 x 48	141	310	Box on pallet
GEM Enclosure, Domestic	2362 x 991 x 874	93 x 39 x 34	68	150	Box on pallet
GEM Enclosure, International	2388 x 1675 x 711	94 x 66 x 28	68	150	Box on pallet
Patient Table	2388 x 737 x 965	94 x 29 x 38	301	663	Pallet
Power Gradient RF Cabinet (PGR) (Domestic)	1701 x 1041 x 2235	67.25 x 41 x 88	1567	3455	Pallet

Table 2-18 MR System Component Shipping Specifications (Table continued)

Component	W x D x H		Weight		Notes
	mm	in.	kg	lb.	
Power Gradient RF Cabinet (PGR) (International)	1727 x 1067 x 2337	68.25 x 42 x 92	1680	3704	Pallet
Penetration Panel Cabinet (PEN) (Domestic)	1117 x 895 x 2051	43.75 x 35.25 x 80.75	345	761	Pallet
Penetration Panel Cabinet (PEN) (International)	1137 x 895 x 2153	44.75 x 35.25 x 84.75	373	823	Pallet
Integrated Cooling Cabinet (ICC)	1338 x 952 x 2180	53 x 38 x 86	716	1575	Crate
Platform Integrated Cooling Cabinet (pICC)	1029 x 959 x 2191	41 x 38 x 87	486	1071	Pallet plus crate
Secondary Penetration Wall	533 x 610 x 1626	21 x 24 x 64	46	101	Pallet with cardboard cover
Operator Workspace Cabinet	610 x 889 x 787	24 x 35 x 31	110	243	Wood pallet with cardboard cover
Operator Workspace Display	686 x 838 x 686	27 x 33 x 27	9	19	Box
Operator Workspace equipment	813 x 813 x 584	32 x 32 x 23	45	100	Box
Operator Workspace Table	1041 x 1372 x 381	41 x 54 x 15	75	165	Box

Table 2-19 MR System Service Component Shipping Specifications

Component	W x D x H		Weight		Notes
	mm	in.	kg	lb.	
Shim Tray Extraction Tool (Box 1)	1016 x 2591 x 648	40 x 102 x 25.5	216	477	Case on pallet (pallet weight 102 kg (225 lb.))
Shim Tray Extraction Tool (Box 2)	1219 x 1473 x 1372	48 x 58 x 54	268	591	Case on pallet (pallet weight 91 kg (200 lb.))

Table 2-20 MR System Optional Component Shipping Specifications

Component	W x D x H		Weight		Notes
	mm	in.	kg	lb.	
8 kW MNS Amplifier Cabinet	600 x 1016 x 1637	23.6 x 40 x 64.4	276.7	610	On cabinet casters, wrapped with plastic

Figure 2-28 Magnet Dimensions (as Shipped) for UA Magnet

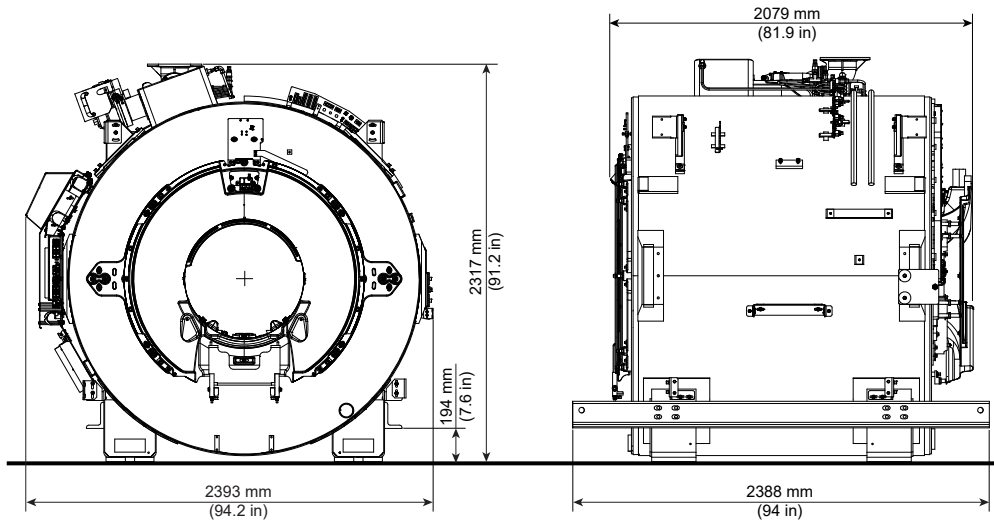


Figure 2-29 Magnet Dimensions (as Shipped) for AR Magnet

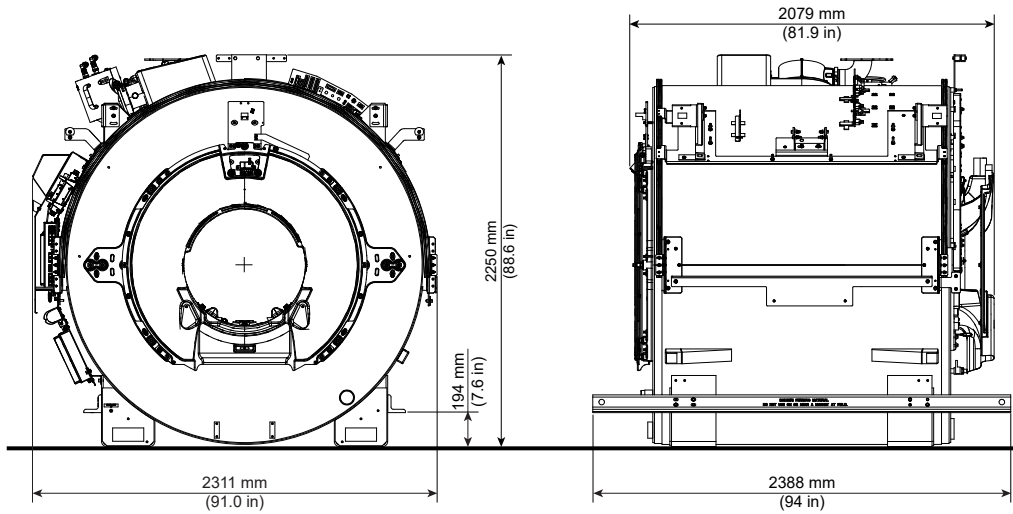


Table 2-21 MR System Component Replacement Shipping Specifications


Component	Component Location	W x D x H		Weight	
		mm	in.	kg	lb.
Split Bridge	Magnet Room	546 x 1963 x 178	21.5 x 77.3 x 7	18	40
Replacement RF Body Coil	Magnet Room	762 x 762 x 1524	30 x 30 x 60	116	255
Replacement XRM Gradient Coil Assembly on a Shipping Cradle/Cart	Magnet Room	910 x 2444 x 1499	35.8 x 96.2 x 59	1449	3194
3.0T RF Amplifier	Equipment room (front of PGR)	671 x 864 x 437	26.4 x 34 x 17.2	181	400
Gradient Power Supply	Equipment room (front of PGR)	508 x 1010 x 565	20 x 39.8 x 22.2	181	400
Gradient Cooling Unit (GCU)	Equipment room (front of ICC)	920 x 605 x 1445	36.2 x 23.8 x 56.9	144.5	318.6
Cabinet Cooling Unit (CCU)	Equipment room (front of ICC)	920 x 705 x 1550	36.2 x 27.8 x 61.0	195.5	431
CCU Pump	Equipment room (front of pICC)	610 x 902 x 607	24 x 35.5 x 23.9	144.5	318.6
Gradient Coil Replacement Tool Kit Crate	At site near magnet room	787 x 2184 x 952	31 x 86 x 37.5	340	750



NOTE

The dimensions and weights listed for the components in [Table 2-21 MR System Component Replacement Shipping Specifications on page 60](#) include packaging.

Figure 2-30 Gradient Coil Cart

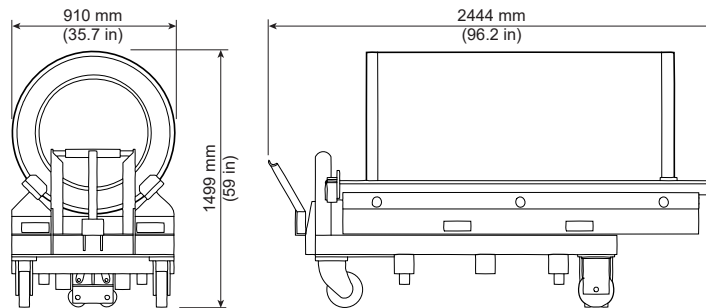

2.11.4 Temperature and Humidity Storage Requirements

MR systems and components must be stored within the environmental requirements listed below.



NOTE

Some equipment is liquid-cooled. After coolant is added, the equipment must be kept from freezing. Phantoms and the coolant itself must also be kept from freezing.

Room	Temperature		Humidity	
	Range °C (°F)	Change °C/Hr (°F/Hr)	Relative % (Non-condensing)	Change %/Hr
Electronic Cabinets and Equipment	-30 — 60 (-22 — 140)	20 (68)	10 — 80	30
Magnet Assembly	-30 — 60 (-22 — 140)	20 (68)	10 — 80	30
Surface Coil	-30 — 50 (-22 — 122)	20 (68)	10 — 95	30

3 Magnet Room


3.1 Magnet Room Introduction



The Magnet Room is best understood as a series of layers, or “rooms within a room.” Each of these rooms has a specific function and associated requirements. All requirements in this chapter must be followed to ensure safe and correct operation of the MR System.

1. The Magnetic shielded room contains the MR Magnet fringe field within a confined space. A site survey is required to determine magnet shield requirements (not all sites require magnetic shielding). Because of the added cost of magnetic shielding, room location should be carefully considered.
2. The Acoustic room is a layer used to help attenuate the noise produced during a scan. An acoustic engineer is strongly recommended to assess the environment.
3. The RF shielded room is critical to the correct MR System operation. RF shielding prevents interaction of external RF radiation with MR System operation and it also prevents MR System RF radiation from interfering with external systems, such as aircraft control. Special care must be used when installing all fixtures penetrating the RF shield (for example, vents, electrical conduit, penetration panels) to ensure the integrity of the RF shield is maintained. Refer to *RF Shielded Room Requirements*, 5850260-1EN.
4. The Finished room includes the wall coverings, ceiling tile, ceiling grid, other fixtures, Magnet (MAG) and Patient Table (PT). When planning the finished room, ensure the following:
 - 4.1. All building codes are met (such as maintaining egress routes).
 - 4.2. Items which may generate or create RF interference (including fluorescent lighting) are not allowed for installation within the Magnet Room.
 - 4.3. Customer is responsible for the selection and installation of all locally required safety devices (for example, smoke detectors, oxygen monitors, and so on).
 - 4.4. Smoke detectors should be located outside of the Magnet Room (for example, within the return air duct) whenever possible. If code does not allow this, use only simple two wire non-addressable smoke detectors in the Magnet Room.
 - 4.5. Ferrous or metallic items which could become projectiles when the magnet is installed (including wall coverings, ceiling tile, ceiling grid, or other fixtures) are not used or are correctly secured.

Figure 3-1 Magnet Room Layers

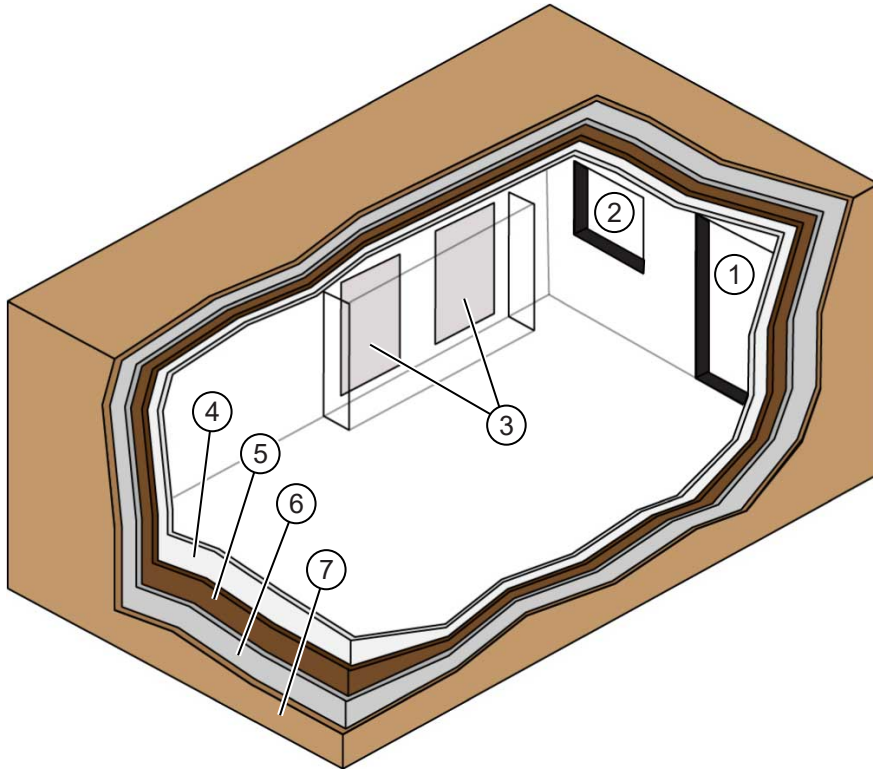

Item	Description	Item	Description
1	Door	2	Window
3	Penetration wall(s)	4	Finished room
5	RF shielding	6	Acoustic barrier
7	Magnetic shielding	-	-



NOTE

The sequence of the room layers can vary based on siting needs.

3.2 Magnet Room Structural Requirements



This section lists the structural requirements that must be considered when performing site evaluation and planning of the Magnet Room.

3.2.1 Overview



1. When preparing a building plan or evaluating a potential site for an MR System, take care to ensure the MR suite will not interact with the surrounding environment (that is, magnetic, acoustic, environmental steel, and vibration).
2. The customer is responsible for vibration testing required to verify suitability of a proposed site. All test results and any questions regarding testing, results, or analysis must be forwarded to the GE HealthCare Project Manager of Installation (PMI).

3.2.2 Environmental Steel Limits



A static magnetic field extends in a three-dimensional space around the magnet isocenter. Environmental steel within the static magnetic field affects the uniformity (or homogeneity) of the field. Field uniformity is critical to both image quality and chemical shift analysis (spectroscopy). An analysis of the environmental steel is required within a 2.5 m (8.2 ft.) spherical radius of the magnet isocenter. Environmental steel includes ferrous pipes, beams, concrete rebar, or any other structural steel in the floors, walls, or ceiling.

The following floor items must be limited per [Table 3-1 Steel Mass Limits to Magnet Isocenter \(3 x 3 m \(10 x 10 ft.\) Area Under Magnet\) on page 65](#).

1. Non-movable steel construction material such as rebar and metal decking
2. Existing or proposed RF/magnetic shielding or shim plates
3. [Table 3-1 Steel Mass Limits to Magnet Isocenter \(3 x 3 m \(10 x 10 ft.\) Area Under Magnet\) on page 65](#) defines the limits of use as a guideline to help the customer understand allowable amounts of ferrous rebar, steel decking, or other components as they design the MR suite and Magnet Room floor structure.
4. The customer must provide detail defining ferrous material below the magnet to the Project Manager so the GE HealthCare MR Siting and Shielding (MRSS) team can review for compliance.

Table 3-1 Steel Mass Limits to Magnet Isocenter (3 x 3 m (10 x 10 ft.) Area Under Magnet)

Limits Of Steel Mass kg/m ² (lb./ft ²)	Distance Below Top Surface Of Floor mm (in.)
0 (0)	0-76 (0-3)
9.8 (2)	76-127 (3-5)
14.7 (3)	127-254 (5-10)
39.2 (8)	254-330 (10-13)
98.0 (20)	330+ (13+)


3.2.3 Vibration Requirements



Excessive vibration can affect MR image quality. Vibration testing must be performed early in the site planning process to ensure vibration is minimized. Both steady state vibration (exhaust fans, air conditioners, pumps, and so on) and transient vibrations (traffic, pedestrians, door slamming, and so on) must be assessed (see [Figure 3-2 Magnet Steady State Vibration Specifications for UA magnet on page 66](#) and [Figure 3-3 Magnet Steady State Vibration Specifications for AR magnet on page 67](#)). Specific requirements for vibration mitigation, include:

1. The Magnet (MAG) cannot be directly isolated from vibration. Any vibration issue must be resolved at the source.
2. MR Suite HVAC must have vibration isolation.
3. A vibration analysis must be performed at the proposed site with the results (and any mitigation) forwarded to the GE HealthCare Project Manager of Installation (PMI). See the [8.2 MR Site Vibration Test Guidelines on page 137](#).
4. A transient vibration test must only be performed after a steady-state test has been performed and all steady-state sources of vibration have been mitigated.
5. Transient vibration levels above the specified limits in the [8.2 MR Site Vibration Test Guidelines on page 137](#) must be given to the PMI for review.
6. Any transient vibration that causes vibration to exceed the steady-state level must be mitigated.
7. The vibration test consultant must account for non-mechanically induced signals such as test equipment instabilities, thermal drift or RF interference.


Figure 3-2 Magnet Steady State Vibration Specifications for UA magnet

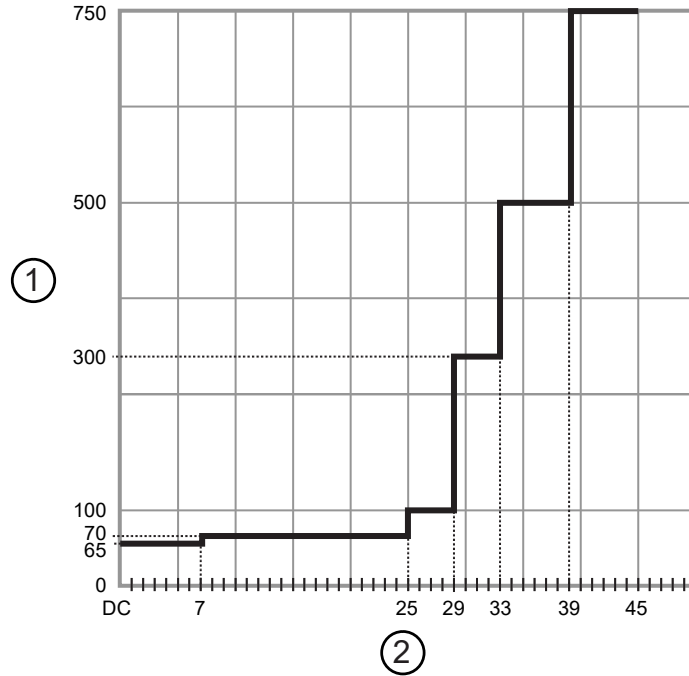

Item	Description
1	Acceleration ($g \times 10^{-6}$) (RMS) above ambient baseline
2	Excitation frequency (Hz)

Acceleration ($g \times 10^{-6}$) (RMS) above ambient baseline	Excitation frequency (Hz)
65	0 - 7
70	7 - 25
100	25 - 29
300	29 - 33
500	33 - 39
750	39 - 45

Figure 3-3 Magnet Steady State Vibration Specifications for AR magnet

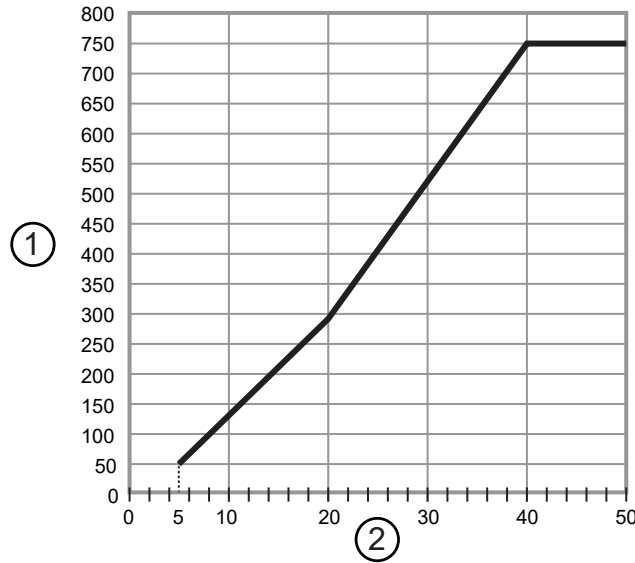

Item	Description
1	Acceleration (g x 10 ⁻⁶) (RMS) above ambient baseline
2	Excitation frequency (Hz)

3.3 Magnetic Shielded Room Requirements



Magnetic shielding prevents interaction between the magnet and nearby sensitive devices. Because of the added cost of magnetic shielding, room location should be carefully considered. All sites, including upgrade sites, must be evaluated for magnetic shielding requirements. Existing magnetic shielding at an upgrade site may not be sufficient for the new system. Contact the GE HealthCare Project Manager of Installation (PMI) to request a site evaluation.

See [MR Suite Magnetic Field Specifications on page 24](#) for detailed magnetic proximity limit information.

1. The GE HealthCare Project Manager of Installation (PMI) works with the customer to coordinate the magnetic shielding site evaluation.
2. The customer is responsible for installation of all magnetic shielding.
3. If rear wall magnetic shield or steel RF wall is closer than 2500mm (98.4 in.) from isocenter, it should be verified by GE HealthCare PMI.

3.4 Penetration Panel Wall Opening Requirements



1. The Equipment Room and the Magnet Room must share at least one common wall to mount the PEN panel and SPW.
2. The penetration panel opening requirements are shown in [Figure 3-4 Penetration Wall Openings on page 69](#).
3. Two GE HealthCare-supplied penetration panel adaptor frames are used to connect the PEN panel and SPW to the wall openings. Connection details are shown in [Figure 3-4 Penetration Wall Openings on page 69](#).
4. The RF vendor must supply 40 fasteners for hole size of $\text{Ø}7.25$ mm (0.285 in) (total quantity for both panels) and install each penetration panel adaptor frame to the wall through the outer ring of holes in the plate.

**NOTE**


The length of the fasteners should be based on the thickness of the RF Wall.


**NOTE**

GE HealthCare-supplied bolts connect the penetration panels to each adaptor frame through the inner ring of holes in the plate.

5. RF shielding in the Magnet Room walls must completely wrap around all edges of the penetration wall openings and continue through the openings to the inside of the Equipment Room to provide full shielding continuity with the adaptor frame. The minimum overlay in the Equipment Room is 60 mm (2.4 in.).
6. A minimum of 102 mm (4 in.) is required between the floor and the bottom of the penetration panel opening to provide clearance for the adaptor frame.
7. The maximum distance between the penetration panels is determined by the available length of the cables between the RF common ground stud and pen panels. Refer to [7.1.2 Available Cable Lengths on page 117](#).
8. Refer to [3.5.4 Penetration Panel Closet on page 74](#) for service hatch requirements.

Figure 3-4 Penetration Wall Openings

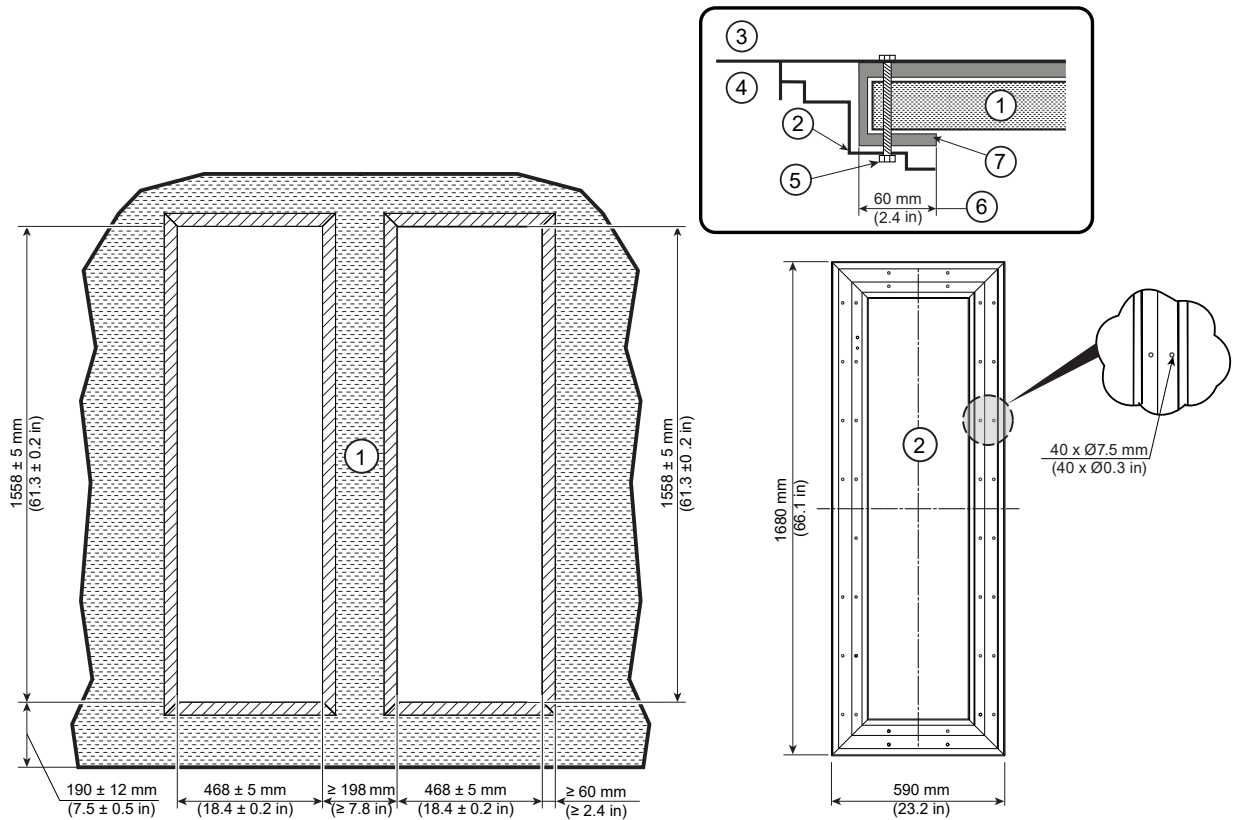

NOTE
 The wall detail cross section is not to scale.

Item	Description	Item	Description
1	RF shield supporting structure (wall or double-sided galvanized panels, and so on)	5	Bolts are provided by screen room vendor and bolt holes are drilled at installation. Bolt length varies by wall thickness.
2	PEN Frame	6	Minimum exposed RF shield material (all around)
3	Magnet Room	7	RF Shield Material
4	Equipment Room	-	-

3.5 Finished Room Requirements

3.5.1 Ferrous Materials in the Magnet Room



1. Non-ferrous (non-metallic) materials or components should be used in the Magnet Room.
2. Ferrous components or material in the Magnet Room that could be removed for servicing, cleaning, or replacement must be secured to prevent the ferrous material from becoming a projectile (ferrous components or material must also be identified as ferrous to prevent untrained personnel from working on the ferrous material while the magnet is energized).

3.5.2 Walls



Refer to *Acoustic Room Details*, 5850262-1EN. Hard, bare wall surfaces may create a harsh Magnet Room acoustic environment due to reflection of sound waves. Finished walls with acoustic properties can reduce reflected noise.

1. GE HealthCare recommends finished walls to protect the RF shielding.
2. Walls and any millwork, cabinets, storage areas, acoustic coverings, and so on, must remain outside the minimum service area.
3. A metallic electrical conduit inside walls and ceilings may be used. Conduit for receptacles must be metallic.

3.5.3 Magnet Preinstallation Markings



For correct cryogen venting, the magnet vent adaptor must align correctly with the ceiling vent when the magnet is installed.

1. The magnet isocenter position must be clearly marked, and the marking must be identifiable throughout construction.



NOTE

If there is no ceiling grid in the room, we recommend to also mark the magnet isocenter location on the ceiling. This can serve as a reference for positioning a new vent pipe, or can be used to reproduce the floor markings if they become lost during construction.

Refer to *Magnet Room Venting Requirements*, 5850263-1EN for the location of the magnet vent.

2. If no ceiling vent pipe exists prior to construction, the location of the magnet isocenter and magnet z-axis orientation must be marked on the Magnet Room floor as shown below.

Figure 3-5 Marking Magnet Isocenter (in a Room Without a Vent Pipe)


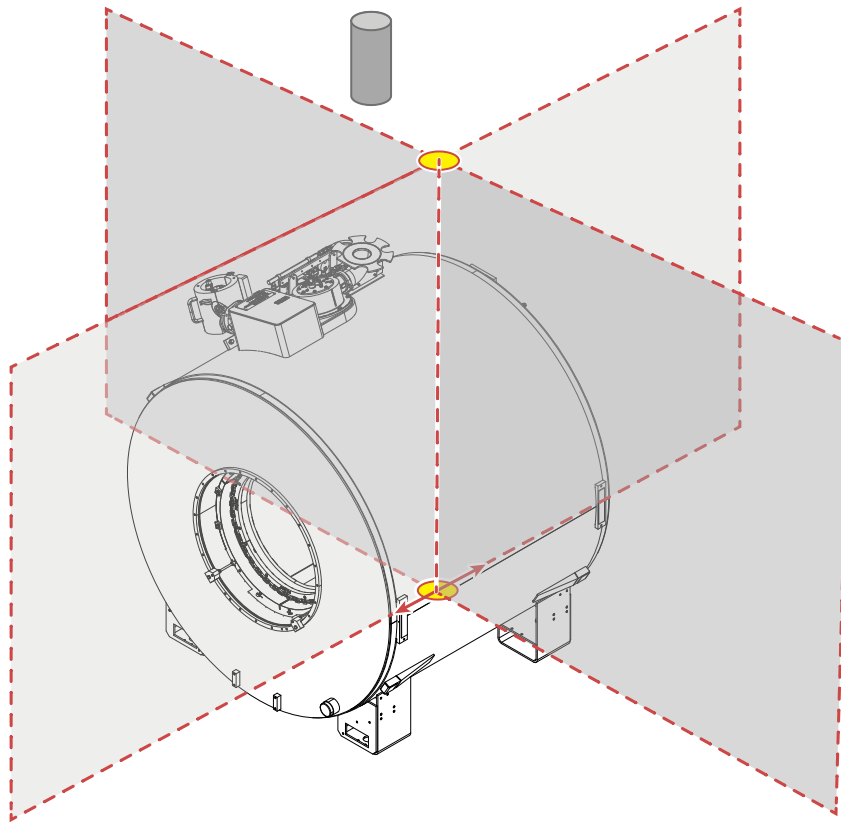


Figure 3-6 Marking the Magnet Isocenter and the Z-Axis Orientation on the Floor

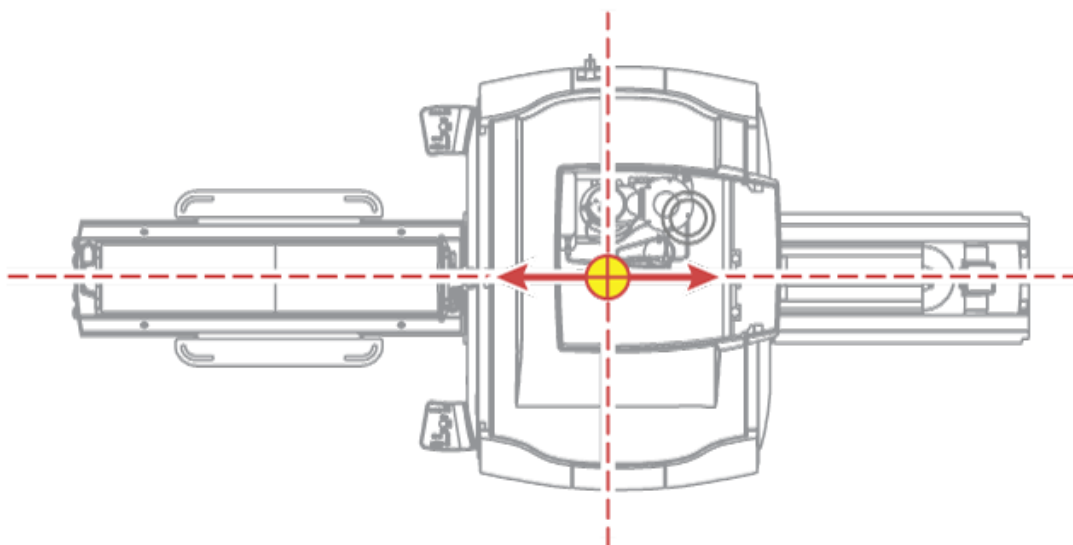


- 3. If a ceiling vent pipe already exists at the start of construction, the magnet isocenter marking must be correctly aligned relative to the position of the vent pipe. Refer to *Magnet Room Venting Requirements*, 5850263-1EN for the location of the magnet vent.

 **NOTE**

We recommend:


- Marking the center of the vent pipe.
- Using a cross line laser or other accurate method to transcribe it to the floor directly below.
- Marking the magnet isocenter relative to the vent pipe marking on the floor. The positioning template can be used as a marking aid.

Figure 3-7 Marking the Center of the Vent Pipe

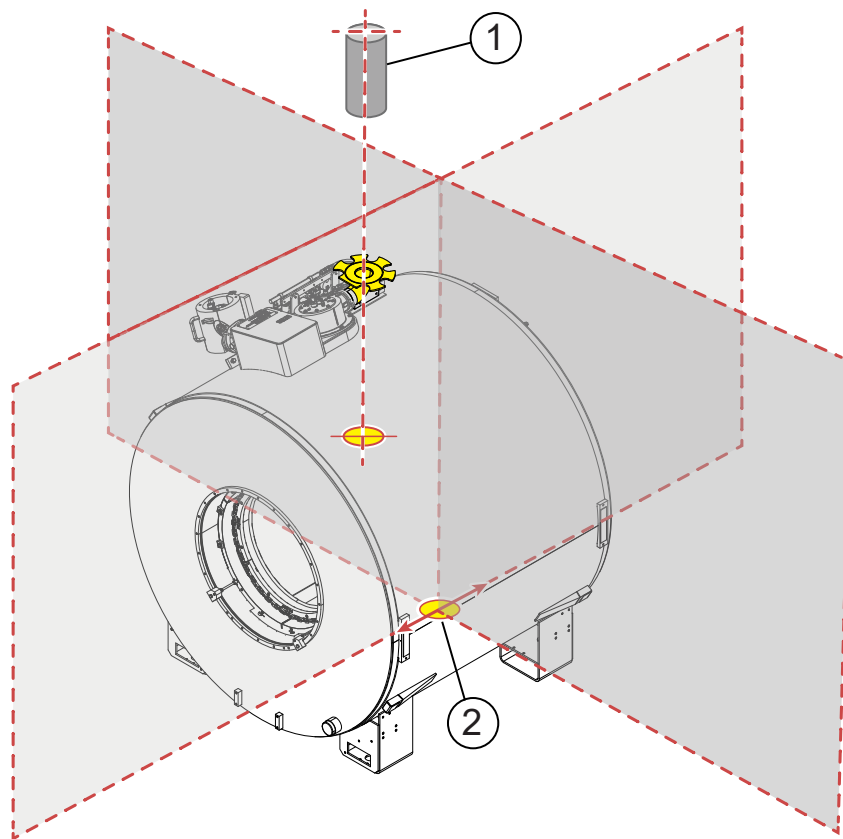

Item	Description
1	Vent center
2	Magnet isocenter

Figure 3-8 Marking the Vent Pipe and Magnet Isocenter on the Floor

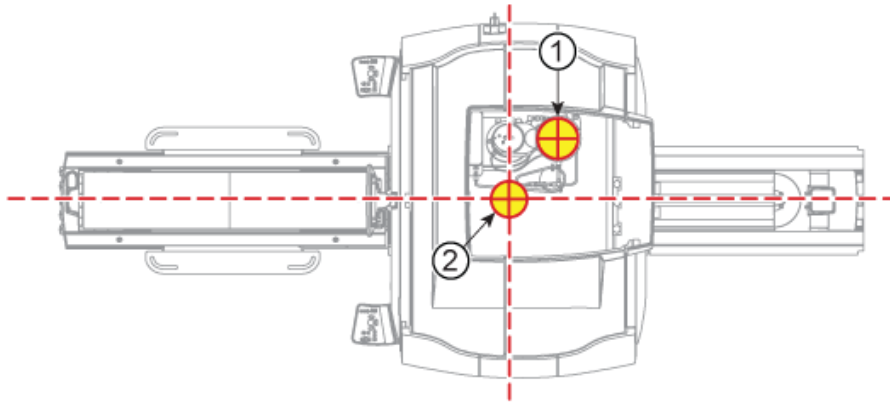

Item	Description
1	Vent center
2	Magnet isocenter

3.5.4 Penetration Panel Closet



1. An enclosure (that is, the penetration panel closet) must be provided to restrict access to the penetration panel(s) and for storage of excess interconnections.
 - 1.1. The penetration panel closet must meet the minimum penetration panel closet outline as shown in [3.5.5 Penetration Panel Closet Specifications on page 75](#).
 - 1.2. The penetration panel closet must have a mechanical locking mechanism to restrict access to the penetration panels.
 - 1.3. The penetration panel closet may be expanded to provide an area for excess cable storage with the following requirements:
 - 1.3.1. Excess cable must not be stored within the minimum closet service area.
 - 1.3.2. Excess cable must not interfere with access or servicing of the Penetration Panel (PP) or SPW.
 - 1.3.3. The area within the penetration panel closet to store the cable should be sized to accept a 711 mm (28 in.) cable loop (2x the minimum bend radius of the largest cable).
 - 1.3.4. Refer to [7.1.6 Storage Requirements for Excess Gradient Cable on page 121](#) for additional excess cable storage requirements.
 - 1.4. The penetration panel closet must allow free air exchange of 680 m³/hour (400 cfm) between the Magnet Room and penetration panel closet for MR System blowers. Airflow may be achieved through door louvers or other openings in the penetration panel closet that meet all other penetration panel closet requirements.
 - 1.5. The primary source of airflow must be from the Magnet Room. Openings into the area above a false ceiling or other storage areas should be minimized.
 - 1.6. The penetration walls may be enclosed by separate closets with the following requirements:
 - 1.6.1. The maximum distance between the penetration panels is determined by the available length of the cables between the RF common ground stud and pen panels. Refer to [7.1.2 Available Cable Lengths on page 117](#).
 - 1.6.2. The separate closets must meet all other service area requirements for each penetration panel.
 - 1.6.3. Airflow as listed above must be provided for both closets.
 - 1.6.4. Both closets must have mechanical locks.


3.5.5 Penetration Panel Closet Specifications

- Minimum Closet Depth: 470 mm (18.5 in)
- Minimum Closet Width: 1669 mm (65.7 in)
- Minimum Closet Height: 2007 mm (79 in)

NOTE
 If the closet depth is less than 914 mm (36 in), the closet must have door(s) configured to clear the service area for the Penetration Panel.

NOTE
 Minimum closet specifications do not account for storage space for excess gradient cable. Refer to [7.1.6 Storage Requirements for Excess Gradient Cable](#) on page 121 for additional excess cable storage requirements.

Figure 3-9 Minimum dimensions of closet in the Magnet Room

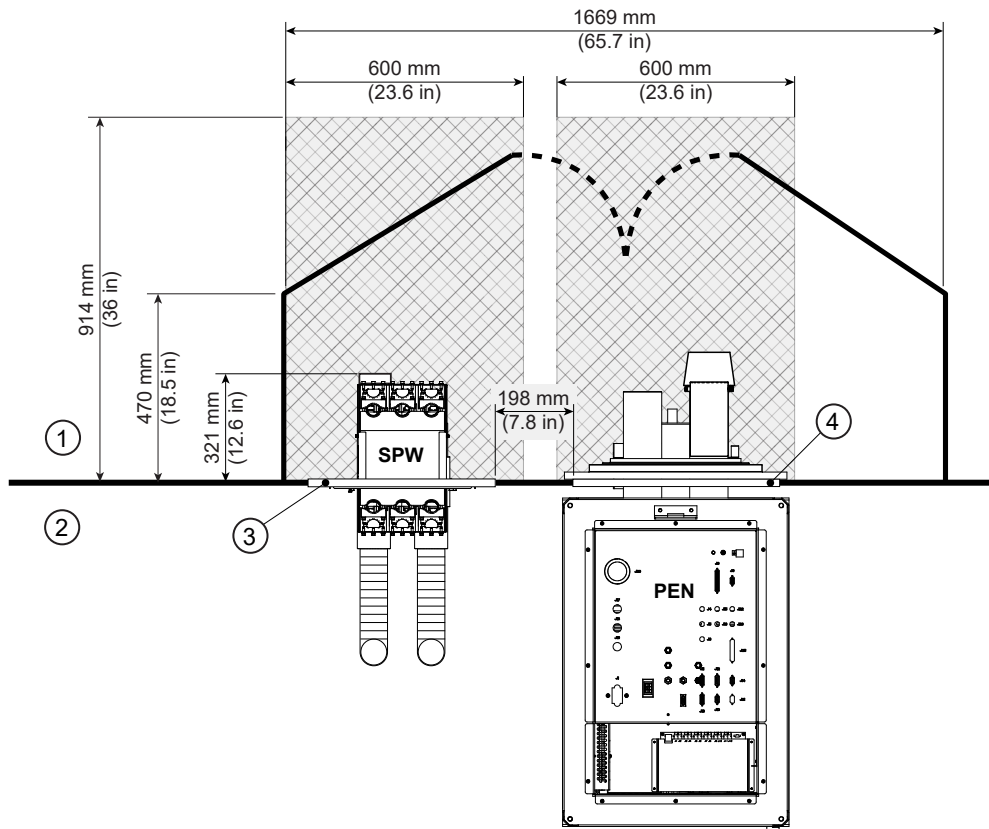

Item	Description		
1	Magnet Room	3	SPW Wall Opening
2	Equipment Room	4	PEN Wall Opening

3.5.6 Doors, Magnet Access Openings, and Patient Viewing Windows



1. The finished opening of the Magnet Room main door must be at least 1092 mm (43 in.) wide to allow for helium dewars and patient tables.
2. Threshold height must not exceed 15 mm (0.6 in.) on both sides of the door with a maximum 10-degree threshold inclination.
3. IEC requires the patient, while in the bore, be in full view of the operator.



NOTE


- GE HealthCare recommends using a window, although other means (for example, camera and video display) may be used as long as all IEC requirements are met.
 - The recommended dimensions for the patient viewing windows are 1219 mm wide x 762 mm high (48 in. wide x 30 in. high).
 - The recommended distance from the bottom edge of the patient viewing window to the finished floor is 1067 mm (42 in.).
4. The magnet delivery requires an opening into the room to allow access for the magnet delivery, rigging, and personnel access.

3.5.7 Finished Ceiling



1. The customer is responsible for the finished ceiling.
2. The finished ceiling grid must be non-ferrous.
3. An optional cable concealment kit is available from GE HealthCare.
4. Ceiling preparation should be completed prior to magnet delivery. See [Figure 3-10 Cable Concealment Kit Ceiling Opening Dimensions \(Top View\)](#) on page 77 for ceiling opening dimensions.

Figure 3-10 Cable Concealment Kit Ceiling Opening Dimensions (Top View)

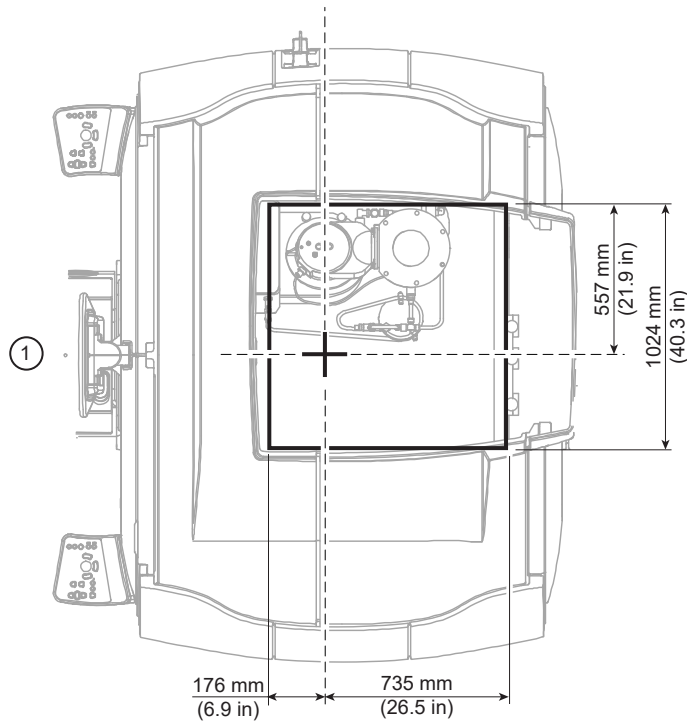

Item	Description
1	Front of magnet

3.5.8 Magnet Room Floors



1. The finished floor must support the weight of all components throughout operation and service life. This includes the magnet, patient table, and gradient coil replacement cart.




NOTE

For gradient coil replacement, field engineers remove the patient table from the Magnet Room before they move the gradient coil replacement cart into the Magnet Room.

2. The finished floor must be water resistant to protect the subfloor and shielding from water damage.
3. The customer is responsible for providing flooring to prevent ESD (Electrostatic Discharge) buildup to 8 kV for protection of the sensitive MR equipment.
4. Magnet, Enclosure, and Patient Table areas must be flat and level within 3 mm (0.125 in.), with the magnet in place, within the shaded area shown in [Figure 3-11 Magnet Room Floor Levelness Area on page 78](#).

Figure 3-11 Magnet Room Floor Levelness Area

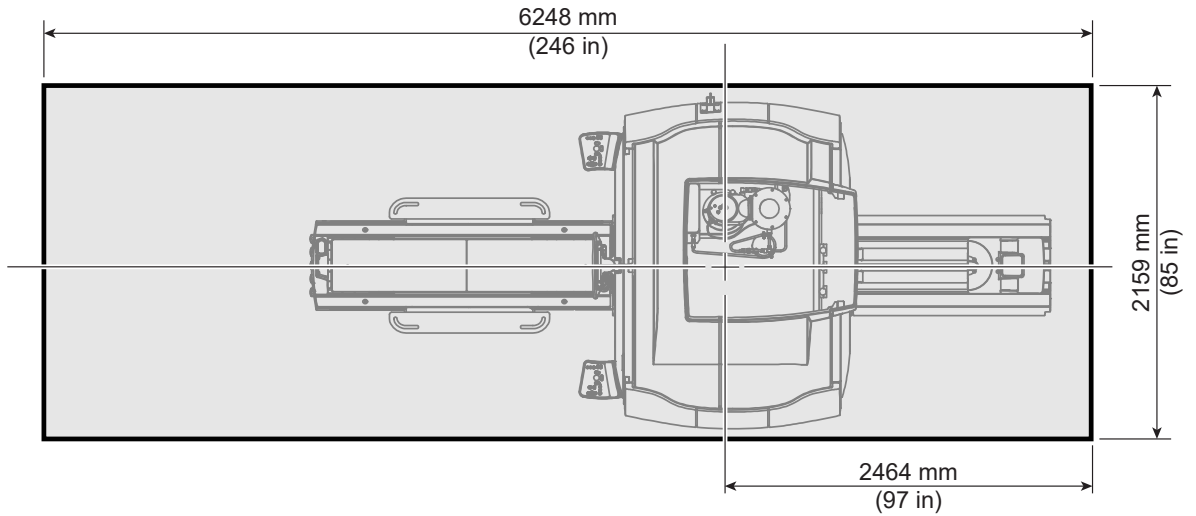


5. **For UA series magnet only:** The VibroAcoustic Dampening kit for systems that do not require Seismic mounting is M7000VA. See below for dimensions that aid in the calculation for floor loading. See [Figure 3-12 For UA series magnet only: Magnet Mounting Detail for M7000VA \(non-seismic\)—Patient End is on the Left Side on page 79](#) for details.
 - Weight: 8 kg (17 lb.) each
 - Size: 390 x 495 mm (15.25 x 19.25 in.) each
 - A 102 mm x 102 mm (4 in. x 4 in.) rebar-free area is necessary under the table dock anchor, in the position shown in the figure below.



IMPORTANT

The table dock anchor hole is drilled only after magnet installation.

Figure 3-12 For UA series magnet only: Magnet Mounting Detail for M7000VA (non-seismic)— Patient End is on the Left Side

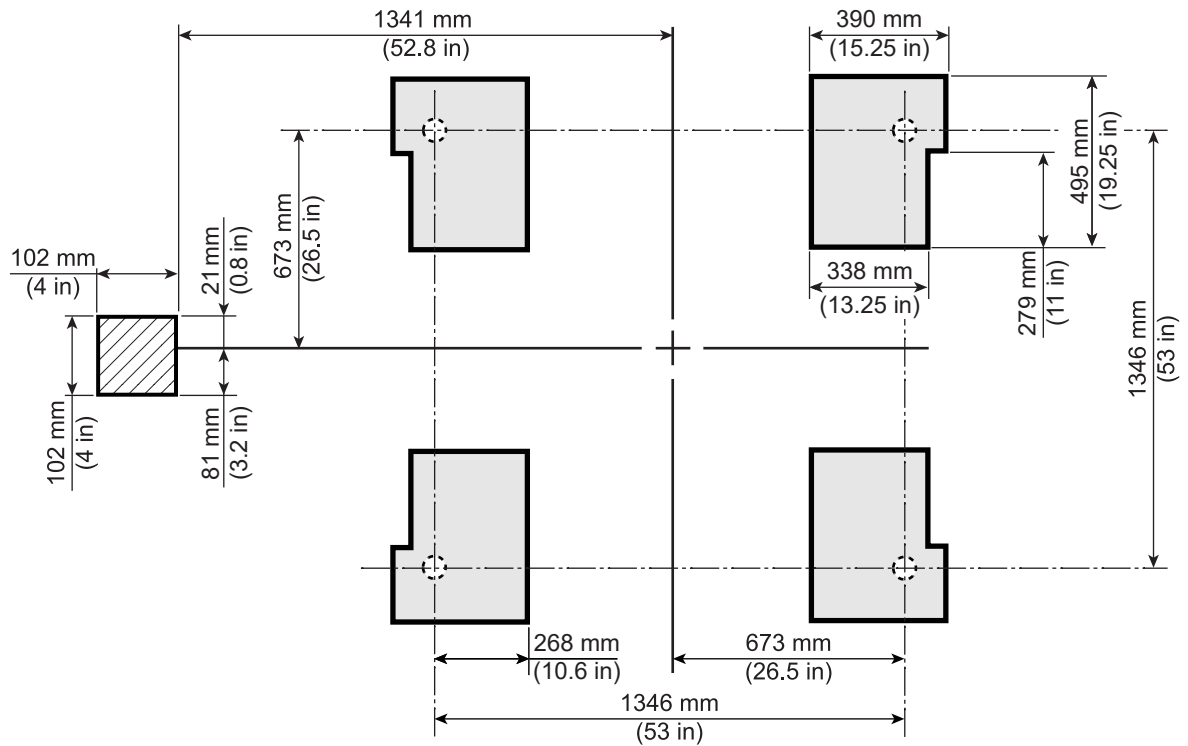

6. **For UA series magnet only:** The VibroAcoustic Dampening kit for systems that require Seismic mounting is M7000VM. See below for dimensions that aid in the calculation for floor loading. See [Figure 3-13 For UA series magnet only: Magnet Mounting Detail for M7000VM \(seismic\)— Patient End is on the Left Side on page 80](#) for details.


- 2 pads per magnet, 116 kg (255 lb.) each
- Nominal size: 610 x 1651 mm (24 x 65 in.) each
- Magnet anchor option 1: Corner seismic anchor points (x4)
- Magnet anchor option 2: Seismic anchor holes (x4) 38 mm (1.5 in.)
- Seismic anchors must be isolated from rebar.
- A 102 mm x 102 mm (4 in. x 4 in.) rebar-free area is necessary under the table dock anchor, in the position shown in the figure below.



IMPORTANT

The table dock anchor hole is drilled only after magnet installation.

Figure 3-13 For UA series magnet only: Magnet Mounting Detail for M7000VM (seismic)— Patient End is on the Left Side

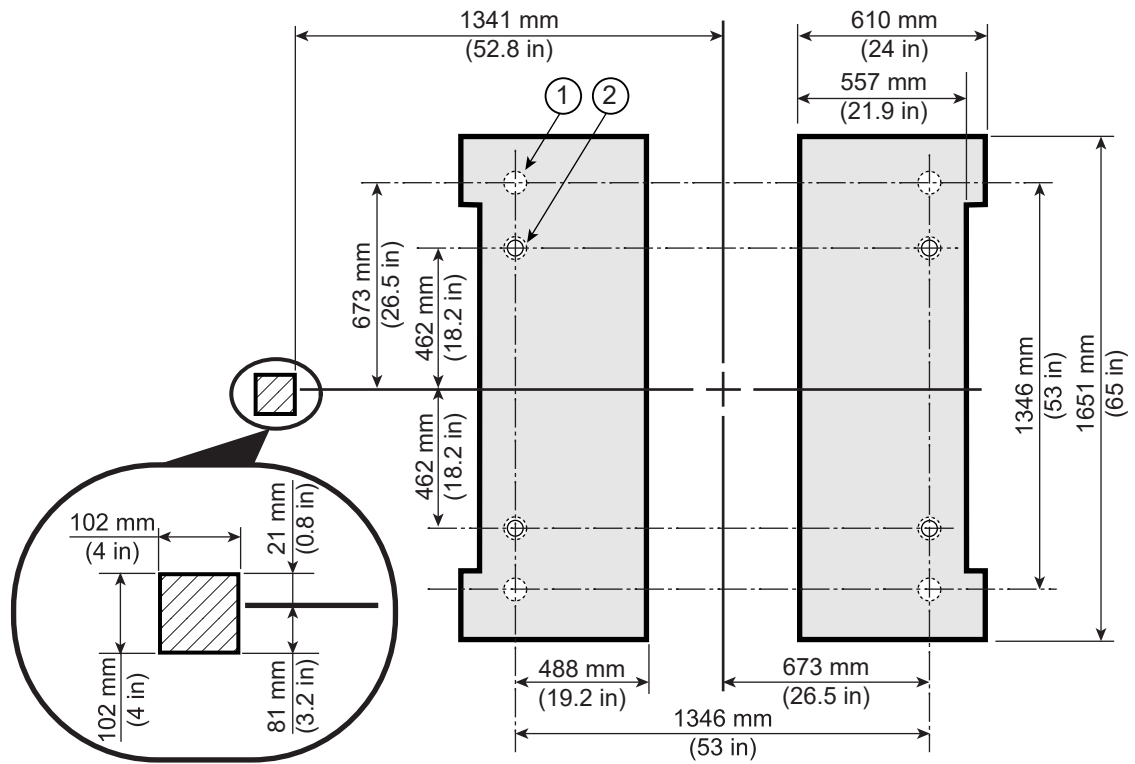

Item	Description	Item	Description
1	Magnet mounting holes	2	Seismic mounting holes

- For AR series magnet only:** The VibroAcoustic Dampening kit for both seismic and non-seismic mounting is included with the magnet shipment. See [Figure 3-14 For AR series magnet only: Magnet Mounting Detail— Patient End is on the Left Side](#) on page 81 for details.


 - Seismic anchors shall have a 76.2 mm (3 in) rebar-free area around the anchor and shall be installed with the assistance of the RF shield vendor.
 - A 102 mm x 102 mm (4 in. x 4 in.) rebar-free area is necessary under the table dock anchor, in the position shown in the illustration below.
 - Seismic anchors must be isolated from rebar.
 - For seismic installations, a mounting hardware kit for the magnet feet (M6001AH) is available from GE HealthCare.



IMPORTANT

The table dock anchor hole is drilled only after magnet installation.

Figure 3-14 For AR series magnet only: Magnet Mounting Detail— Patient End is on the Left Side

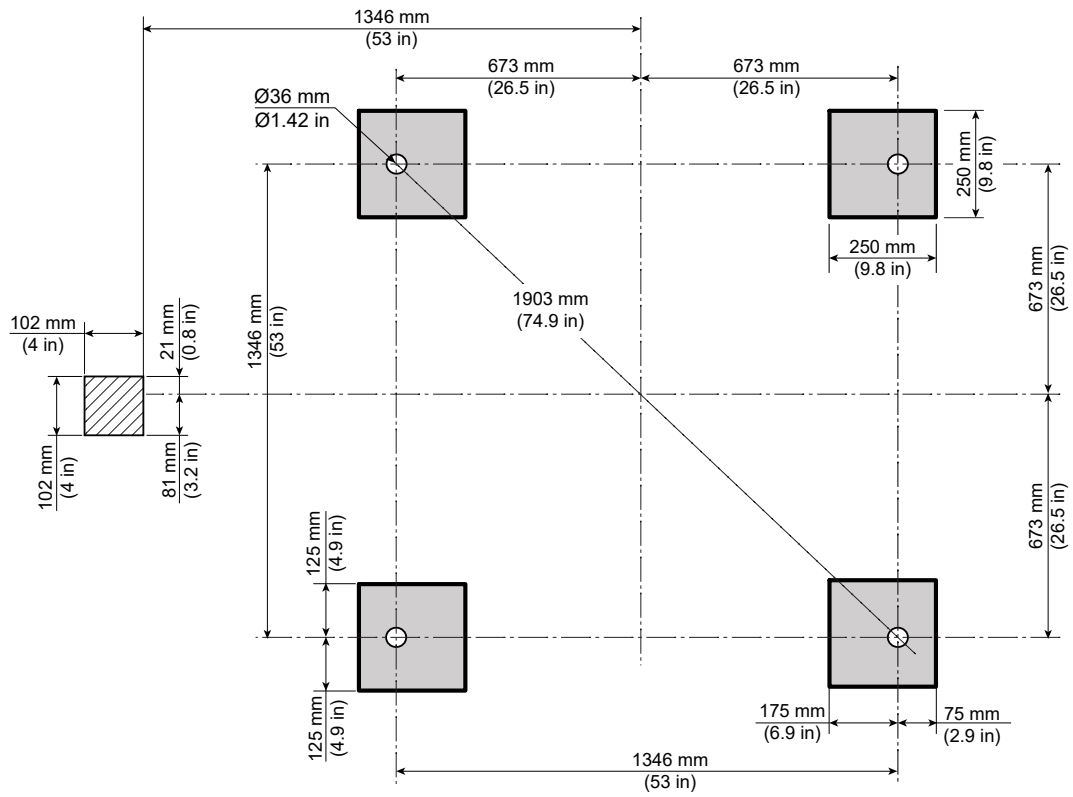

8. RF shield seams, joints, or overlaps must not be located under the VibroAcoustic mats.

3.5.9 Storage Cabinets



NOTE

GE HealthCare no longer provides a storage solution for system phantoms.

1. The customer shall provide storage for phantoms in the magnet room (for example, a cart, shelving unit or cabinet). Storage needs to be large enough to accommodate system phantoms listed in Table 1-1 of *Customer Site Storage Requirements*, 5182674 (available in the Customer Documentation Portal).
2. The storage solution can not interfere with the magnet room minimum service area.

3.6 Magnet Room Equipment Specifications




(Applies to all subsections within this section)

3.6.1 Magnet (MAG) Assembly Specifications


1. Rear Pedestal weight: 96 kg (212 lb.)
2. UA Magnet weight, with cryogenics at maximum capacity: 7443 kg (16410 lb.)
3. AR Magnet weight, with cryogenics at maximum capacity: 5665 kg (12489 lb.)

Figure 3-15 Magnet (MAG) Dimensions and Rear Pedestal for UA Magnet

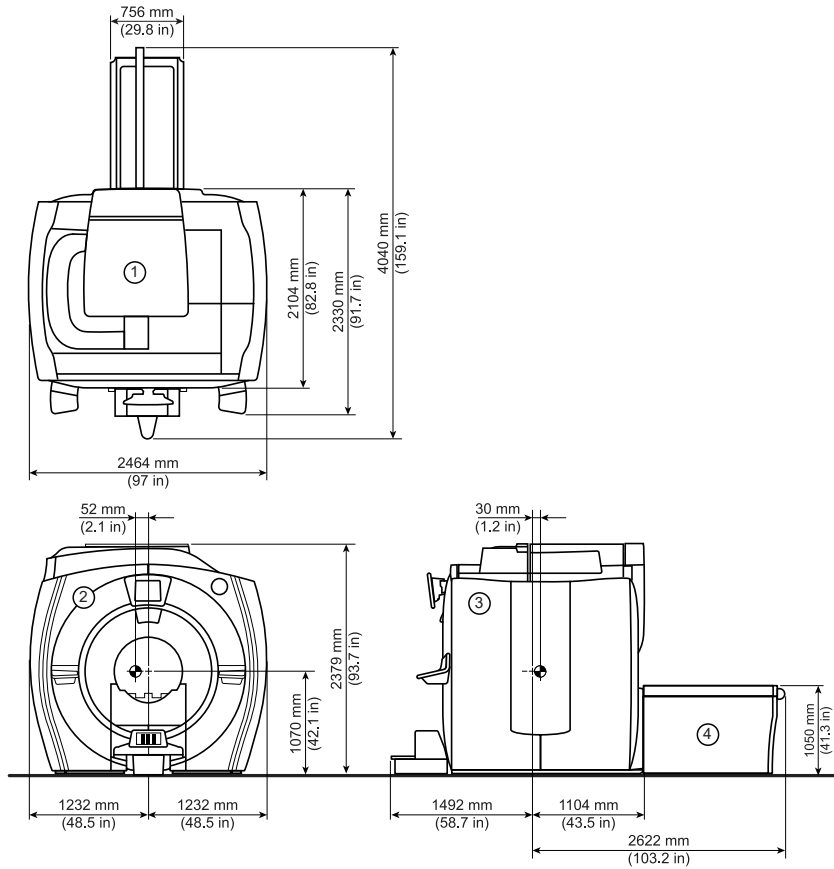

Item	Description	Item	Description
1	Top view	3	Side view
2	Front view	4	Rear Bridge

Figure 3-16 Magnet (MAG) Dimensions and Rear Pedestal for AR Magnet

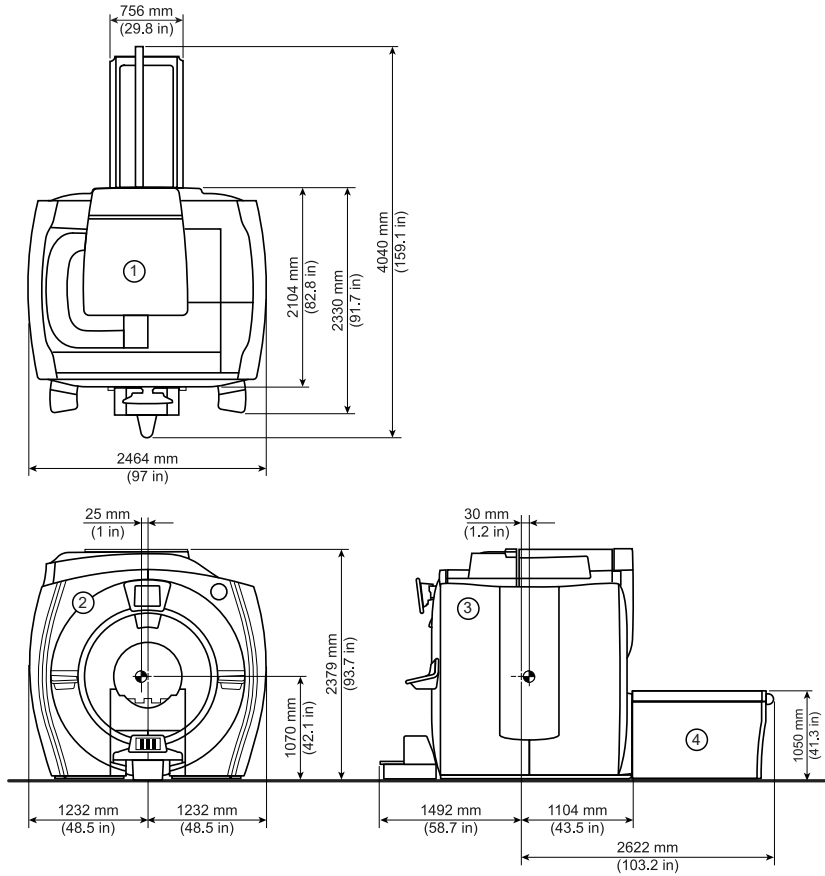


Item	Description	Item	Description
1	Top view	3	Side view
2	Front view	4	Rear Bridge

3.6.2 Patient Table (PT) Specifications

1. eXpress Table weight, empty: 210 kg (463 lb.)
2. eXpress Table weight with maximum patient weight of 227 kg (500 lb.): 437 kg (963 lb.)

Figure 3-17 GEM Patient Table (PT)

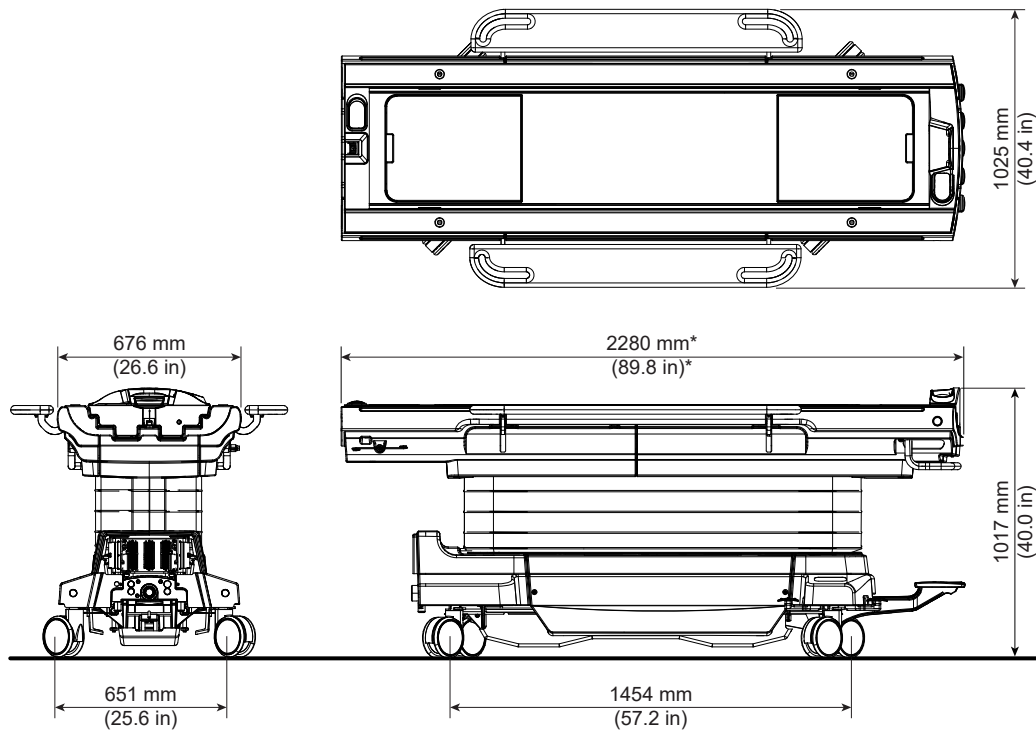


* When the table handle is rotated into its extended position, the length of the table is 2370 mm (93.3 in.)

3.6.3 Optional Surgical Suite Transmobile Table

1. Weight, empty: 164 kg (362 lb.)
2. Weight with maximum patient weight of 225 kg (496 lb.): 389 kg (858 lb.)

Figure 3-18 Optional Surgical Suite Transmobile Table

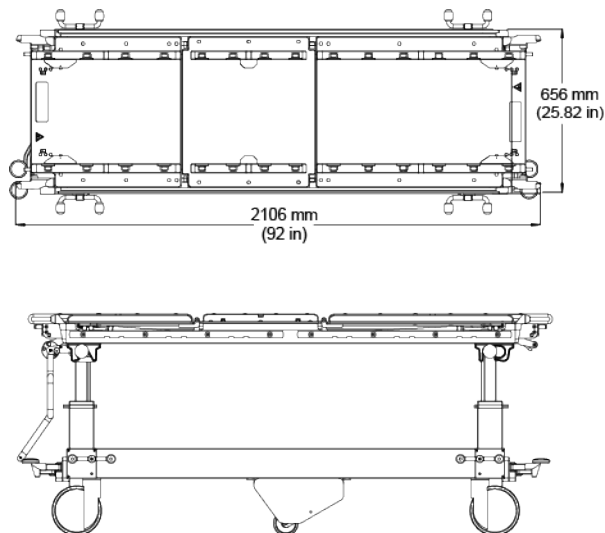

3.6.4 Magnet Rundown Unit (MRU) Specifications and Requirements


1. Location: The bottom edge of the MRU must be mounted 1524 ± 25 mm (60 ± 1 in.) above the Magnet Room floor near the front of the magnet enclosure.
2. Weight: 3.2 kg (7 lb.)
3. Magnetic Field Limit: 20 mT (200 G)
4. The MRU is installed by the facility contractor.



NOTE


An optional remote MRU may be installed outside the Magnet Room. The remote MRU does not require facility power. For more information, refer to MRU vendor manual.

Figure 3-19 Magnet Rundown Unit (MRU)

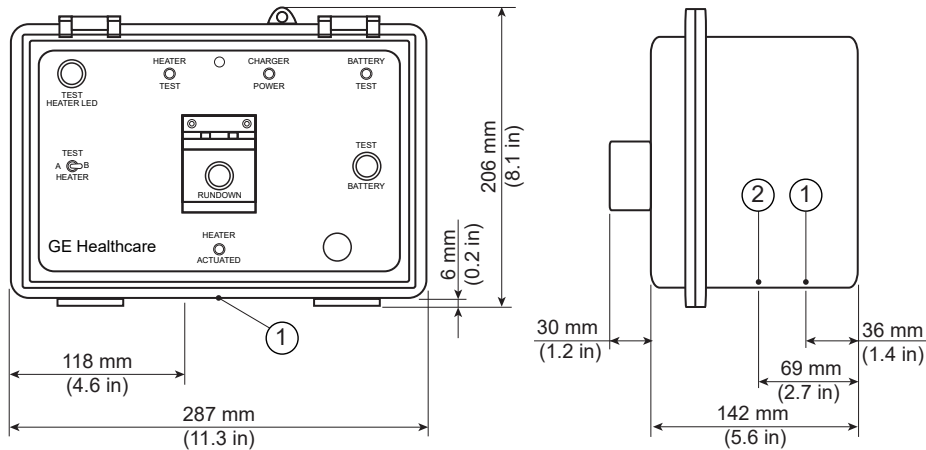

Item	Description
1	Cable access
2	Power access

Figure 3-20 MRU Mounting Pattern

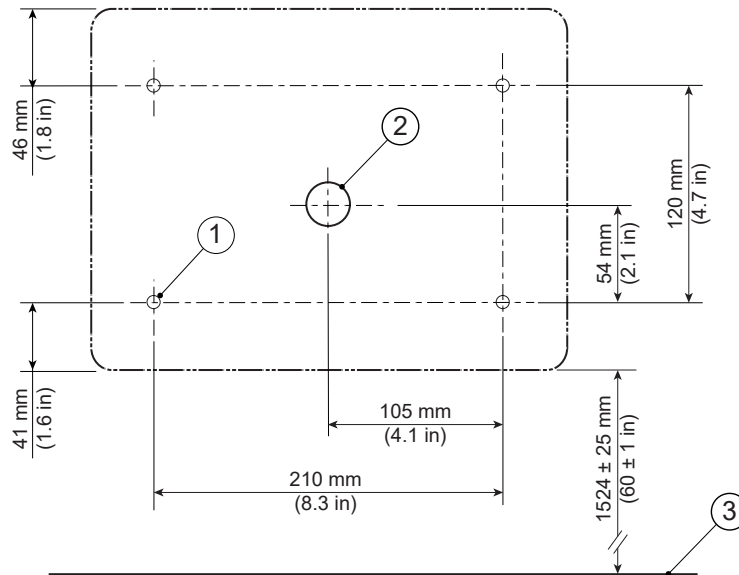

Item	Description	Item	Description
1	7 mm (0.275 in.) diameter mounting hole	3	Finished floor
2	26 mm (1.025 in.) diameter cable access	-	-

3.6.5 Oxygen Monitor Option Sensor Specifications

See [Oxygen Monitor \(OXY\) Option](#) on page 113.

3.7 Magnet Room Lighting Requirements



1. All lighting fixtures and associated components must meet all RF shielded room and RF grounding requirements (for example, track lighting is not recommended due to possible RF noise).
2. All removable lighting fixtures and associated components must be non-magnetic.
3. All lighting must use direct current (the DC must have less than 5% ripple).
4. At least 300 lux must be provided at the front of the magnet for patient access and above the magnet for servicing.
5. Fluorescent lighting must not be used in the Magnet Room.
6. Lighting must be adjusted using a discrete switch or a variable DC lighting controller.
7. SCR dimmers or rheostats must not be used.
8. DC LED lighting may be used if the DC power converter and RF sources are all located outside the Magnet Room RF Shield.



NOTE

LED lighting could cause image quality issues due to RF interference. Make sure a MR-compatible LED lighting solution is chosen.

9. Battery chargers (for example, used for emergency lighting) must be located outside the Magnet Room.
10. LED Lighting or short filament length incandescent bulbs are recommended.
11. Linear lamps are not recommended due to the high burnout rate.

4 Equipment Room

4.1 Equipment Room Overview



(Applies to all sections within this chapter)

1. The vertical distance between the coolant connection points of the ICC and the Gradient Coil must be less than 5 meters (196.8 in.).
2. The ICC and PGR Cabinets must be located on the same floor.



NOTE

The ICC contains the cryocooler compressor.

The following illustration shows minimum equipment room service clearances. Refer to [MR Suite Minimum Room Size Requirements on page 18](#) for a list of considerations not included in the minimum area dimensions.



NOTE

Colored areas indicate service/installation areas. These areas can overlap as necessary as shown below. See individual component descriptions and room requirements for service area details. Magnet Monitor (MON) can be mounted on the side of the wall near the ICC. Optional equipment is not shown; additional space may be required for options.



NOTE

A 381 mm (15 in.) clearance on the cabinet sides for PGR, ICC, PEN is needed for cabinet dolly removal. This clearance must extend from the cabinet front for dolly removal if the adjacent cabinet or wall is fixed.



NOTE

Refer to [Table 2-21 MR System Component Replacement Shipping Specifications on page 60](#) for the dimensions of the replacement parts. The parts must be able to be positioned in front of the noted cabinet for replacement procedures after the system has been installed.

Figure 4-1 Typical Minimum Equipment Room with Service Clearances (with ICC)


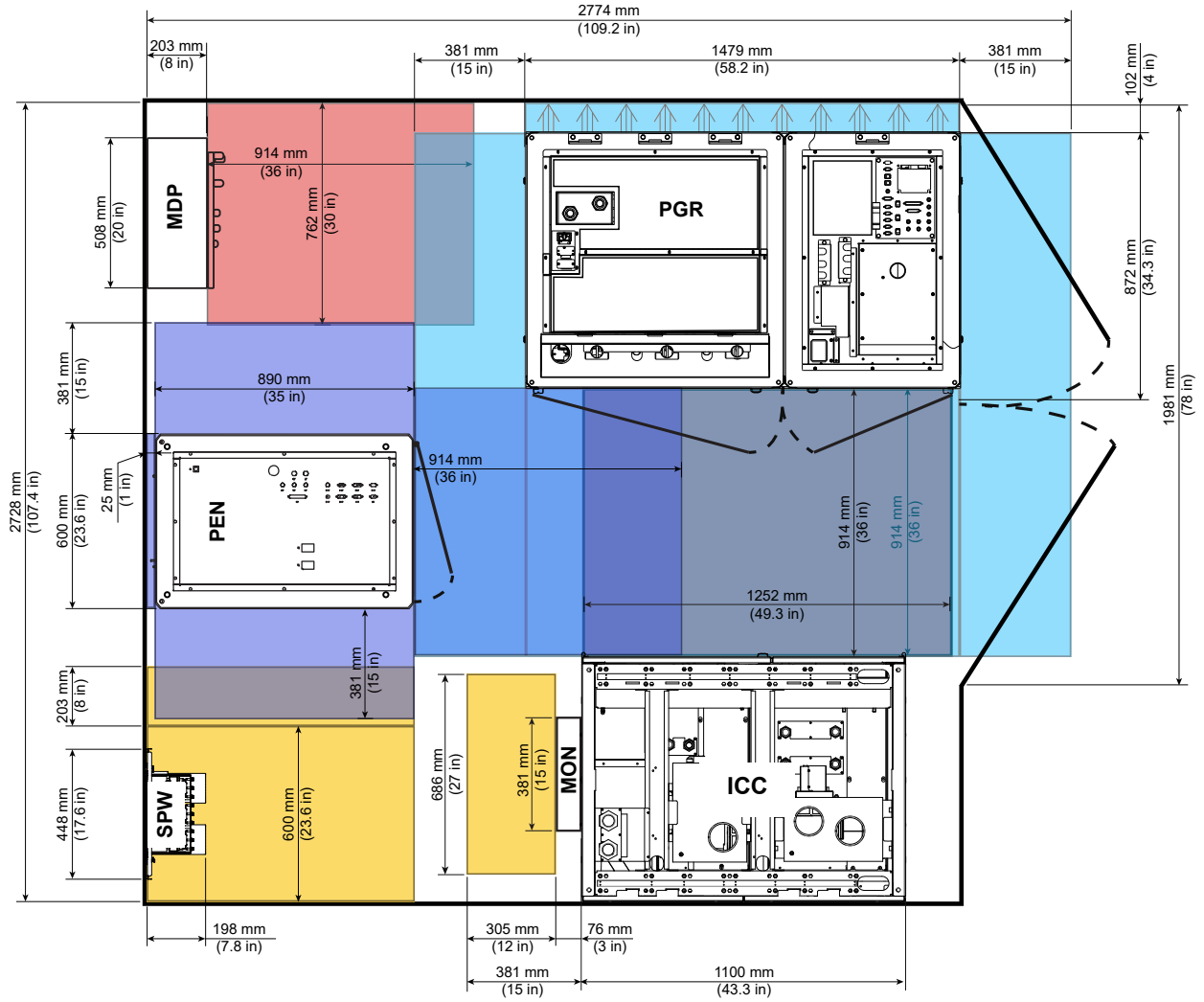


Figure 4-2 Typical Minimum Equipment Room with Service Clearances (with pICC)

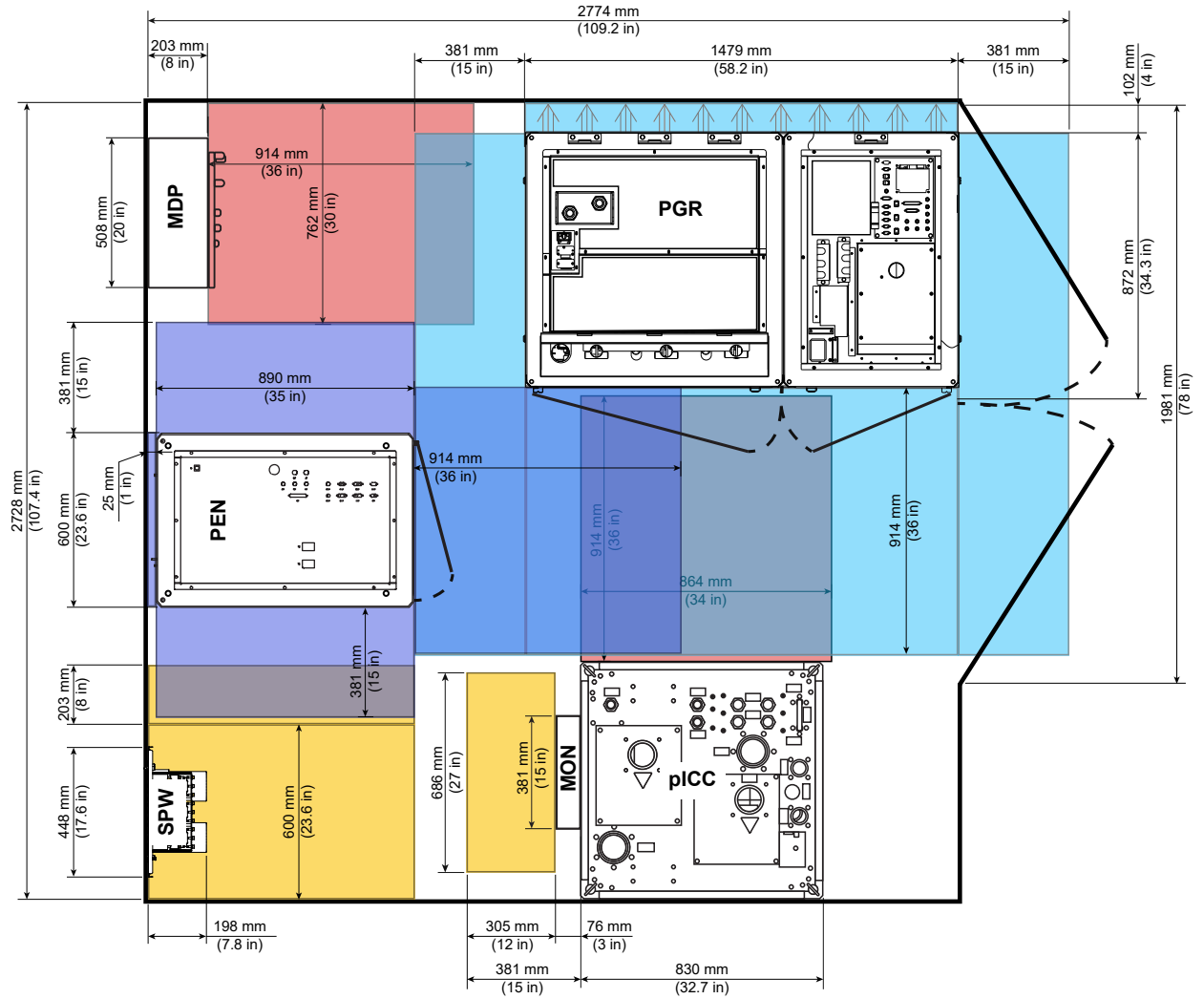

4.2 Main Disconnect Panel (MDP) Requirements and Specifications

4.2.1 Requirements

1. It is recommended to install the following items to support a power monitor:
 - 1.1. A T100 network connection with RJ45 connector near the MDP
 - 1.2. An electrical outlet
2. The cable must be Cat 5 or better.
3. The network connection must not be routed through the Ethernet switch in the Global Operator Cabinet (GOC).


4.2.2 Specifications

The Main Disconnect Panel (MDP) is provided with the MR System. Only exempt countries can supply their own MDP (see [2.10.3 Customer-supplied Main Disconnect Panel \(MDP\) Requirements \(exempt countries only*\)](#) on page 52).

M7100ZA and M7100ZB

1. Weight: 57 kg (125 lb.)
2. Magnetic Field Limit: 5 mT (50 G)

Figure 4-3 GE HealthCare supplied Main Disconnect Panel (MDP) M7100ZA and M7100ZB

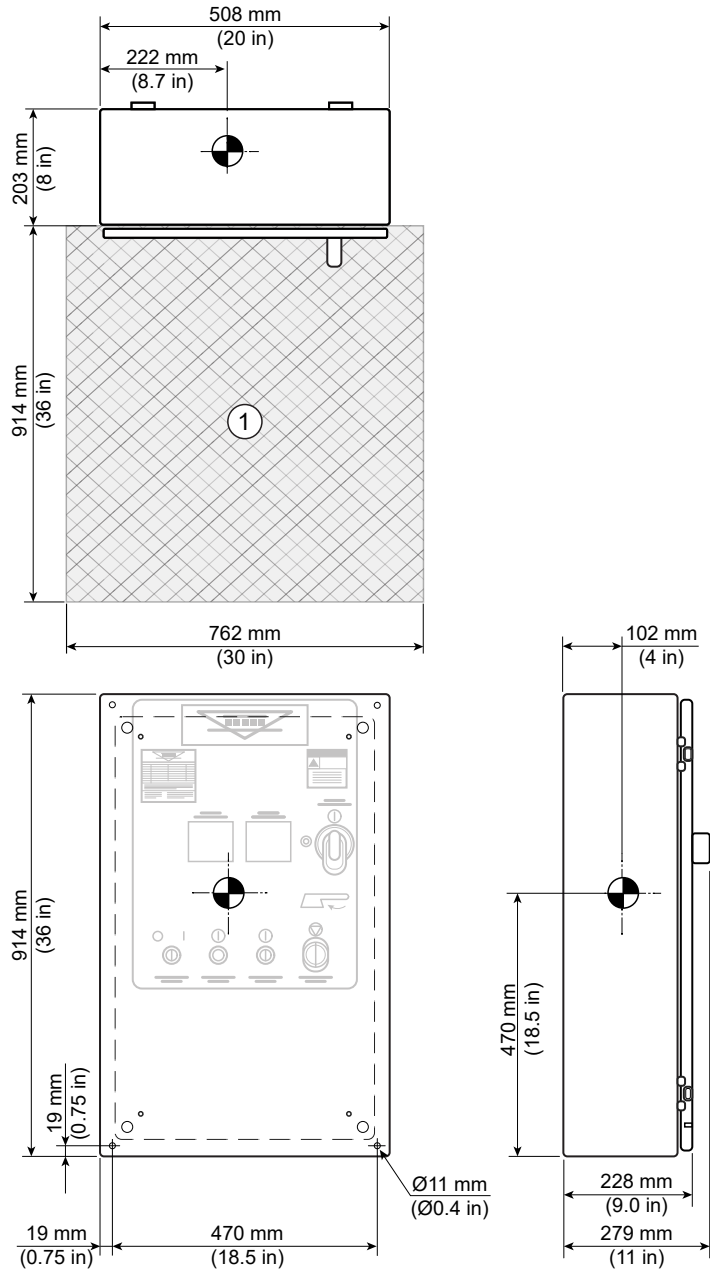


Item	Description
1	Service Clearance

4.3 Power, Gradient, RF Cabinet (PGR) Specifications

1. Weight: 1510 kg (3329 lb.)
2. Magnetic Field Limit: 5 mT (50 G)

Figure 4-4 Power, Gradient, RF (PGR) Cabinet

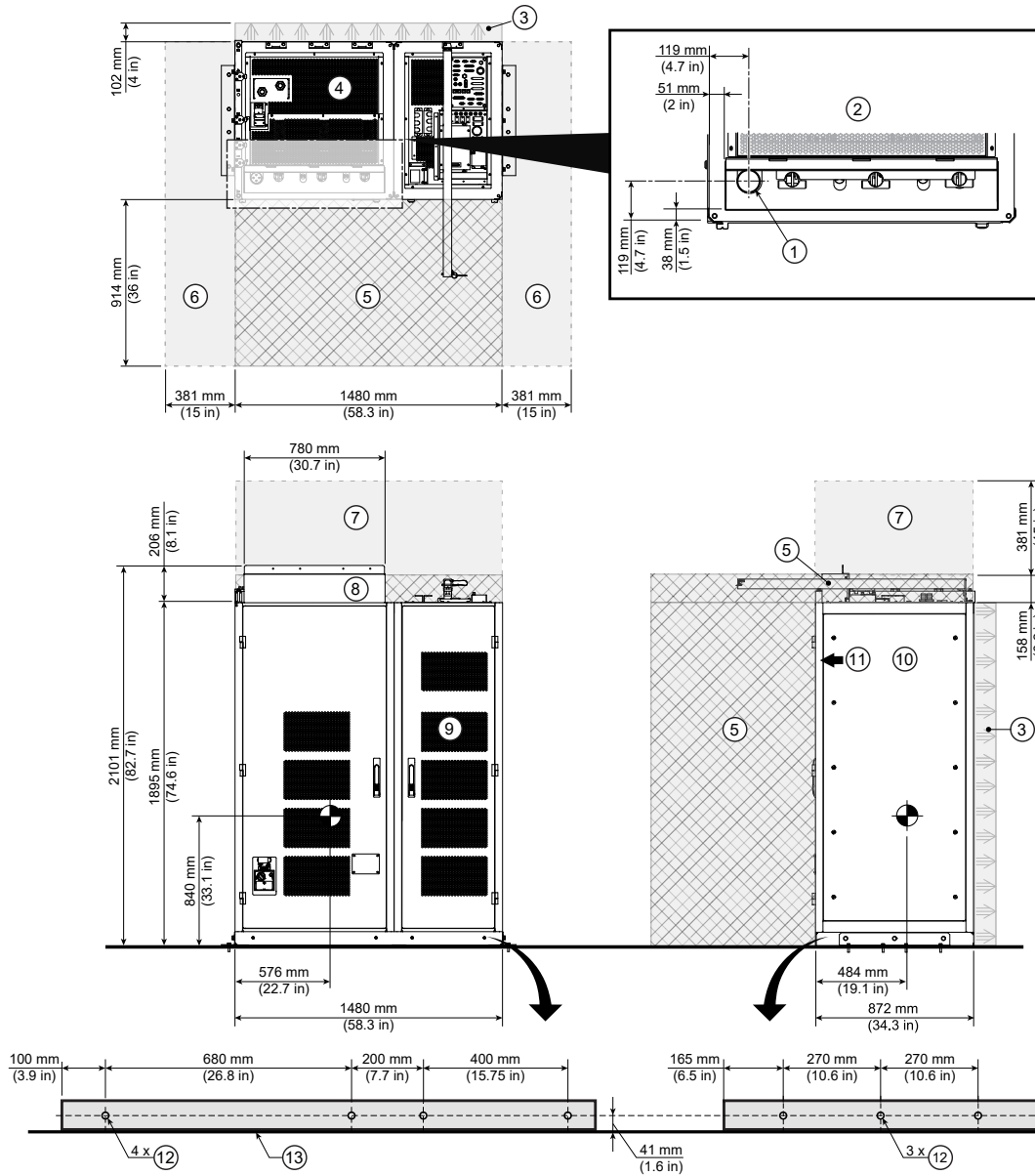


Item	Description	Item	Description
1	Conduit Hole for Incoming Power from MDP	8	Cable Strain Relief
2	Strain Relief Detail (Top View)	9	Front View
3	Airflow Clearance	10	Side View
4	Top View	11	Cabinet Front
5	Service Clearance	12	Seismic Anchor Mounting Holes to accommodate M12 fasteners (Same on front and both sides of cabinet)
6	Dolly Installation Clearance	13	Finished floor
7	Cable Clearance		

4.4 Penetration Cabinet (PEN) Specifications

The Penetration Cabinet Penetration Panel provides interconnects from the PEN Cabinet through the Magnet Room RF Shield.

1. Weight: 313 kg (691 lb.)
2. PEN Panel Magnetic Field Limit: 20 mT (200 G) for the entire PEN Panel
PEN Cabinet Magnetic Field Limit: 5 mT (50 G)
3. The PEN cabinet must be positioned directly in front of the PEN Panel. Refer to [3.4 Penetration Panel Wall Opening Requirements on page 67](#) for PEN panel mounting and location requirements.

Figure 4-5 PEN Cabinet

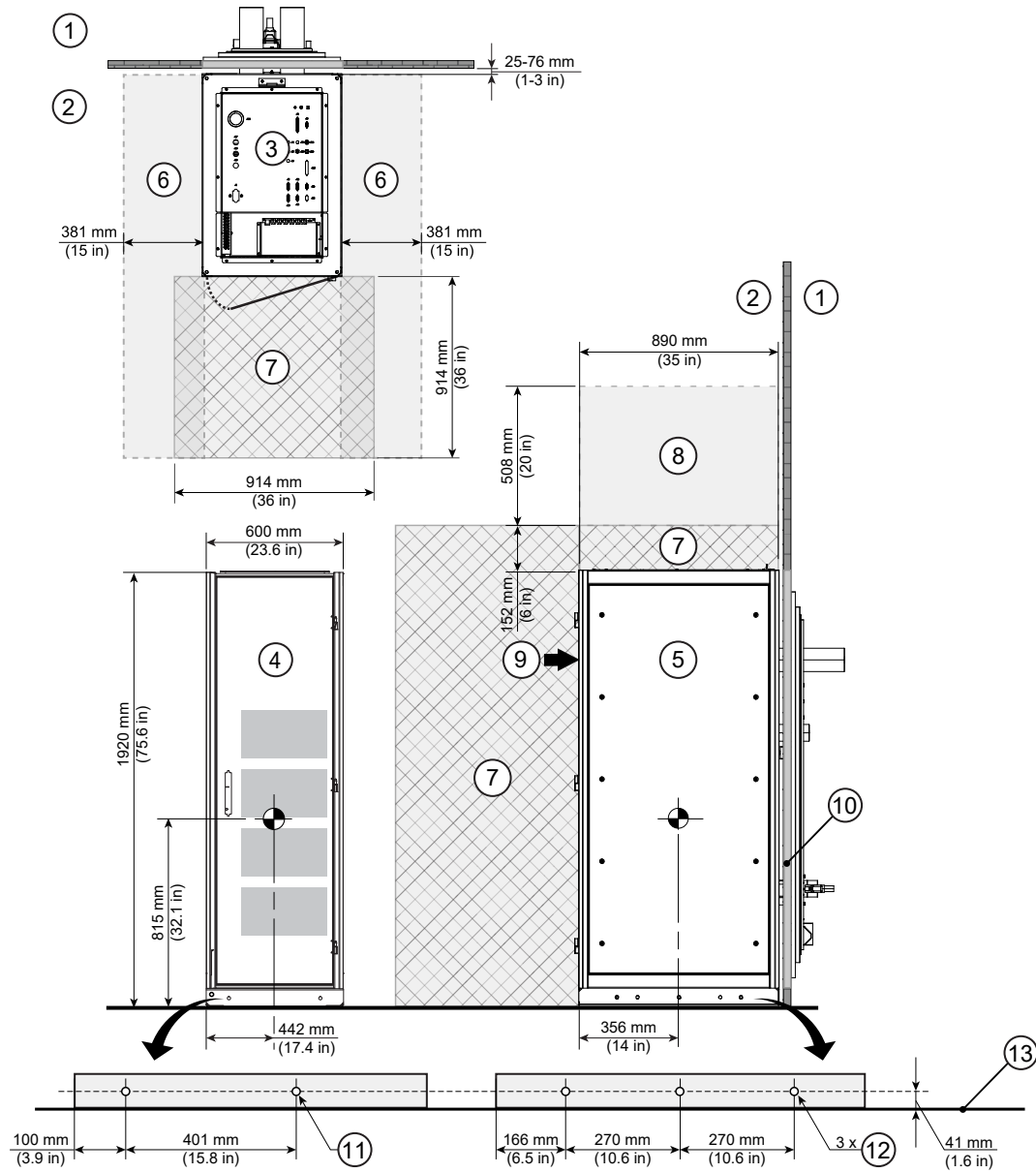


Item	Description	Item	Description
1	Magnet Room	8	Cable Drop Area
2	Equipment Room	9	Cabinet Front
3	Top View	10	Penetration Panel (See 3.4 Penetration Panel Wall Opening Requirements on page 67)
4	Front View	11	All M12 x 1.75 mm Front Mounting Holes
5	Side View	12	Seismic Anchor Mounting Holes (Same on both sides of cabinet)
6	Dolly Installation Clearance	13	Finished floor
7	Service Clearance	-	-

4.5 Integrated Cooling Cabinet (ICC) Specifications


1. ICC Weight with coolant, but without Cryocooler Compressor: 615 kg (1353 lb.)
2. pICC Weight with coolant, but without Cryocooler Compressor: 459 kg (1011 lb.)
3. F-50SH Cryocooler Compressor weight: Approx 125 kg (276 lb.)
4. Magnetic Field Limit: 5 mT (50 G) (Penetration Panel side)
5. The threaded seismic anchor mounting holes accommodate M12 bolts.
6. Seismic mounting of the ICC must be performed at the sides of the cabinet.

Figure 4-6 Integrated Cooling Cabinet (ICC)

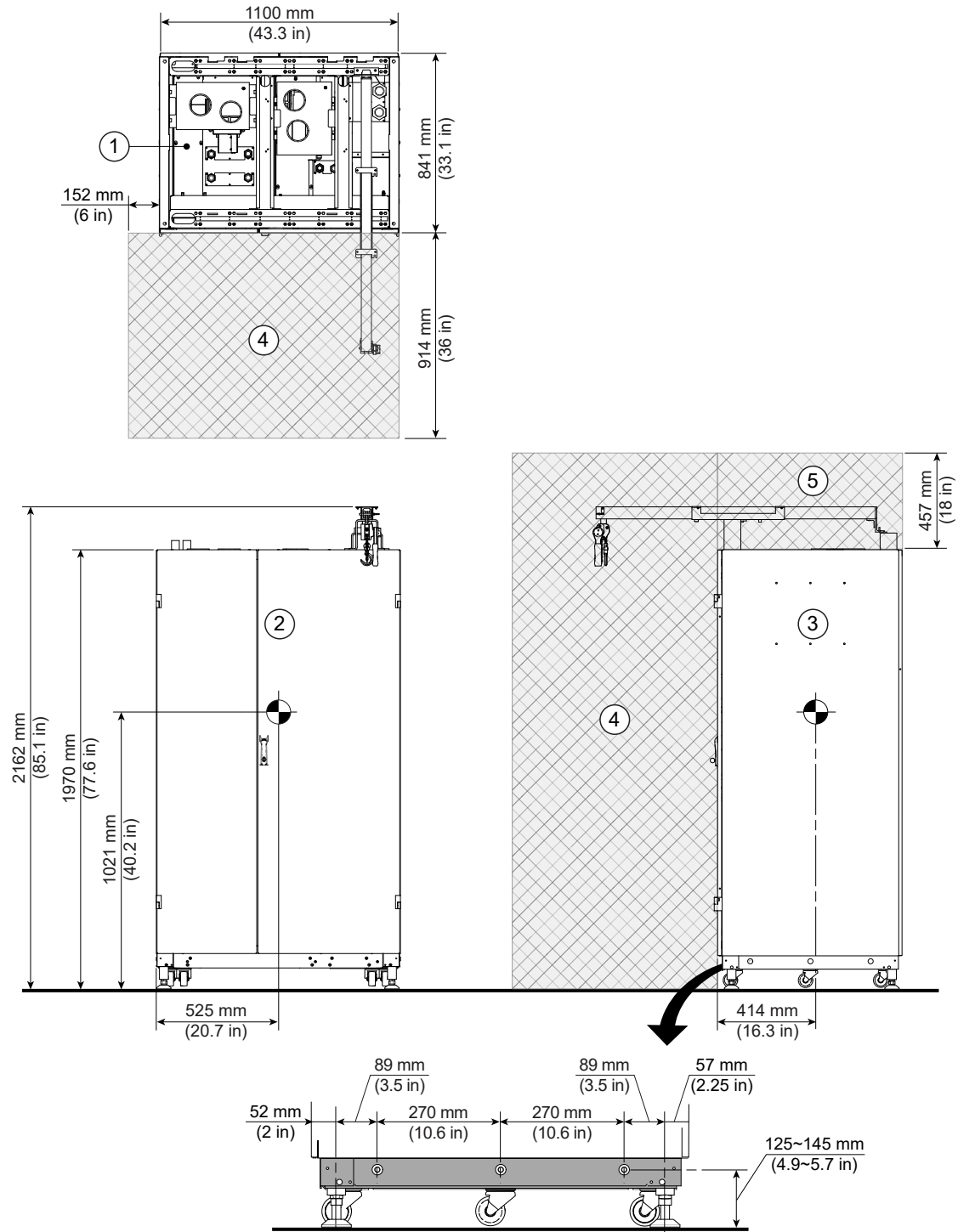


Item	Description	Item	Description
1	Top View	4	Service Area
2	Front View	5	Air Flow and Cables
3	Side View	-	-

Figure 4-7 Platform Integrated Cooling Cabinet (pICC)

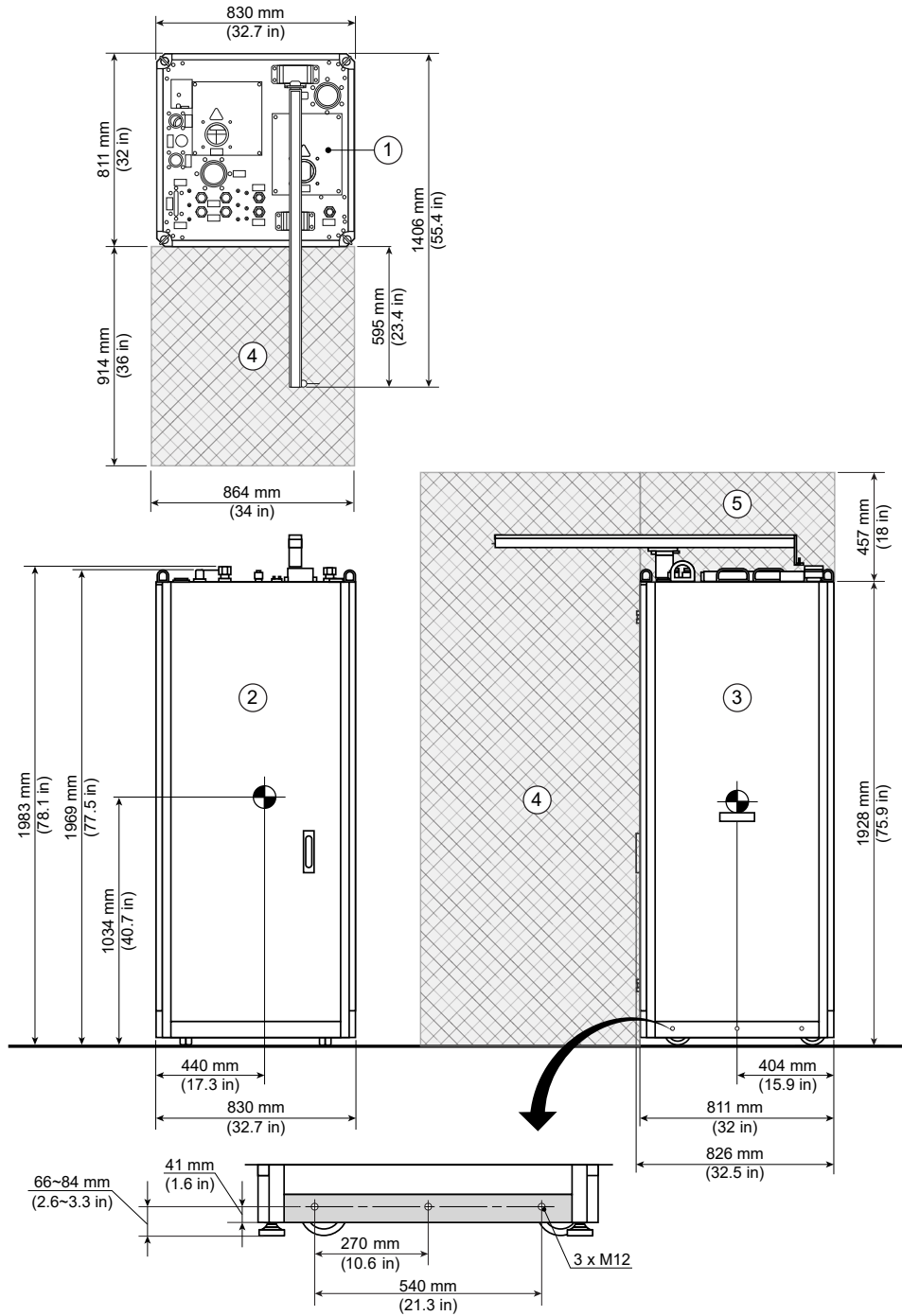

Item	Description	Item	Description
1	Top View	4	Service Area
2	Front View	5	Air Flow and Cables
3	Side View	-	-

Figure 4-8 pICC Facility Water Supply and Return Connections Location

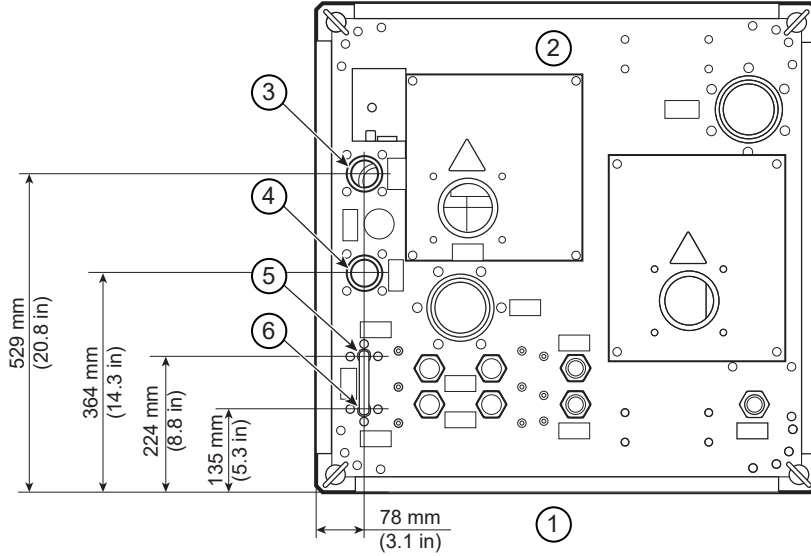


Item	Description	Item	Description
1	Cabinet Front	4	Facility Supply
2	Top View	5	City Water Return for Cryocooler compressor backup
3	Facility Return	6	City Water Supply for Cryocooler compressor backup

4.6 Secondary Penetration Wall (SPW) Specifications

The Secondary Penetration Wall (SPW) provides interconnects from the Equipment Room through the Magnet Room RF Shield.

1. Magnetic Field Limit: 20 mT (200 G)
2. Refer to [3.4 Penetration Panel Wall Opening Requirements](#) on page 67 for PEN panel mounting and location requirements.
3. The service area for the SPW is 914 mm (36 in.) from both sides the RF wall, and from floor to cable trays.

Figure 4-9 Secondary Penetration Wall (SPW)

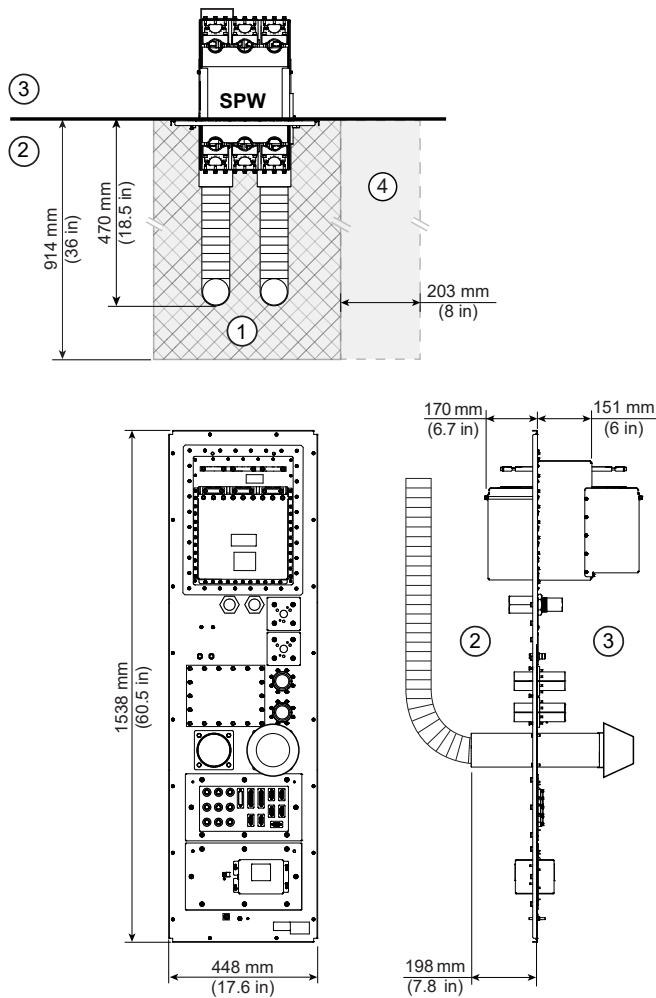

Item	Description		
1	Service clearance (shaded area)	3	Magnet Room
2	Equipment Room	4	Additional service clearance needed. The clearance can be on either side of the SPW.

4.7 Magnet Monitor (MON) Requirements and Specifications


4.7.1 Requirements

1. Customer must supply T100 network connection with RJ45 connector to the Magnet Monitor (MON). Network connectivity must be active prior to magnet delivery.
2. The cable must be Cat 5 or better.
3. The network connection must not be routed through the Ethernet switch in the Global Operator Cabinet (GOC).

4.7.2 Specifications


1. Mounting location: On either side of the Integrated Cooling Cabinet (ICC) or on a wall that is close to the ICC
2. Weight: 4.5 kg (10 lb.)
3. Magnetic Field Limit: 20 mT (200 gauss)
4. Power cord length: 1829 mm (72 in.)

Figure 4-10 Magnet Monitor (MON)

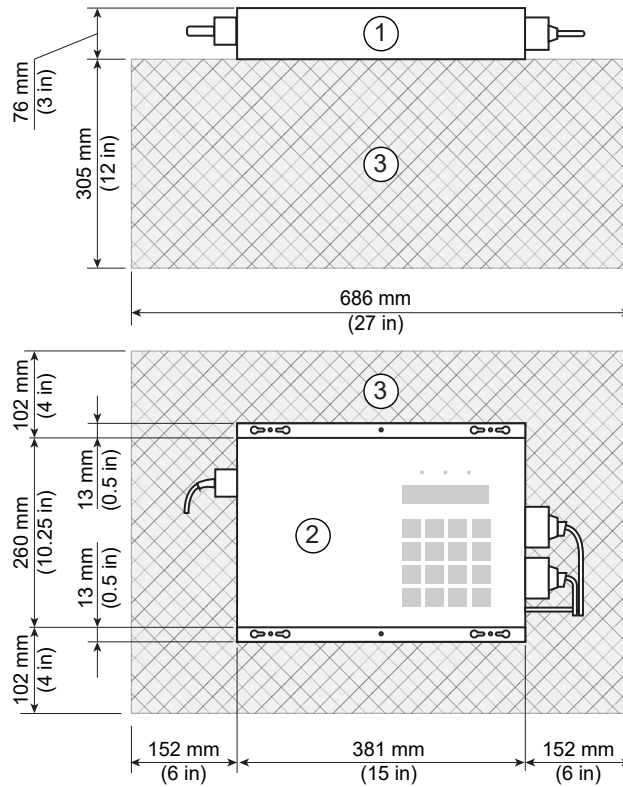

Item	Description	Item	Description
1	Top View	3	Service area
2	Front View	-	-

Figure 4-11 Magnet Monitor (MON) Mounting Patterns

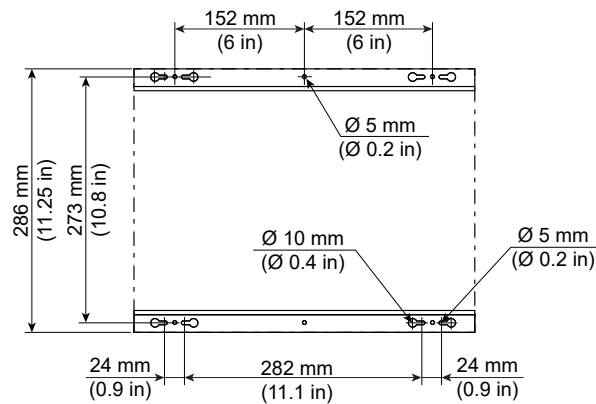

4.8 Magnetic Resonance Elastography (MRE) Specifications (Optional Equipment)


4.8.1 Requirements

1. The customer must work with the RF shield vendor to provide a waveguide for the 25 mm (1 in.) diameter tube.
2. MRE Resoundant Acoustic Driver location is limited to the length of the 25 mm (1 in.) tube (see the available cable lengths in [7.1.5 Magnetic Resonance Elastography \(MRE\) Option on page 121](#)).

4.8.2 Specifications

1. Weight: 24.22 kg (53.4 lb.)
2. Magnetic Field Limit: 5 mT (50 G)
3. Power Cord Length:
 - 60 Hz: 6096 mm (240 in.)
 - 50 Hz: 7620 mm (300 in.)


Figure 4-12 Magnetic Resonance Elastography (MRE) Resoundant Acoustic Driver

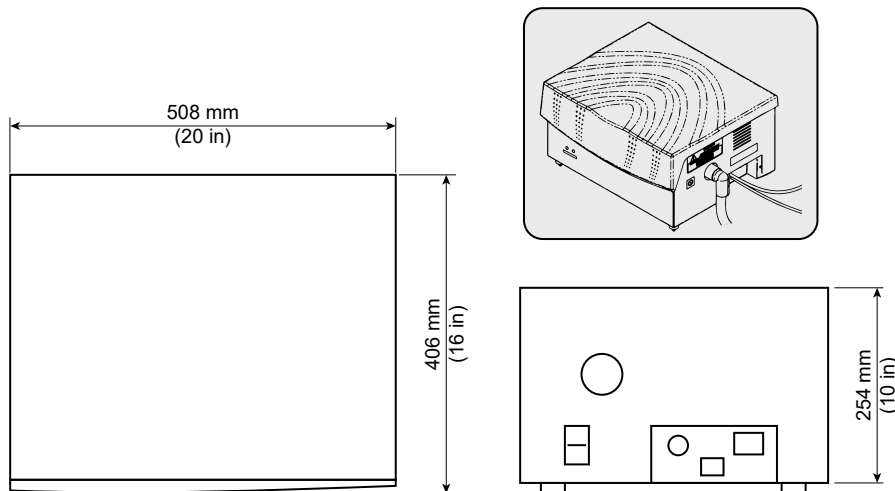

4.9 Multi-Nuclear Spectroscopy (MNS) Specifications (Optional Equipment)

1. Weight: 295 kg (650 lb.)
2. Magnetic Field Limit: 20 mT (200 G)
3. The threaded seismic anchor mounting holes accommodate 3/8"-16 bolts.

Figure 4-13 Seismic Mounting Detail for MNS Amplifier Cabinet for 8 kW MNS Option

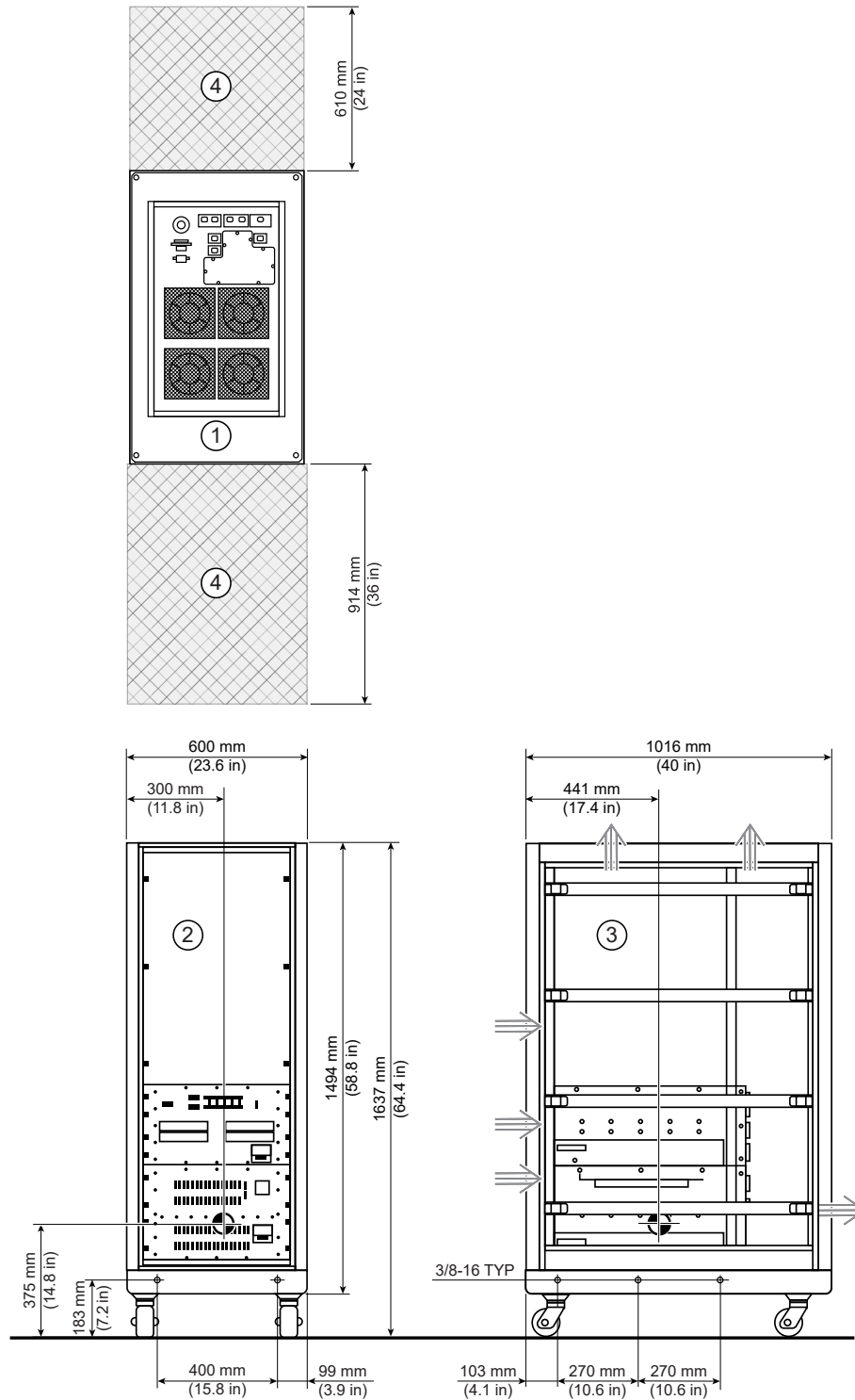

Item	Description	Item	Description
1	Top View	3	Side View
2	Front View	4	Service area

5 Control Room


5.1 Operator Workspace Equipment Specifications



(Applies to all sections within this chapter)


5.1.1 Operator Workspace Assembly

Figure 5-1 Operator Workspace Assembly

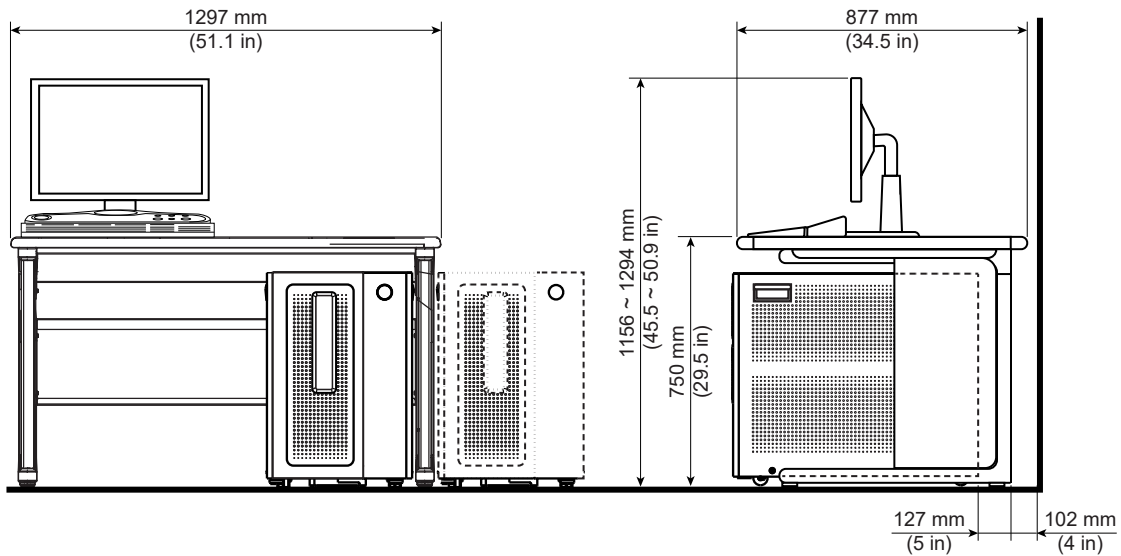

5.1.2 Operator Workspace (OW) (Optional Equipment)

Figure 5-2 Operator Workspace (OW) Table (Top and Front View)

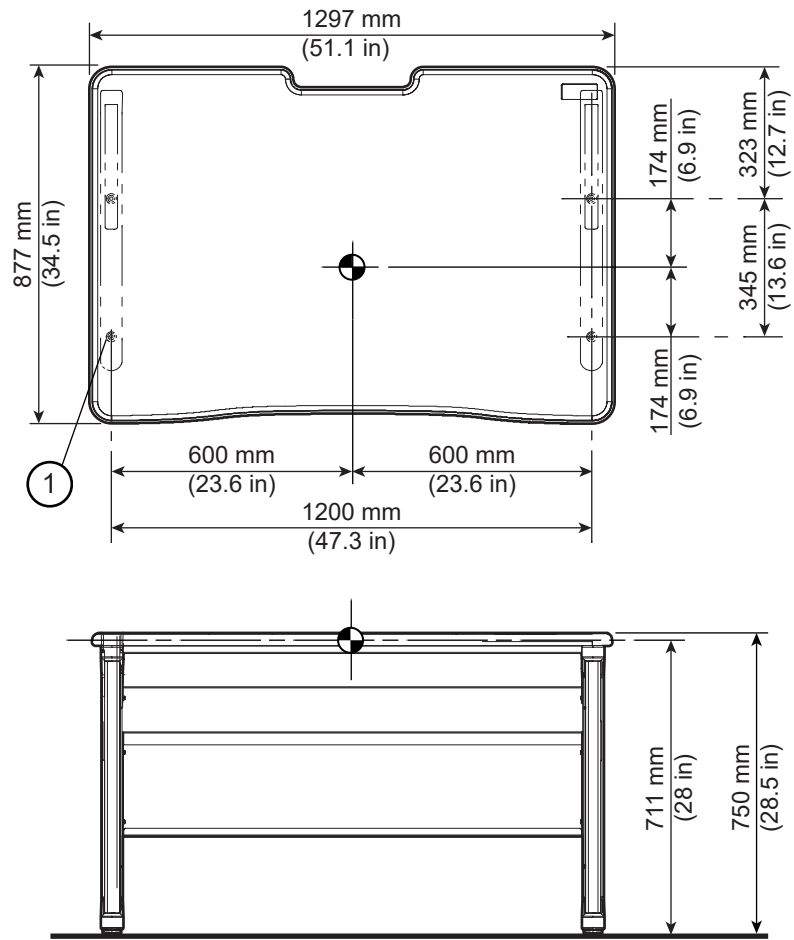


Item	Description
1	Four (4) 15.9 mm (5/8 in.) thru mounting holes for 9.5 mm (3/8 in.) seismic anchors

5.1.3 Global Operator Cabinet (GOC)

1. Weight:
 - For Dell T5820:** 60.4 kg (133.1 lb.)
 - For HP Z4G5:** 55.0 kg (121.3 lb.)
2. Magnetic Field Limit: 5 mT (50 G)
3. Anchor size: M10 (3/8 in.)

Figure 5-3 Global Operator Cabinet (GOC)

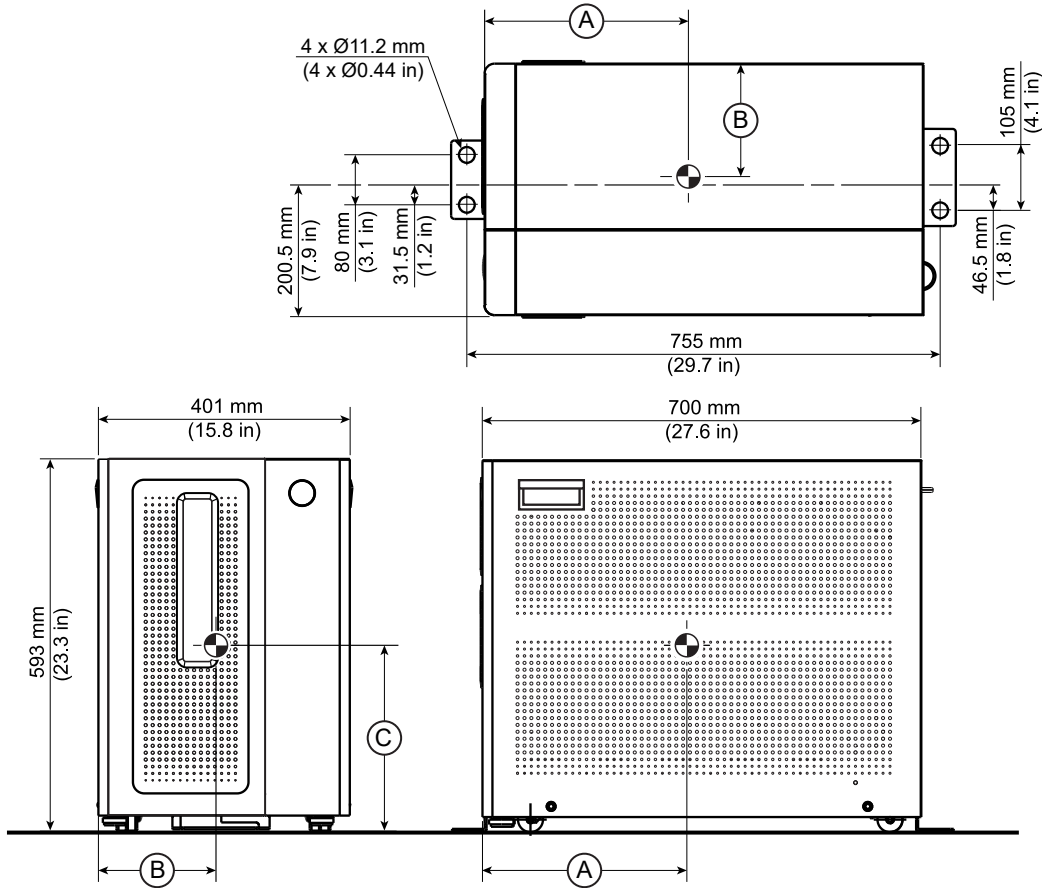


	A	B	C
Dell T5820	324 mm (12.7 in.)	192 mm (7.6 in.)	279 mm (11 in.)
HP Z4G5	375 mm (14.8 in.)	205 mm (8.1 in.)	330 mm (13 in.)

5.1.4 Host Display

Weight and dimensions for the Host Display are approximate and might vary depending on the display model.


- 1. Weight: 6.7kg (14.8 lb.)
- 2. Magnetic Field Limit: 5 mT (50 G)

Figure 5-4 Host Display

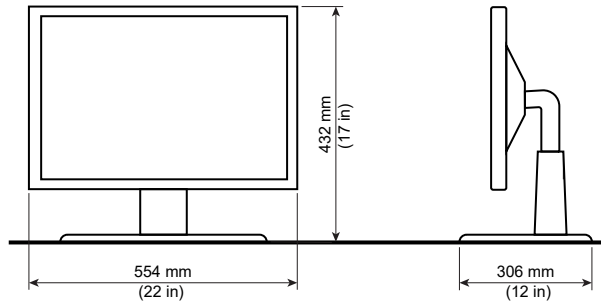

5.1.5 Host Keyboard

Figure 5-5 Host Keyboard

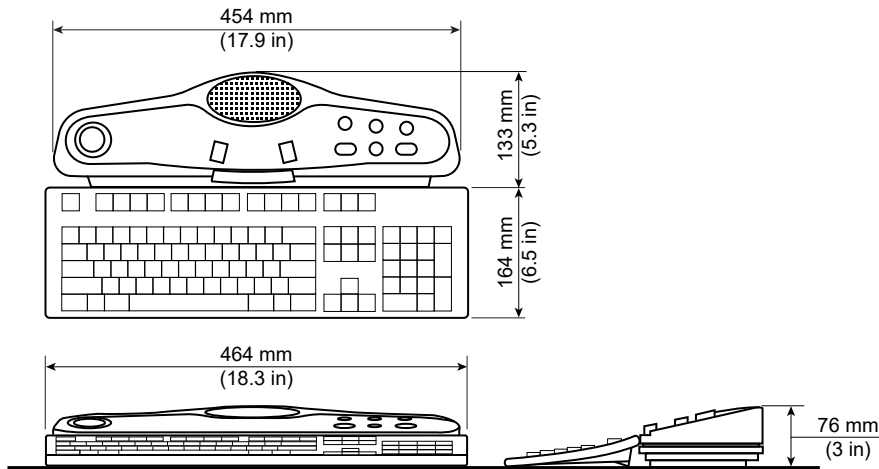


5.1.6 Pneumatic Patient Alert

The Pneumatic Patient Alert system allows the patient to contact the operator. The Control Box audible and visual alarm will be activated by the patient squeeze bulb which is located on the Magnet Enclosure and connected by pneumatic tubing through the Penetration Panel to the Control Box.

- 1. Weight 0.2 kg (0.5 lb.)
- 2. Magnetic Field Limit: 5 mT (50 G)
- 3. The Control Box must be placed or mounted within reach of the operator and within 1.5 m (5 ft.) of an electrical outlet.

Figure 5-6 Pneumatic Patient Alert (PA) with Mounting Pattern

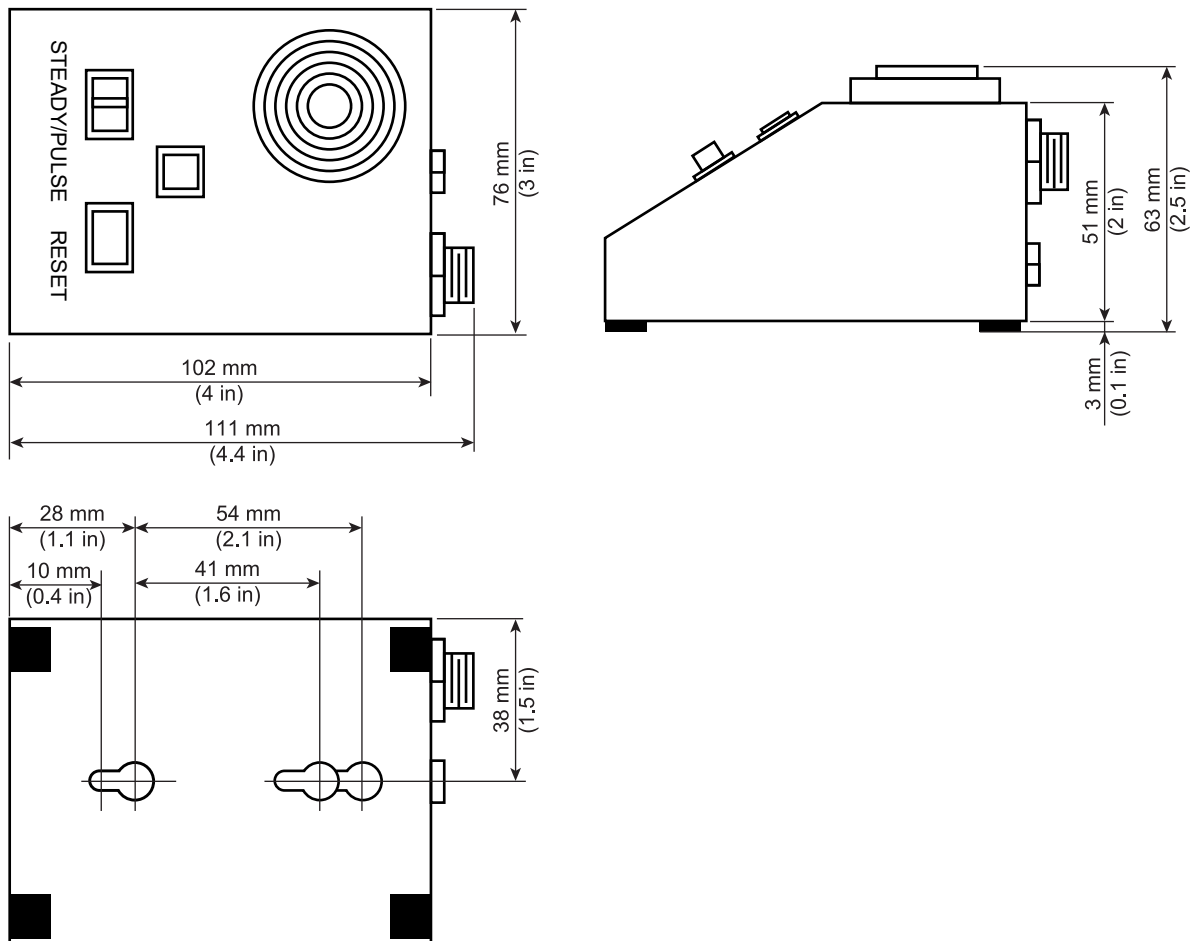


5.2 Oxygen Monitor (OXY) Specifications (Optional Equipment)

The optional Oxygen Monitor system consists of the Oxygen Monitor, the Remote Oxygen Sensor Module, and interconnects through the Secondary Penetration Wall (SPW). The Oxygen Monitor alarm is located near the Operator Workspace and is activated by the Remote Oxygen Sensor Module in the Magnet Room.

1. Oxygen Monitor Weight: 3.6 kg (8 lb.)
2. Oxygen Sensor Module Weight: 0.9 kg (2 lb.)
3. Magnetic Field Limit: 5 mT (50 G)

Figure 5-7 Oxygen Monitor and Remote Sensor

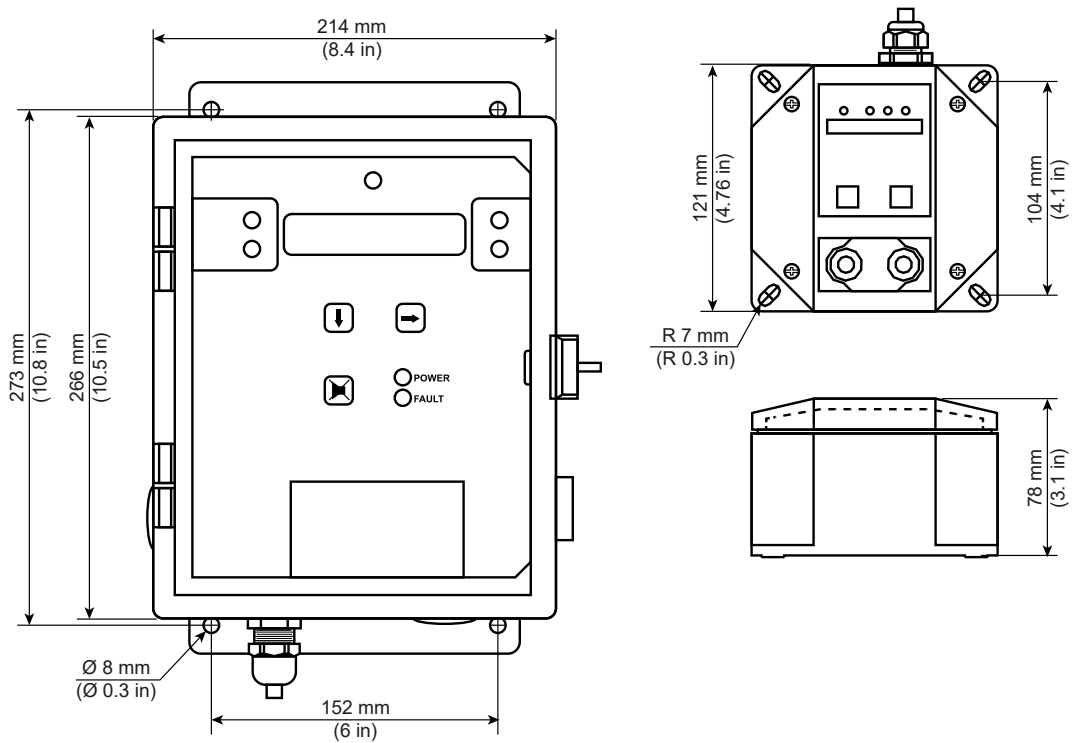


5.3 Physiological Acquisition Transceiver (PAT) Charging Station Specification (Optional Equipment)

The charging station supports the docking and charging of the Physiological Acquisition Transceiver (PAT). It is intended for installing outside the MR scanner magnet room, typically in the Control room. The charging station provides hooks for draping the ECG cable and respiratory sensor cable to facilitate cable management.

1. The charging cradle can be placed on the desktop or mounted on the wall.
2. A clean, level surface tabletop is required for desktop placement method.
3. Charging station wall mounted configuration maximum weight: 4.2kg (9.3 lbs)
4. Magnetic Field Limit: 5 mT (50 G)
5. The charging cradle must be placed or mounted within 3000 mm (118 in.) of an electrical outlet.

Figure 5-8 Charging Cradle Desktop and Wall Mounting Configuration

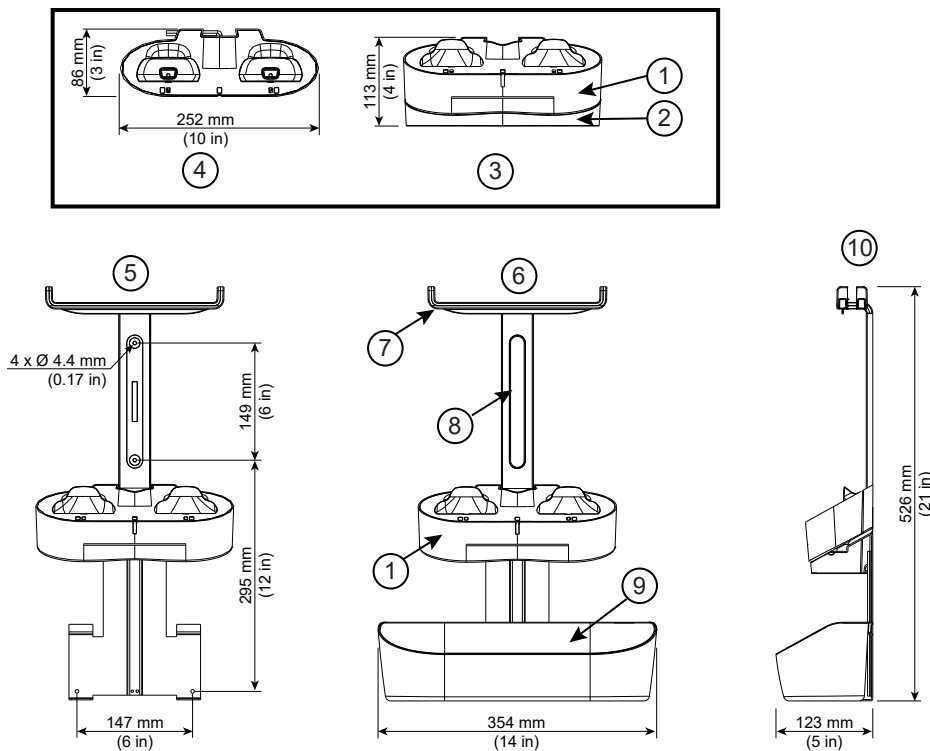

Item	Description	Item	Description
1	Charging cradle	6	Front View - wall mounting configuration with storage basket
2	Desktop base adapter	7	Storage hooks
3	Front View - desktop configuration	8	Wall mount support
4	Top View - desktop configuration	9	Storage basket
5	Front View - wall mounting configuration	10	Side View - wall mounting configuration with storage basket

6 Digital Service and Connectivity Requirements

6.1 InSite RSvP (Remote Service Platform) Requirements



(Applies to all sections within this chapter)

6.1.1 InSite RSvP Connectivity Requirements

Following are the requirements for InSite RSvP connectivity:

1. The customer shall provide a physical connection or a route to an existing enterprise LAN.
2. The cable must be Cat 5 or better.
3. The customer shall provide outbound internet access for the device using HTTPS protocol over port 443.
4. The customer's network administrators shall provide DNS IP Address or Proxy IP address and authentication information (if applicable for the proxy server).
5. The customer's network administrators shall whitelist the following URLs:
 - Enterprise production:
 - <https://insite.gehealthcare.com:443>
 - <https://as1-insite.gehealthcare.com>
 - <https://as2-insite.gehealthcare.com>
 - Flexera URL: <https://gehealthcare-ns.flexnetoperations.com>
 - Flexera Software Download URL: <https://download.flexnetoperations.com>
 - For EU regions, whitelist the following:
 - <https://as1-insite-eu.gehealthcare.com>
 - <https://insite-eu.gehealthcare.com>

7 MR System Interconnects

7.1 MR System Interconnects Specifications



(Applies to all subsections within this section)

7.1.1 Component Designator Definitions

GE HealthCare uses Component Designators to identify system components. All subsystem cabinets and other components are referred to by their component designators in the Interconnect Data diagrams and tables.

Table 7-1 MR System Component Designators

Component Designator	Description
DS, DS1	Door Switch
E01, E02, and so on	Emergency-Off (E-Off) Buttons
ICC	Integrated Cooling Cabinet
MAG	Magnet and Enclosure (all magnet enclosure components in Magnet Room)
MDP	Main Disconnect Panel
MON	Magnet Monitor
MRU	Magnet Rundown Unit
OW	Operator Workstation
PA1	Pneumatic Patient Alert Control Box
PDU	Power Distribution Unit (PDU) is a module in the PGR cabinet
PED	Magnet Rear Pedestal
PEN	Penetration Cabinet
PGR	Power Gradient RF Cabinet
PT	Patient Transport Table
SPW	Secondary Pen Wall

Table 7-2 MR System Options Component Designators

Component Designator	Description
MNS	MNS Cabinet
MRE	Magnetic Resonance Elastography
OXY	Oxygen Monitor
PAT CS	Physiological Acquisition Transceiver Charging Station

7.1.2 Available Cable Lengths

Three configurations of cable lengths are available for order, and multiple selections for gradient cables. Gradient cables are no longer field terminated, and it is critical to select an appropriate length to avoid too much excess. To determine required cable lengths, find the total distance of the cable path between the specified equipment by measuring the following:

- distance from the top of the specified equipment up to the cable tray
- horizontal distance across the cable tray
- distance from the cable tray to the top of the other specified equipment

Compare the total distance to the lengths specified in [Table 7-4 Available Cable and Hose Lengths on page 117](#). If your total distance is more than the "Short" distance, but less than the "Long" distance, the Site Option is "Long." If both of your Site Options are "Long", or if either option is more than the "Long" distance, you must reconfigure the layout of the room.

Table 7-3 Order Configuration Options

Configuration	Equipment Room — Site Option	Magnet Room — Site Option
A	Short	Short
B	Long	Short
C	Short	Long

Table 7-4 Available Cable and Hose Lengths

Length Identifier (shown in the figure below)	Point A	Point B	Site Option: Short	Site Option: Long
			mm (in.)	
Equipment Room				
L1	ICC, Top panel (CRY)	MON	8200 (322.8)	
L2	MON	SPW, top edge	16800 (661.4)	
L3 ¹	MDP	Top of ICC (CRY power)	7500 (295.3)	
Not shown	ICC, top panel	Customer-Supplier Network	28000 (1102.4)	
L4	PEN, top panel	OW, rear panel	29000 (1141.7)	
L5	PEN, top panel	MON	16400 (645.7)	
Not shown	PEN, top edge	RF common ground stud	2673 (105.2)	
Not shown	SPW, top edge	RF common ground stud	1276 (50.2)	
L6	PGR, top panel	OW, rear panel	28800 (1133.9)	
L7	PGR, top panel	ICC, top panel	8800 (346.5)	
L8	PGR, top panel	Magnet Room door switch	24400 (960.6)	
L9	PGR, top panel	PEN, top panel	9400 (370.1)	16400 (645.7)

Table 7-4 Available Cable and Hose Lengths (Table continued)

Length Identifier (shown in the figure below)	Point A	Point B	Site Option: Short	Site Option: Long
			mm (in.)	
L10 (This length is for all cables and hoses except for gradient cables. For gradient cable available lengths, see Table 7-5 Available Cable Lengths for Gradient Cable (Equipment Room, SPW I/O Panel to Top of PGR) on page 119.)	PGR, top panel	SPW, top edge	9200 (362.2)	16200 (637.8)
Not shown	PGR, top panel	RF common ground stud	9200 (362.2)	16200 (637.8)
L11	ICC, Top panel (CRY)	MAG (cold head)	15000 (590.5)	
Not shown	SPW, top edge	E-Off switch, Control Room or Equipment Room	9800 (385.8)	16800 (661.4)
L12	SPW, top edge	OW, rear panel	29200 (1149.6)	
L17	ICC, top panel	SPW, Top edge	9000 (354.3)	11500 (452.7)
L18	ICC, top panel	MAG (Gradient cooling manifold)	28000 (1102.3)	
L19	ICC, top panel	PEN, top panel	9000 (354.3)	11500 (452.7)
Magnet Room				
L13	PEN (magnet room side) top edge	MAG, gradient cable clamp block	9000 (354.3)	13000 (511.8)
L14 (This length is for all cables and hoses except for gradient cables. For gradient cable available lengths, refer to Table 7-6 Available Cable Lengths for Gradient Cable (Magnet Room, Magnet Gradient Cable Strain Relief Plate to Top of SPW I/O Panel) on page 119.)	SPW, top edge	MAG, gradient cable clamp block	7700 (303.2)	11700 (460.6)
Not shown	SPW, top edge	E-Off switch, Magnet Room	29800 (1173.2)	
L15	OW, rear panel	MAG, gradient cable clamp block	38000 (1496.1)	42000 (1653.5)
L16	MRU	MAG, gradient cable clamp block	24700 (972.4)	

Table 7-4 Available Cable and Hose Lengths (Table continued)


Length Identifier (shown in the figure below)	Point A	Point B	Site Option: Short	Site Option: Long
			mm (in.)	
Not shown	RF common ground stud	MAG, gradient cable clamp block	10000 (393.7)	14000 (551.2)
 NOTE ¹ A customer-supplied substitute can be used if the supplied run is shorter than required.				

Table 7-5 Available Cable Lengths for Gradient Cable (Equipment Room, SPW I/O Panel to Top of PGR)

Length Identifier (shown in the figure below)	Available lengths mm (in.)
L10	3700 (145.7)
	5700 (224.4)
	7700 (303.1)
	9700 (381.9)
	11700 (460.6)
	13700 (539.4)
	15700 (618.1)


Table 7-6 Available Cable Lengths for Gradient Cable (Magnet Room, Magnet Gradient Cable Strain Relief Plate to Top of SPW I/O Panel)

Length Identifier (shown in the figure below)	Available lengths mm (in.)
L14	4000 (157.5)
	6000 (236.2)
	8000 (315.0)
	10000 (393.7)
	12000 (472.4)

**NOTE**

Figure 7-1 Available Cable/Hose Lengths on page 120 does not show cable trays. This figure is for reference only, to show available cable lengths and not to be used to design cable routing.

Figure 7-1 Available Cable/Hose Lengths

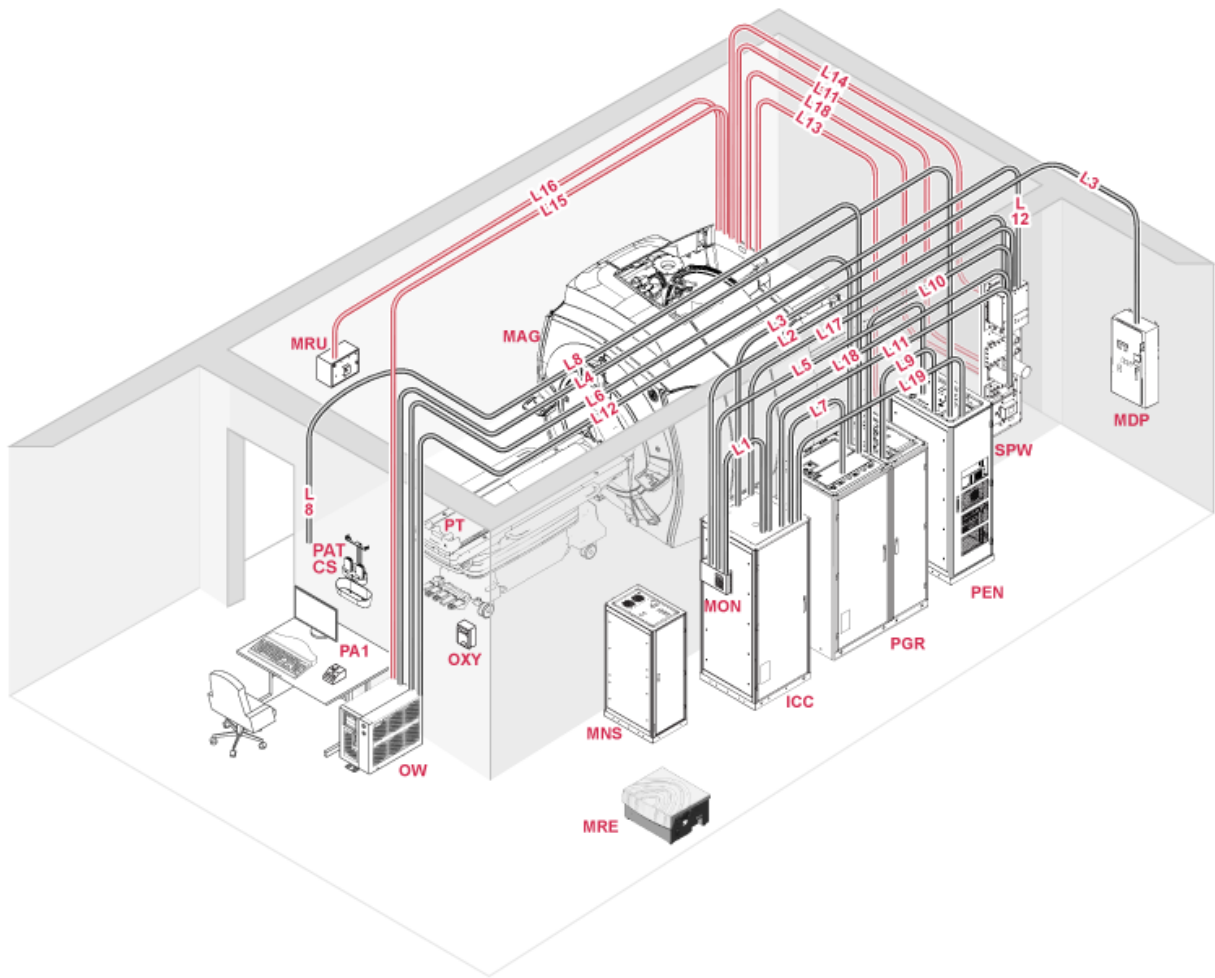

NOTE

MNS, MRE, PAT Charging Station and OXY shown above are optional components of the system.

7.1.3 Multi-Nuclear Spectroscopy (MNS) Option

Table 7-7 MNS Option Available Cable Lengths

Point A	Point B	Site Option A: Short ER, Short SR	Site Option B: Long ER, Short SR	Site Option C: Short ER, Long SR
		mm (in.)		
PGR	MNS Cabinet, top panel	9850 (387.8)		
MNS Cabinet, top panel	SPW	7850 (309.1)	14000 (551.2)	7850 (309.1)
PGR	SPW	8200 (322.8)	14350 (565.0)	8200 (322.8)



NOTE

All cable lengths include a 600 mm (24 in.) service loop.

7.1.4 BrainWave Option

Table 7-8 BrainWave Option Usable Cable Lengths

Cable	From	To	Cable Length mm (in.)
9-pin D-sub to JF-50	SPW, top edge	Top of GOC	34187 (1346)
9-pin D-sub to BNC	PEN, top edge	Top of GOC	27432 (1080)
9-pin D-sub to Coax	SPW, top edge	MAG, gradient cable clamp block	13867 (546)

7.1.5 Magnetic Resonance Elastography (MRE) Option


Table 7-9 MRE Option Available Cable Lengths

Cable	Point A	Point B	mm (in.)
25 mm (1 in.) Tubing	Resoundant Active Driver	Magnet (Isocenter)	Nominal 7315 (288) Maximum 10058 (396)
BNC	Resoundant Active Driver	Penetration Panel	15240 (600)
Ethernet	Resoundant Active Driver	Ethernet Hub in PGR	15240 (600)
Power	Resoundant Active Driver	Customer-Supplied Outlet	60 Hz: 6096 (240) 50 Hz: 7620 (300)

7.1.6 Storage Requirements for Excess Gradient Cable

1. Excess gradient cable must be stored on a wall hook inside the penetration panel closet in the Magnet Room and on a hook anchored into the Equipment Room wall. GE HealthCare supplies one hook per cable.
2. Supports and anchorage must be able to hold up to 22.7 kg (50 lb.) of weight per cable.
3. Excess gradient cable can be stored in a single loop (no more than one loop per cable) with a minimum bend radius of 330 mm (13 in.).
4. Supports for the 1 AWG (X and Y axis) cables shall be anchored no more than 2591 mm (102 in.) above the floor and no more than 2133 mm (84 in.) above the floor for the 2/0 AWG (Z axis) cable.
5. The Z-gradient cable is heavier than the X or Y gradient cable and must be positioned lower, as shown in [Figure 7-4 Gradient Cable Storage Horizontal Layout \(Example\) on page 124](#).


Figure 7-2 Gradient Cable Mounting Hook Mounting Detail

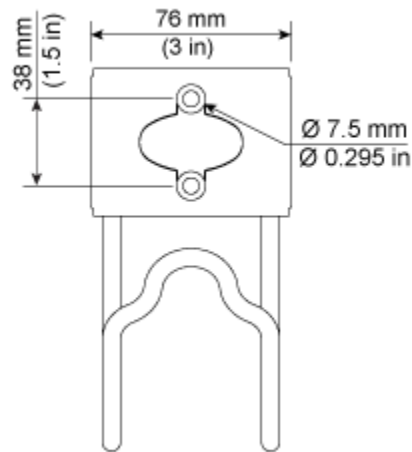

NOTE


The examples below are representative of possible cable layout or configuration. Layout may vary depending on size, shape and configuration of the penetration panel closet. Configure cables to best suit the install configuration without violating the requirements above.

Figure 7-3 Gradient Cable Storage Vertical Layout (Example)

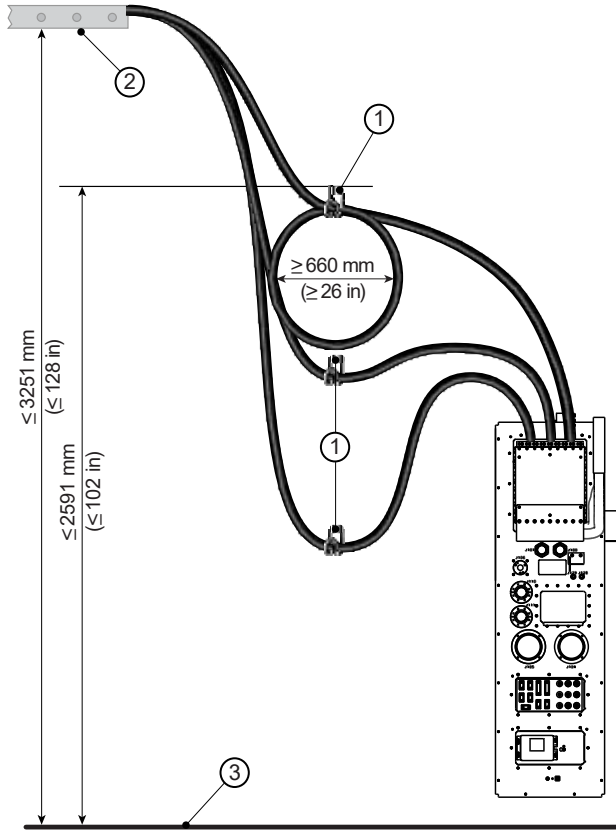

Item	Description	Item	Description
1	GE HealthCare-supplied gradient cable mounting hooks	3	Finished floor
2	Cable tray		

Figure 7-4 Gradient Cable Storage Horizontal Layout (Example)

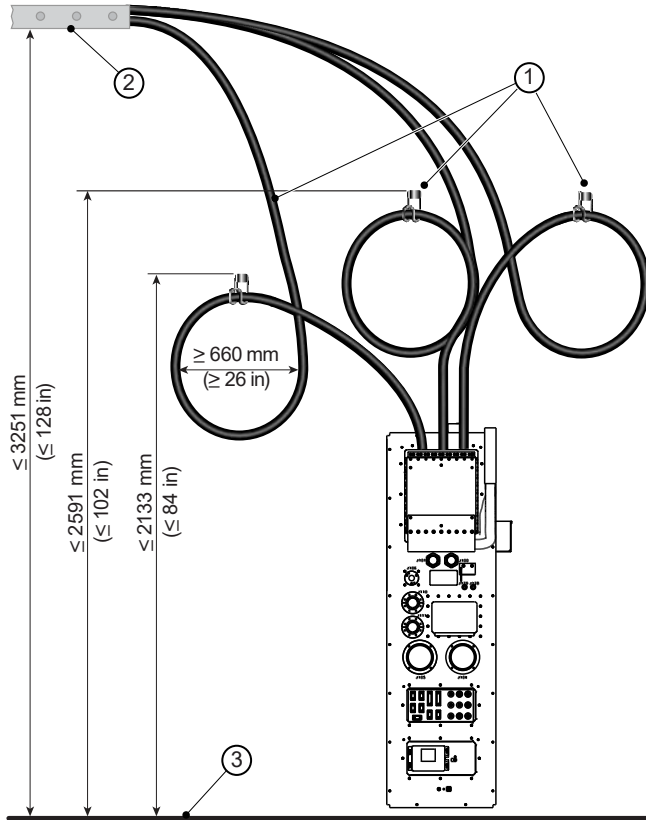

Item	Description	Item	Description
1	GE HealthCare-supplied gradient cable mounting hooks	3	Finished floor
2	Cable tray		

7.2 MR System Interconnects Routing Requirements



(Applies to all subsections within this section)

7.2.1 General Requirements

1. The customer is responsible for the purchase and installation of all cable support mechanisms.
2. Any type of nonferrous cable support can be used, such as a commercially available ladder or wire rack style cable trays, if the cable trays meet all MR System requirements.
3. The distance between cable supports must be less than 305 mm (12 in.). For example, the distance between rungs on a ladder tray, or the distance from the end of a cable tray to a final nonferrous cable support must be 305 mm (12 in.) or less.
4. The cable supports must have the minimum cable bend radius per MR System Cable Specifications. For example, the vertical and horizontal bends of the gradient cables must have a bend radius of 330 mm (13 in.).
5. Cable supports can be stacked or side-to-side.
6. If trays are stacked, the air, water, and cryogen lines must be run in the lower support (see [Figure 7-6 Cable Groupings on page 127](#)).
7. Each cable tray must support a weight of at least 74.8 kg/m (50 lb./ft.).



NOTE

If stacked, each cable tray must support the weight of both cable trays: that is 149.6 kg/m (100 lb./ft.).

8. Cables must be accessible on at least one side of the cable support and require a minimum of 254 mm (10 in.) from the top of the tray. (For example, see [Figure 7-7 Cable Tray Requirements Example \(Side-By-Side\) on page 128](#) or [Figure 7-8 Cable Tray Requirements Example \(Stacked\) on page 129](#).) If this is not possible due to obstructions, see the exceptions below.
9. All individual cable supports (for example, rungs) require lateral support to maintain the positions specified in the illustrations both during installation and after installation of the interconnects (for example, cables).


Exceptions for Obstructions:

1. The top of the cable tray must not touch an obstruction. A minimum of 254 mm (10 in.) of clearance is required on either side of the obstruction. See [Figure 7-5 Obstruction Example on page 126](#).
2. A minimum of 178 mm (7 in.) of clearance is required from the top of the tray rung to the lowest point of any obstruction.

**NOTE**

The illustration below shows how to route cable trays around HVAC ducts, light fixtures, medical gases, structural beams, and other obstructions. If local code permits, a tray with a bend can route the cables underneath an obstruction. At the rear of the magnet, the cable tray must be installed at the minimum height. The part of the cable tray that is under the obstruction can be installed at a lower height. The tray bend must provide the minimum cable bend radius.

Figure 7-5 Obstruction Example

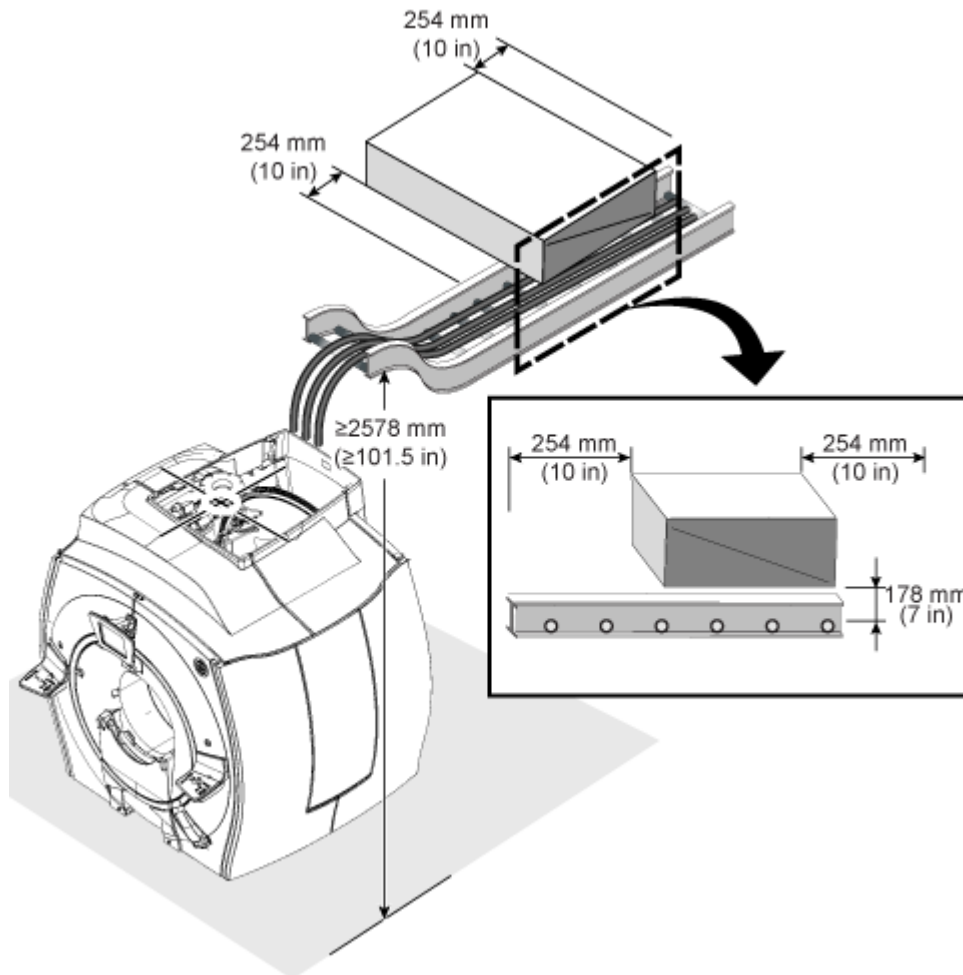

7.2.2 Magnet Room Requirements

1. Two cable trays must be used, each at least 450 mm or 18 in. wide.
2. Installation and routing of cable trays must be coordinated with the RF shield vendor.
3. Side-to-side trays in the Magnet Room must not touch to prevent RF broadband noise caused by metal-to-metal sidewall contact.
4. Ceiling grid work, medical gas lines, lighting fixtures, and so on, must not touch MR System cabling or cable supports.
5. Excess cable length in the Magnet Room must be stored in either:
 - 5.1. Penetration Panel closet. If utilizing the Penetration Panel Closet for cable storage, the supports and anchors must be able to hold up to 22.7 kg (50 lb.).

- 5.2. Magnet Room cable trays (excess cable must be at least 915 mm (36 in.) from the magnet end of the tray)


7.2.3 Recommended Cable Groupings

We recommend the grouping order shown below. The cable groups must not touch. Local code may require dividers, channeling, bundling, and so on.

- Group a — Gradient and RF common ground cables

In the Magnet Room, the gradient cable group must be laid in a single layer. The center of cable group “a” is 89 mm (3.5 in.) from the inside edge of the tray, in line with the magnet center. See [Figure 7-6 Cable Groupings on page 127](#).
- Group b — Coax, RF, and AC power cables with jacket rating of 600V and above
- Group c — Data and fiber optic and 300V coax or RF clock cables
- Group d — 300V signal, 300V power, and 300V power and signal cables
- Group e — Air, water, and cryogen lines

Figure 7-6 Cable Groupings

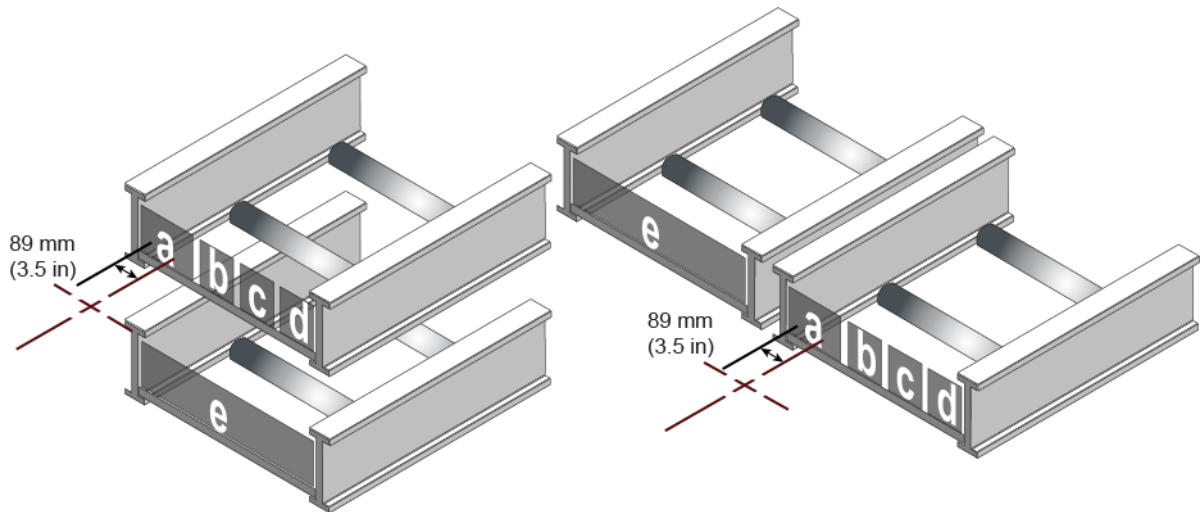

7.2.4 Cable Tray Requirements and Examples

1. The gradient cable support must end at the back of the magnet 1245 ± 12 mm (49 ± 0.5 in.) from the geometric isocenter, and must be parallel to the rear face of the magnet. See [Figure 7-7 Cable Tray Requirements Example \(Side-By-Side\) on page 128](#), [Figure 7-8 Cable Tray Requirements Example \(Stacked\) on page 129](#), and [Figure 7-9 Cable Tray Requirements Example \(90° Magnet Interface\) on page 130](#).
2. Supports for all other cables and hoses must end at the back of the magnet 864 ± 12 mm (34 ± 0.5 in.) from the geometric isocenter, and must be parallel to the rear face of the magnet.
3. Cable supports must have a minimum height of 2578 mm (101.5 in.) at the back of the magnet measured from the floor. See [Figure 7-7 Cable Tray Requirements Example \(Side-By-Side\) on page 128](#).



NOTE

Supports may be lower at other points along the route to clear obstructions as long as all other requirements are met.

- 4. The top of the cable tray must be less than 3251 mm (128 in.) above the finished floor. See [Figure 7-7 Cable Tray Requirements Example \(Side-By-Side\)](#) on page 128 .

Figure 7-7 Cable Tray Requirements Example (Side-By-Side)

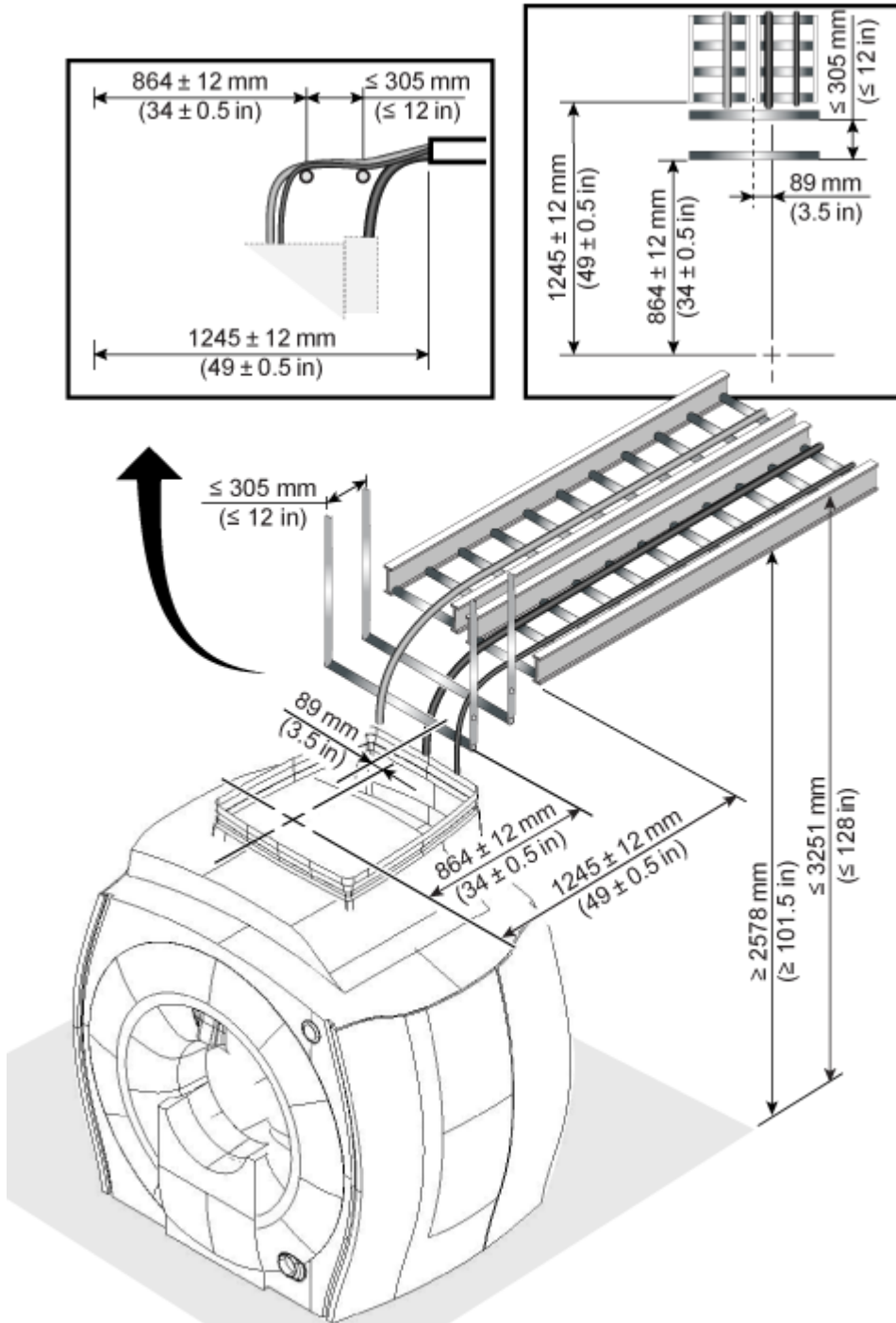


Figure 7-8 Cable Tray Requirements Example (Stacked)


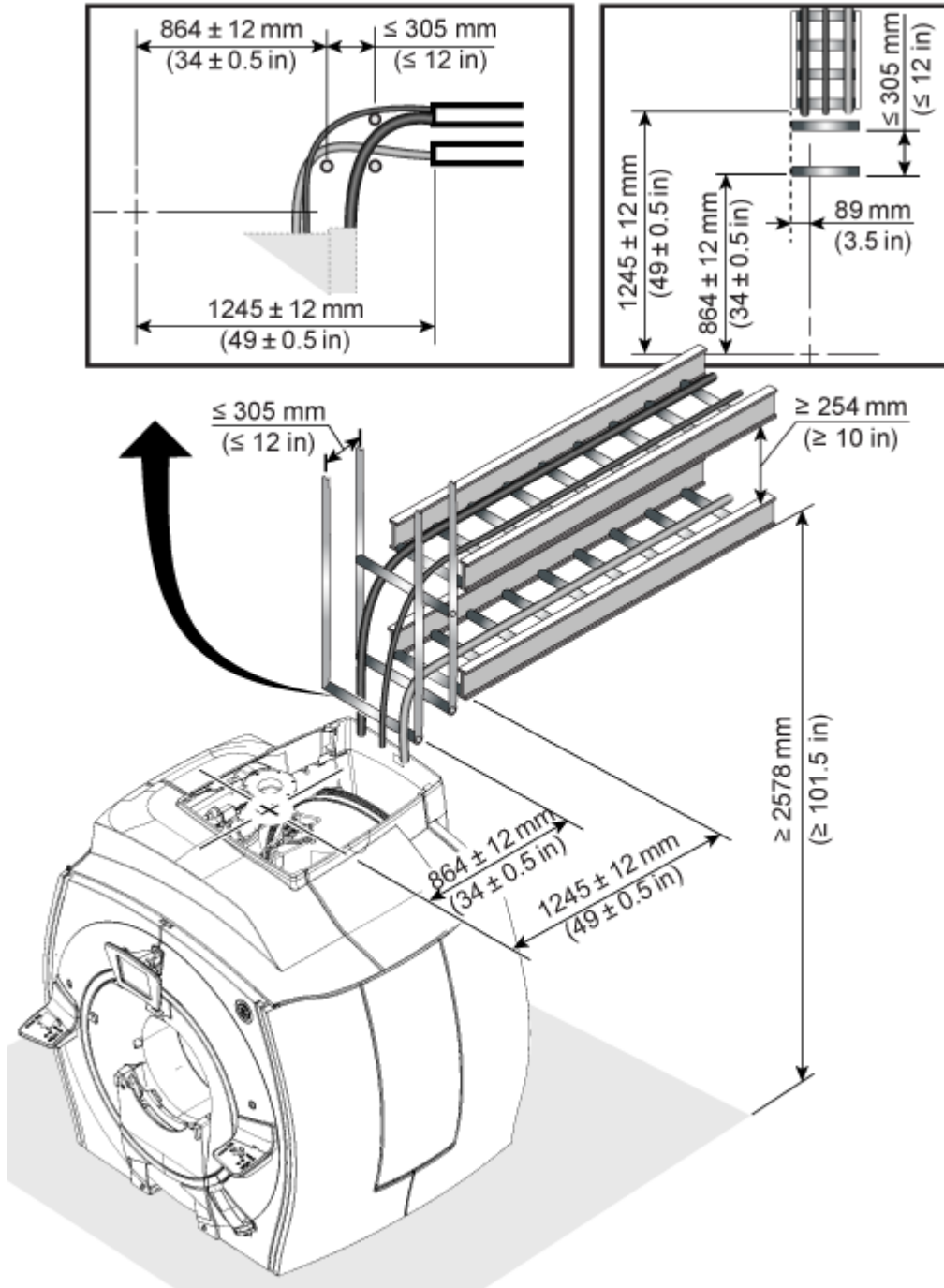


Figure 7-9 Cable Tray Requirements Example (90° Magnet Interface)

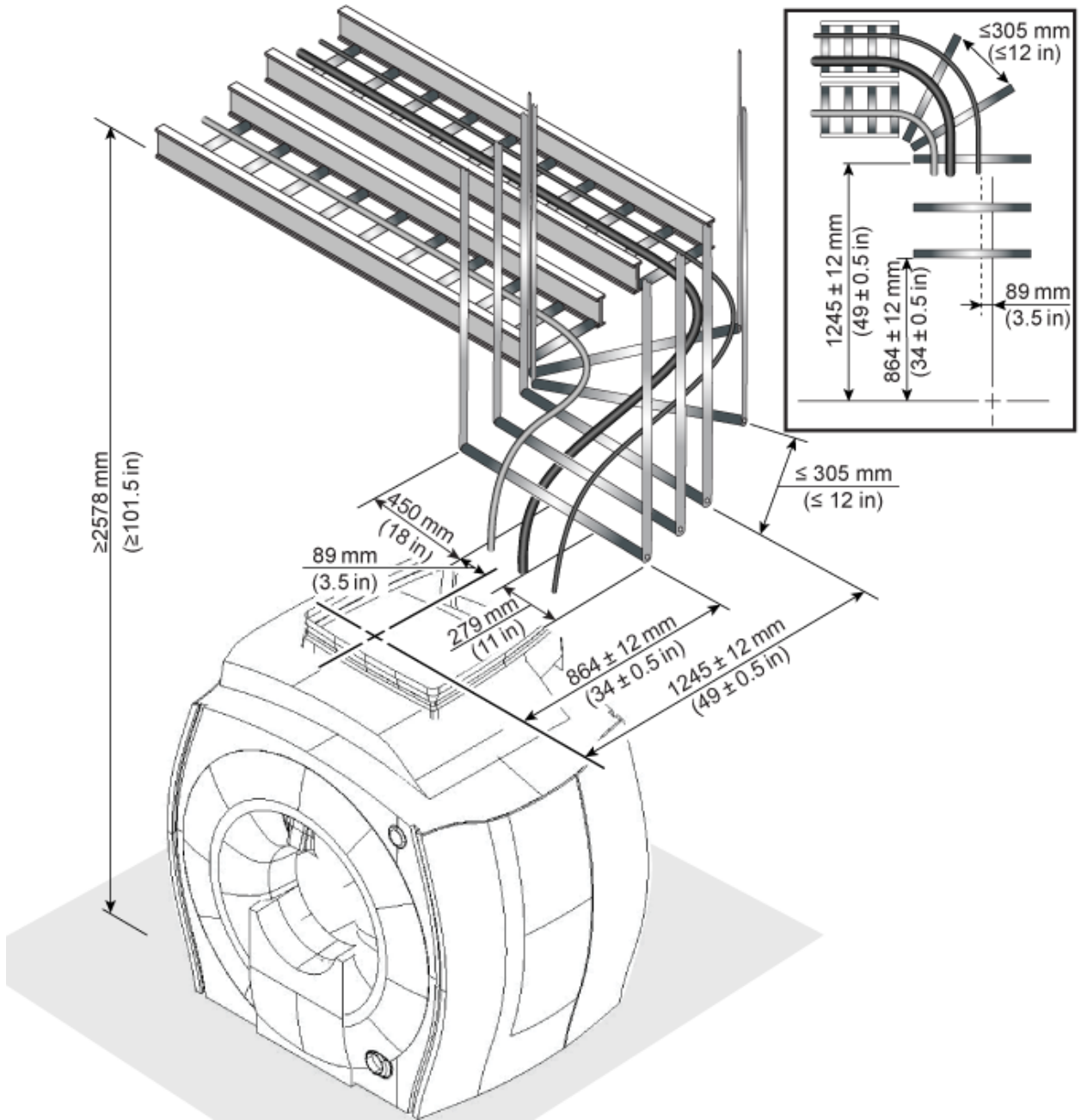

7.2.5 Penetration Panel Closet Requirements

1. The cable support must be positioned to prevent cables from resting on the top of the Penetration Panel Closet. The support must end 560 mm (22 in.) in front of the Penetration Panel.
2. The gradient cable support height and distance from the SPW must support a minimum cable bend radius of 330 mm (13 in.) to accommodate gradient cables.

7.2.6 Equipment Room Requirements

1. All equipment interconnects must route overhead, except for helium hoses.
2. Cables or hoses must drop through the bottom or off the end of the cable support directly to the top of the cabinets (see individual components in Chapter 4 for height requirements).
3. Cable supports do not have minimum height restrictions except above Equipment Room components (to maintain the 330 mm (13 in.) bend radius drop to the PGR cabinet).
4. The end of the gradient cable support must be parallel to the Secondary Penetration Wall (SPW).
5. The gradient cable support height and distance from the SPW must support a minimum cable bend radius of 330 mm (13 in.) to accommodate the gradient cable bend radius.
6. Excess cable length must be stored in the Equipment Room.
7. For multiple MR System installations, cables from different MR Systems must not share the same cable support.
8. The table below lists the minimum width for cable trays between Equipment Room (and Operator Workspace) components.

Table 7-10 Minimum Cable Tray Width

From	To	Minimum Cable Tray Width mm (in.)	For
PEN	PGR	150 mm (6 in.)	Electrical
PEN	ICC	150 mm (6 in.)	Electrical
		300 mm (12 in.)	Air
PEN	OW	150 mm (6 in.)	Electrical
SPW	PGR	300 mm (12 in.)	Electrical
SPW	ICC	150 mm (6 in.)	Electrical
		300 mm (12 in.)	Air and Water
SPW	OW	150 mm (6 in.)	Electrical
ICC	PGR	150 mm (6 in.)	Electrical
		76 mm (3 in.)	Water
PGR	OW	150 mm (6 in.)	Electrical

7.3 Facility-Supplied System Interconnects Specifications



On installation sites in China, make sure that the power cables and ground cables provided by customers have China Compulsory Certification (CCC). This information is supplied to the customer in China Power Cable Requirements, 5159493. (Go to the Customer Documentation Portal or contact the PMI.)

The following table lists the required facility-supplied system interconnects. Refer to [Figure 7-10 Facility-Supplied System Interconnects](#) on page 133 for additional information.

Table 7-11 Facility-Supplied System Interconnects

Group	Between Units		Comments	Requirements
	From	To		
C01	Facility Power	MDP	Facility Power and Ground	MR Suite Electrical Requirements on page 47
C02	MDP	ICC		
C03	MDP	PGR		
	Facility Cooling Water	ICC	Cooling Water Supply	MR System Facility Water Requirements on page 41
	Facility Cooling Water	ICC	Cooling Water Return	
C04	Facility Network	MON	Facility must provide separate network access for the Magnet Monitor (MON), the Global Operator Cabinet (GOC) and the MDP. The MON connection must be available at all times.	Magnet Monitor (MON) Requirements and Specifications on page 103
	Facility Network	GOC		InSite RSvP (Remote Service Platform) Requirements on page 115
	Facility Network	Receptacle near MDP		4.2 Main Disconnect Panel (MDP) Requirements and Specifications on page 92
C05	MDP	E-Off Switch	Facility must supply cable from MDP to the E-Off Switch in the Equipment Room.	MR Suite Electrical Requirements on page 47
	Facility Power	Outlet near MON	Facility outlet for MON power	Magnet Monitor (MON) Requirements and Specifications on page 103
C06	GE HealthCare supplied Cable from SPW	E-Off Switch in Control Room or Equipment Room	Facility must supply additional wiring between the GE HealthCare supplied cable and the E-Off switch if the length needed is greater than the available length listed in MR System Interconnects Specifications on page 116 .	MR Suite Electrical Requirements on page 47


Table 7-11 Facility-Supplied System Interconnects (Table continued)

Group	Between Units		Comments	Requirements
	From	To		
	Facility Power	MRU	Facility power to MRU	Magnet Room Equipment Specifications on page 83
	Facility Power	Service Receptacle	Receptacle to be installed in Magnet Room, using appropriate filter.	MR Suite Electrical Requirements on page 47
	Facility Power	Pneumatic Patient Alert	Facility Power to Pneumatic Patient Alert	MR Suite Electrical Requirements on page 47

Table 7-12 Optional Facility-Supplied System Interconnects

Group	Between Units	Comments	Requirements
C01	Facility Power	Outlet near MRE	Facility outlet for MRE power (optional, not shown)
C01	Facility Power	Control room	Hardwired connection for OXY power (optional, not shown)
C01	Facility Power	Outlet near PAT Charging Station	Facility outlet for PAT Charging Station (optional, not shown)

Figure 7-10 Facility-Supplied System Interconnects

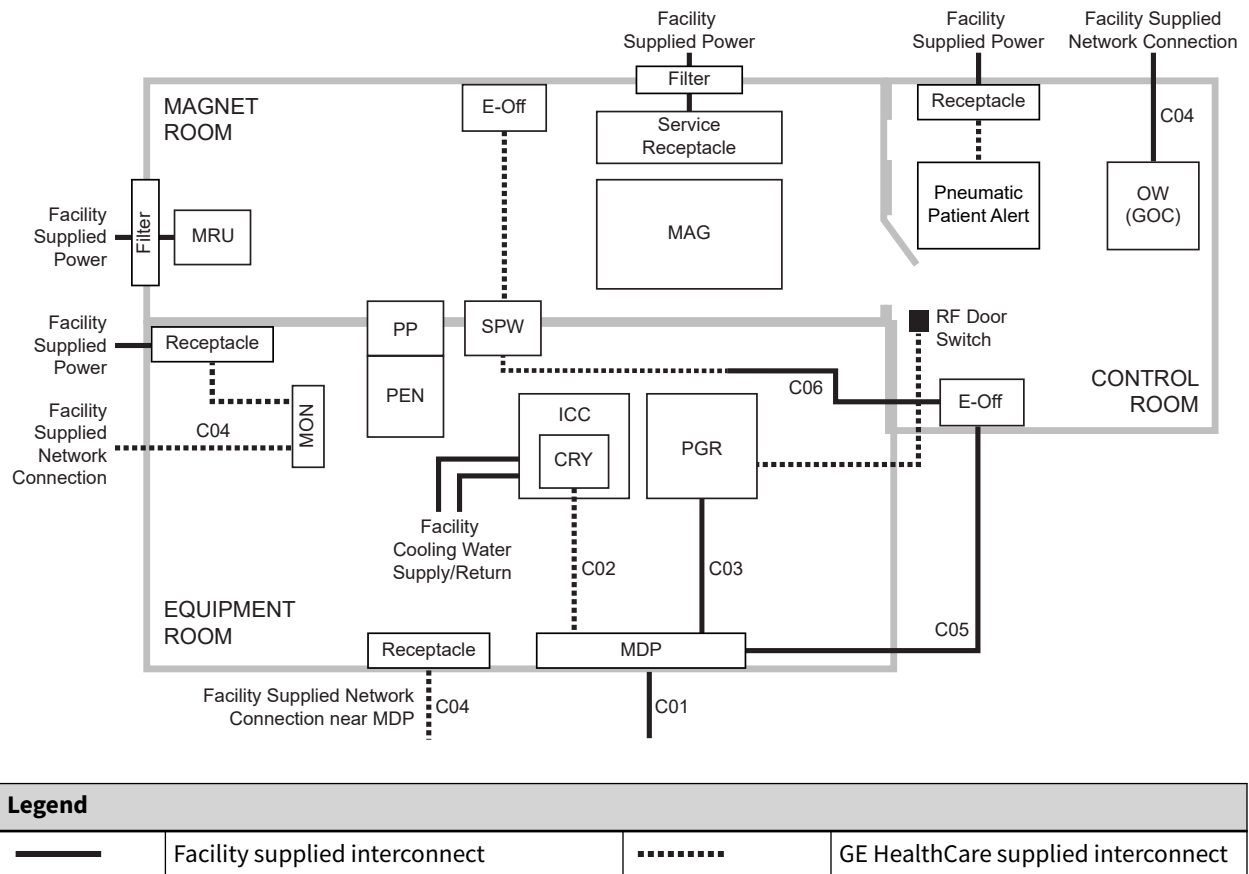

**NOTE**

- GE HealthCare recommends installing the RF Door switch on the outside wall of the Magnet Room.
- The illustration is not to scale and component positioning/interconnect runs are typical.
- Only GE HealthCare equipment interconnects are shown. Additional facility interconnects are required for non-GE HealthCare equipment (for example, Magnet Room DC Lighting).
- The E-Off button placement and cable routing shown indicates one possible configuration. Final E-Off button placement and cable routing is the responsibility of the customer. Refer to [2.10.1 General Electrical Requirements on page 47](#) for requirements for E-Off button placement.

The RF Shielded Room Vendor is responsible for installing the RF door switches. Refer to *RF Shielded Room Requirements*, 5850260-1EN.

8 Appendix

8.1 Glossary



Cryogen

A substance for producing low temperatures. Liquid helium is the cryogen used to cool the magnet to approximately 4 K (-269°C or -452°F).

Dewar

A container with an evacuated space between two highly reflective walls used to keep low temperature substances at near-constant temperatures. Liquid helium is usually stored and shipped in dewars.

Ferrous Material

Any substance containing iron which is strongly attracted by a magnetic field.

Gauss (G)

A unit of magnetic flux density. The earth's magnetic field strength is approximately one half Gauss to one Gauss depending on location. The internationally accepted unit is the tesla (1 tesla = 10000 G and 1 millitesla = 10 G).

Homogeneity

Uniformity. The homogeneity of the static magnetic field is an important quality of the magnet.

Isocenter

Center of the imaging volume ideally located at the magnet center.

Isogauss Line

A line on a field plot connecting identical magnetic field strength points.

Magnetic Field

A condition in a region of space established by the presence of a magnet and characterized by the presence of a detectable magnetic force at every point in the region. A magnetic field exists in the space around a magnet (or current carrying conductor) and can produce a magnetizing force on a body within it.

Magnetic Resonance (MR)

The absorption or emission of electromagnetic energy by nuclei in a static magnetic field, after excitation by a suitable radio frequency field.

Magnetic Shielding

Using material (for example, steel) to redistribute a magnetic field, usually to reduce fringe fields.

Quench

Condition when a superconducting magnet becomes resistive thus rapidly boiling off liquid helium. The magnetic field reduces rapidly after a quench.

Radio Frequency (RF)

Frequency intermediate between audio frequency and infrared frequencies. Used in magnetic resonance systems to excite nuclei.

Radio Frequency Shielding

Using material (for example, copper, aluminum, or steel) to reduce interference from external radio frequencies. A radio frequency shielded room usually encloses the entire Magnet Room.

Resonance

A large amplitude vibration caused by a relative small periodic stimulus of the same or nearly the same period as the natural vibration period of the system. In magnetic resonance imaging, the radio frequency pulses are the periodic stimuli which are at the same vibration period as the hydrogen nuclei being imaged.

Superconducting Magnet

A magnet whose magnetic field originates from current flowing through a superconductor. Such a magnet is enclosed in a cryostat.

Superconductor

A substance whose electrical resistance essentially disappears at temperatures near zero Kelvin. A commonly used superconductor in magnetic resonance imaging system magnets is niobium-titanium embedded in a copper matrix.

Tesla

The internationally accepted unit of magnetic flux density. One tesla is equal to 10000 Gauss. One millitesla is equal to 10 Gauss.

Waveguide

A hollow linear structure used in components such as the penetration wall to route cables and hoses, while limiting and controlling electromagnetic waves from entering the Magnet Room.

8.2 MR Site Vibration Test Guidelines



(Applies to all subsections within this section)

8.2.1 Test Measurements

1. Vibration measurements must be in the range of 10^{-6} g. Test equipment must have the required sensitivity to these levels.
2. Instrumentation must have a low tolerance to temperature effects since many times the low frequency thermal drift may influence the measurements.
3. All measured data must be acquired real time. Recording of vibration data will not allow for a correct site survey, specifically when studying transient vibration and when searching for specific vibration sources.
4. All analyses must be narrow-band Fast Fourier Transforms (FFT) over the frequency bands listed in [Table 8-1 Frequency Bands for FFT on page 137](#).
5. Time histories of the vibration must be recorded as acceleration levels vs. time. The resolution of the time history must be adjusted to clearly capture the transient event. The analyzer set-up will be site dependent and, in special cases, vibration response dependent. It is the responsibility of the vibration consultant to study the transient environment, capture data to confirm that transient activity exceeds the trigger level, then expand the time history data to exhibit the structural response.

Table 8-1 Frequency Bands for FFT

Frequency Band	Frequency Resolution
0.2 to 50 Hz	$\Delta f = 0.125$ Hz

8.2.2 Equipment (Spectral Analyzer) Set-Up

1. Frequency average should be a minimum of 20 linear averages (Do not use peak hold or 1/3 octave analysis).
2. Average and store should be a minimum of 20 plots steady state and 20 plots transient to support the consistency of the site vibrations.
3. Hanning windows must be applied to the entire spectra.
4. Spectrum analyzers capable of these measurements are readily available for purchase or rental. Models, such as the HP 3560A, Nicolet Phaszer, B&K Pulse, and HP 35670, are all capable of making the site vibration measurements. Accelerometers must have the capability to measure from 0.2 Hz beyond 50 Hz. Time histories can be recorded using any of the analyzers listed above.



NOTE

The equipment mentioned is for example only. It is the responsibility of the Engineering Test Firm to provide equipment that will allow measurements compliant with this guideline.

8.2.3 Data Collection

8.2.3.1 Ambient Baseline Condition

1. All of the measurements listed above must be made in a “quiet” environment—that is, areas where excessive traffic, subway trains, and so on, do not exist. A vibration measurement must also be made during periods without traffic or during periods of light traffic. Measurements must define the lowest levels of vibration possible at the site.
2. The source of any steady state vibration, whose level exceeds the magnet specifications found in [Magnet Room Structural Requirements on page 64](#), must be identified. A second measurement should be made with all of the identified contributors powered down if possible. In situations where it is not possible to power down equipment, vibration data must be collected to identify the specific source of the vibration concern. The majority of steady state vibration problems can be negated by isolating the vibration source.

8.2.3.2 Normal Condition

1. All of the vibration measurements listed above must be repeated during periods of “normal” environmental conditions, including the Fast Fourier Transforms (FFTs) and time histories. The transient measurements must be provided to define the dynamic disturbances the MR System may be exposed to. Transient analysis is required for a true assessment of the site.
2. Special attention must be paid to the site assessment during the entire analysis. Since transient vibration is not easily addressed once the MR suite is fully constructed, the test consultant must fully understand the needs for this analysis. The source of any transient vibration must be identified and supported with vibration plots. If the source of any transient vibration is not locatable, it is recommended that the customer have an alternate location identified and the vibration studied.
3. Transient vibration can be difficult to assess if the details are not understood. The **0.0005 g, zero-to-peak trigger level** is a starting point to understanding the vibration stability. The transient vibration peak amplitude, structural (time variant) response, decay rate and an estimate of the number of events per unit of time would constitute a complete transient analysis. All transient failures must be supported by time history plots. The plots must clearly show the structural response, the frequency of the signature and the decay rate. From this data, GE HealthCare can help determine compliance with the vibration guidelines.
4. The test consultant must provide site data to show the design recommendations for all sites/building structures meet the magnet specifications found in [Magnet Room Structural Requirements on page 64](#).

8.2.4 Presentation/Interpretation of Results

1. The recommended format for site vibration data collection, presentation, and analysis is demonstrated in the examples in [Figure 8-1 Example of Site Environmental Vibration Levels on page 139](#), [Figure 8-3 Acceleration Time History on page 140](#), and [Figure 8-4 Acceleration Time History \(Zoomed In on Transient Event\) on page 141](#). Presentation of the data in any other format (linear units only) may result in incorrect interpretation and diagnosis of the site. Additional data collection or presentation methods are at the option of the vibration testing service.
2. All plots must be properly annotated with:

- 2.1. Instrumentation setup including number of averages, frequency resolution, and so on
- 2.2. Test location
- 2.3. Test conditions:
 - 2.3.1. Steady state
 - 2.3.2. Transient
 - 2.3.3. Heel drop
 - 2.3.4. Normal environment
 - 2.3.5. Typical traffic
 - 2.3.6. Any other conditions necessary to demonstrate understanding of potential sources of vibration
- 3. The customer's vibration testing service is responsible for interpreting the results and determining if that site meets GE HealthCare specifications.
- 4. If the vibration levels are too high, additional data acquisition may be necessary to:
 - 4.1. Determine the source of the vibration
 - 4.2. Propose a solution to the problem
 - 4.3. Find an alternate site location
- 5. Any questions regarding test equipment requirements, test parameters, or general questions should be discussed with the GE HealthCare Project Manager of Installation (PMI).

Figure 8-1 Example of Site Environmental Vibration Levels

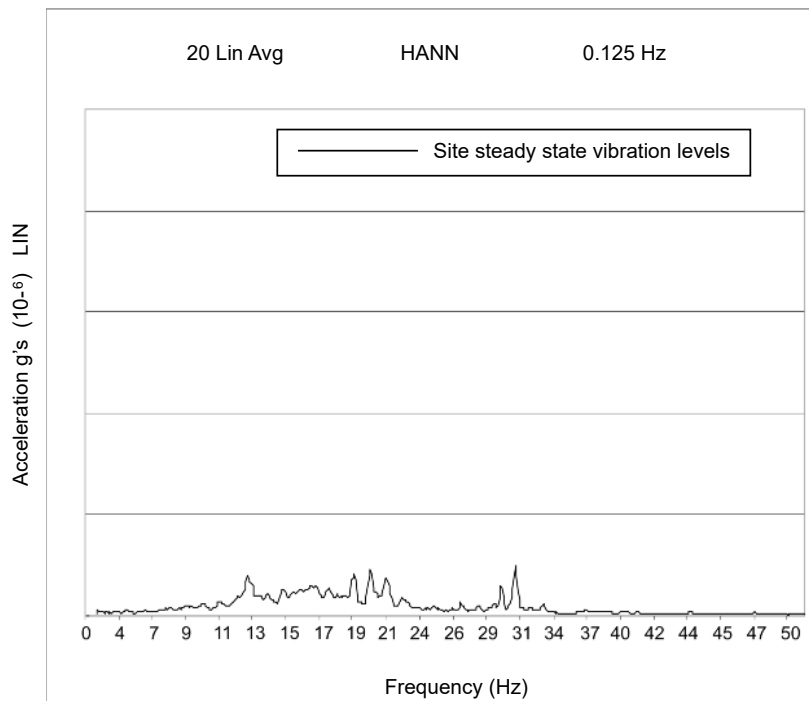


Figure 8-2 Example Site Environmental Vibration

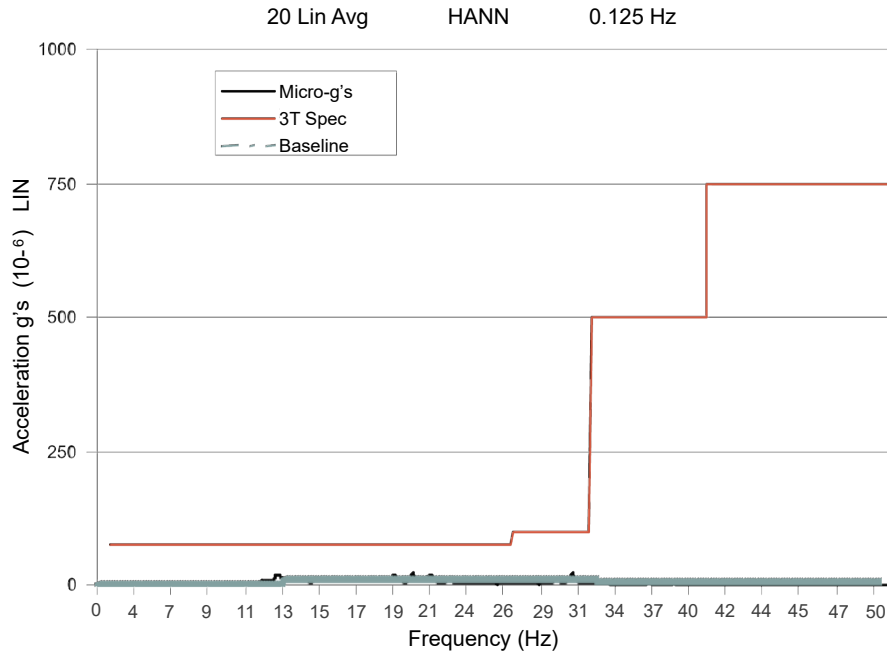


Figure 8-3 Acceleration Time History


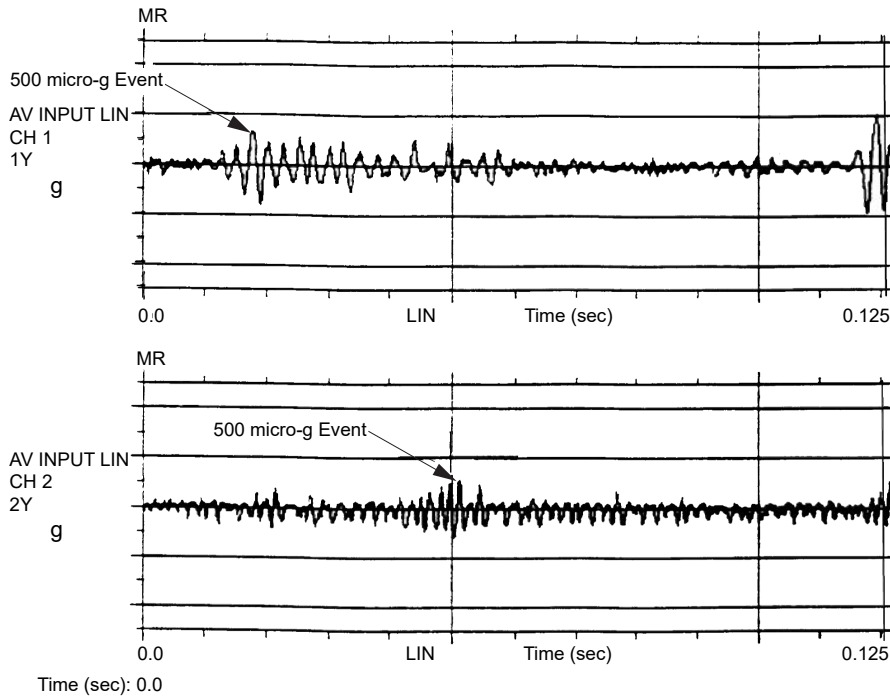


Figure 8-4 Acceleration Time History (Zoomed In on Transient Event)



8.3 Sample Calculation AC Power Equipment Minimum Distance

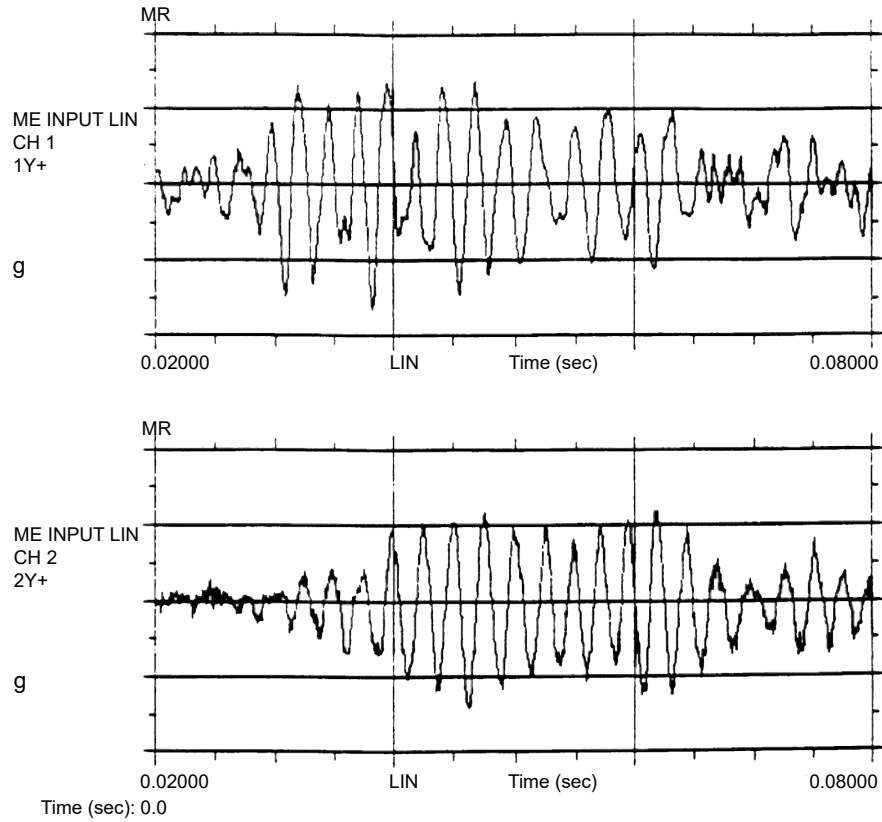

This is a sample calculation to determine minimum distance from a feeder, transformer, or other AC electrical source, using the formula found in 2.6.3 Electrical Current on page 35 to determine minimum distance from a feeder, transformer, or other AC electrical source.

$$I \text{ (amps)} = 8.55x^2 \text{ (meters)} \div S \text{ (meters)}$$


Note that the formula has three variables. If you have two of them, you can calculate the third. In this example, we calculate the minimum distance **x** from the source—in this case, a main electrical feeder carrying 450 amps of current in a 127 mm (5 in.) conduit.

Rearranging:	$x^2 = \frac{I \times S}{8.55}$ $x = \sqrt{\frac{I \times S}{8.55}}$
--------------	--

where:	
x	Minimum distance (in meters) from the feeder lines to isocenter of the magnet
I	Maximum allowable RMS single phase current (in amps) or maximum allowable RMS line current (in amps) in three phase feeder lines
S	Separation (in meters) between single phase conductors or greatest separation between three phase conductors

The separation **S** is the spacing between the conductors, and when all 3 conductors are run in a single conduit, **S** is simply the diameter of the conduit.

$$S = 5 \text{ inches} = 0.127 \text{ meters}$$

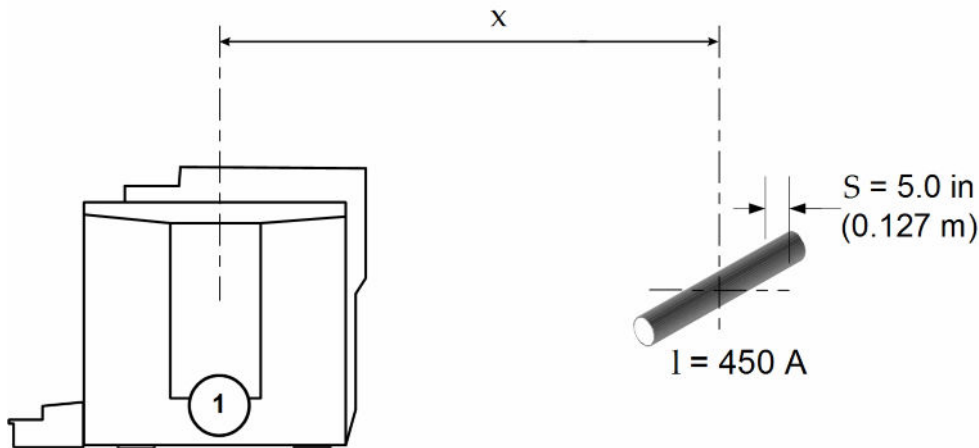

Item	Description
1	Magnet

$$x = \sqrt{\frac{450 \text{ A} \times 0.127 \text{ m}}{8.55}}$$

The conduit should be 2.6 m (8.5 ft.) from the magnet isocenter.

In other situations, the spacing **S** may be the spacing between HV feeders, the distance between transformer lugs, or the spacing between conduits when the phase conductors are run in separate conduits.

What if it is too close?

If this is an existing condition, you should request an *EMI study* to quantify the magnitude and direction of the AC disturbances. The calculation is worst-case and does not take into account the vector direction of the AC interference. The magnet is only sensitive to AC disturbances that are directed horizontally (magnet z-axis). Also, the calculation does not account for any magnetic shielding effect of steel conduit.

8.4 Selecting Anchor Size



The following is an example to illustrate the selection of correct anchors to install a magnet in a building with 13.8 MPa (2000 psi) concrete. For this example the area is not under seismic requirements.

1. Determine the magnet clamping force (for the Magnet: 11100 N + 900 N = 12000 N (2500 lb. + 200 lb. = 2700 lb.)).
2. Refer to the examples of anchor vendor catalogs below to select the anchor diameter and embedment that meets the clamping force (tension) determined in Step 1.

Diameter : ≥ 15.875 mm (0.625 in.) ≤ 31.75 mm (1.25 in.)

For 203 mm (8 in.) embedment select 19 mm (0.75 in.) diameter

For 114.3 mm (4.5 in.) embedment select 25.4 mm (1 in.) diameter

or

Diameter : Min. M16 Max. M32

For 130 mm embedment select M20 diameter

For 114 mm embedment select M24 diameter

3. The vendor instructions and torque to the maximum recommended level for the anchor selected in Step 2 must be provided to the RF shielded room vendor for correct installation of the anchor and equipment.

Table 8-2 Allowable Anchor Loads in Concrete (English Units)

Anchor Diameter mm (in.)	Embedment Depth mm (in.)	13.8 MPa (2000 psi)		20.7 MPa (3000 psi)		27.6 MPa (4000 psi)		41.4 MPa (6000 psi)	
		Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)
15.9 (5/8)	70 (2 3/4)	5.6 (1250)	12.5 (2800)	7.1 (1600)	13.7 (3070)	8.1 (1810)	14.8 (3300)	8.5 (1920)	12.5 (3330)
	102 (4)	8.3 (1870)	14.8 (3330)	10.7 (2400)	14.8 (3330)	13.0 (2930)	14.8 (3330)	14.2 (3200)	12.5 (3330)
	178 (7)	11.2 (2500)	14.8 (3330)	13.4 (3010)	14.8 (3330)	16.2 (3650)	14.8 (3330)	16.2 (3650)	12.5 (3330)
19.1 (3/4)	83 (3 1/4)	6.9 (1550)	12.8 (2880)	8.7 (1950)	14.7 (3310)	10.5 (2350)	16.6 (3730)	11.6 (2610)	21.4 (4800)
	121 (4 3/4)	11.2 (2510)	20.1 (4510)	14.5 (3250)	20.7 (4650)	17.2 (3870)	21.4 (4800)	20.8 (4670)	21.4 (4800)
	203 (8)	13.0 (2930)	21.4 (4800)	17.2 (3870)	21.4 (4800)	20.2 (4530)	21.4 (4800)	22.8 (5120)	21.4 (4800)

Table 8-2 Allowable Anchor Loads in Concrete (English Units) (Table continued)

Anchor Diameter mm (in.)	Embedment Depth mm (in.)	13.8 MPa (2000 psi)		20.7 MPa (3000 psi)		27.6 MPa (4000 psi)		41.4 MPa (6000 psi)	
		Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)
25.4 (1)	114 (4 1/2)	13.9 (3120)	27.0 (6080)	17.2 (3870)	30.1 (6770)	20.5 (4610)	33.2 (7470)	21.4 (4800)	33.2 (7470)
	152 (6)	19.6 (4400)	33.2 (7470)	28.5 (6400)	33.2 (7470)	32.0 (7200)	33.2 (7470)	32.6 (7330)	33.2 (7470)
	229 (9)	24.9 (5600)	33.2 (7470)	35.59 (8000)	33.2 (7470)	41.77 (9390)	33.2 (7470)	41.8 (9390)	33.2 (7470)



NOTE  All bolded values in this table fail to meet the clamping force (tension), and are therefore not acceptable anchors.

Table 8-3 Allowable Anchor Loads in Concrete (Metric Units)

Anchor Diameter	Embedment Depth mm (in.)	13.8 MPa (2000 psi)		20.7 MPa (3000 psi)		27.6 MPa (4000 psi)		41.4 MPa (6000 psi)	
		Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)	Tension kN (lb.)	Shear kN (lb.)
M16	105 (4 1/8)	11.2 (2500)	25.1 (5650)	20.9 (4705)	39.9 (8965)	24.2 (5450)	45.0 (10125)	30.7 (6900)	46.9 (10550)
M20	130 (5 1/8)	25.1 (5650)	52.9 (11900)	30.7 (6910)	58.7 (13195)	36.4 (8175)	64.5 (14490)	44.5 (10005)	64.5 (14490)
M24	155 (6 1/8)	30.0 (6735)	61.2 (13760)	36.9 (8300)	70.5 (15855)	43.9 (9860)	29.8 (17950)	57.7 (12980)	95.6 (21490)


NOTE  All bolded values in this table fail to meet the clamping force (tension), and are therefore not acceptable anchors.

8.5 Sample control schematic for customer-supplied MDP



This section provides an example of a control schematic for the auto-restart and Emergency Power Off (EPO) functions that meets the minimum GE HealthCare PIM requirements. This schematic does not include control, protection, wiring or devices that may be required due to local safety and regulatory requirements. Only the minimum equipment, devices and wiring is shown to meet the performance requirements of GE HealthCare equipment. The final MDP design must be compliant with applicable local codes and regulations.

Figure 8-5 Customer MDP control schematic

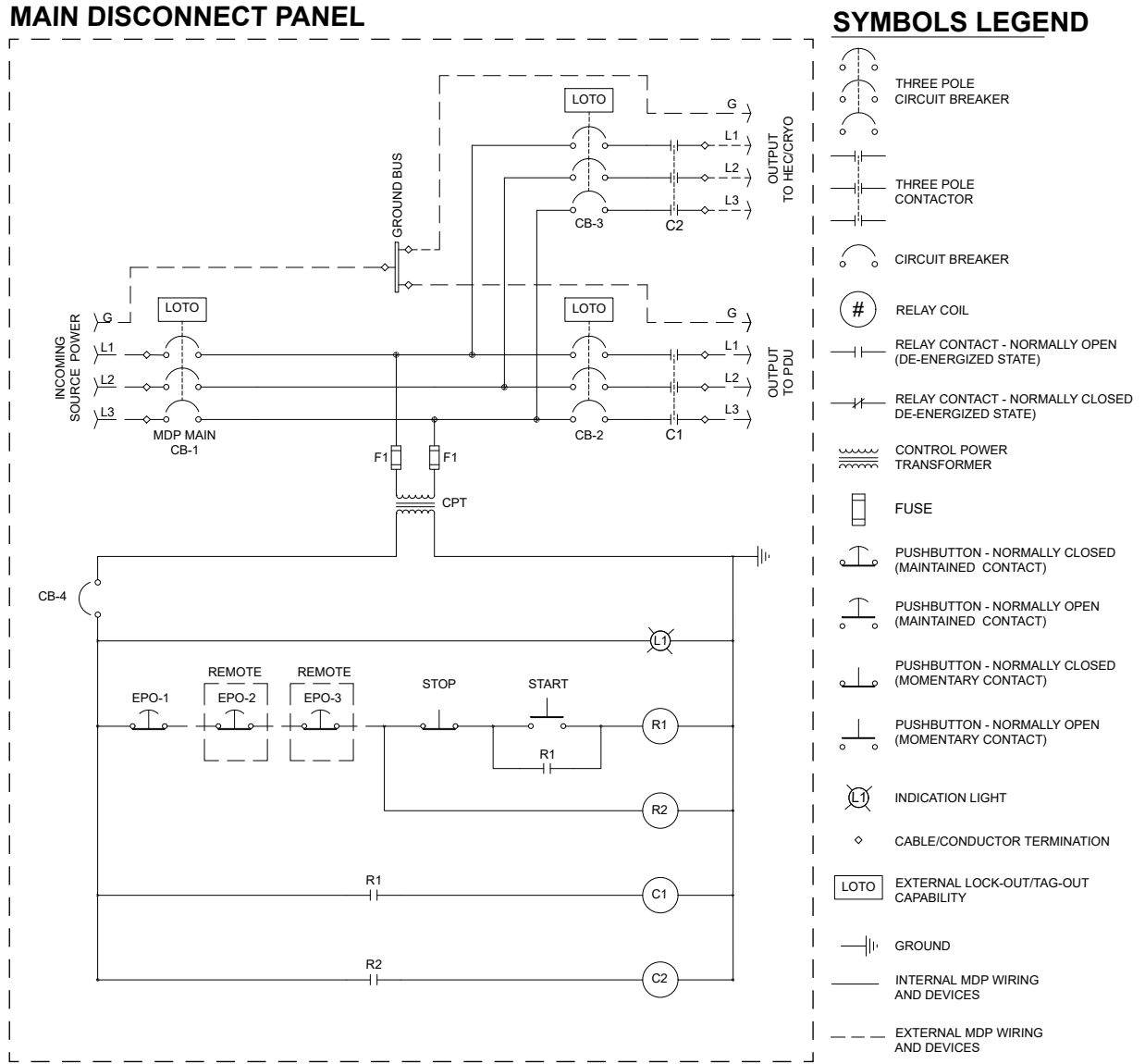

Revision History

Revision 15, October 2024, English Document review and approval per DOC0841068	
Location	Description
Chapter 2	<ul style="list-style-type: none"> Table 2-17: Updated weight for PGR Table 2-18: Updated weight for PGR and AR magnet
Chapter 3	<ul style="list-style-type: none"> Figure 3-9: Removed service area on equipment room side
Chapter 4	<ul style="list-style-type: none"> Section 4.3: Updated weight Figure 4-4: Updated CoG dimensions Figure 4-6: Added front cabinet service area to side view Figure 4-7: Added front cabinet service area to side view
Chapter 5	<ul style="list-style-type: none"> Section 5.1.3 Added weight for HP Z4G5 Figure 5-3: Added callouts and callout table Figure 5-5: Updated dimensions
Chapter 7	<ul style="list-style-type: none"> Figure 7-10: Added filter to MRU
Revision 14, February 2024, English Document review and approval per DOC0841068	
Location	Description
Chapter 2	<ul style="list-style-type: none"> Section 2.2.2: Added ICC, removed HEC and CRY Table 2-1: Updated room size for Equipment Room and Total System Area Section 2.7.2: Updated step 1 (removed 1.5T products) Table 2-5: Added ICC, removed HEC and Blower Box Section 2.9.1: Added entire section for ICC Section 2.9.2: Added ICC references, removed HEC Section 2.10.1, Step 1.1: Added MDP models M7100ZA and M7100ZB, removed M7000ZA and M7000ZB Section 2.10.1, Step 5: Updated for ICC connection Table 2-10: Removed "HEC Continuous Power" Section 2.10.2: Updated entire section for MDP models M7100ZA and M7100ZB Section 2.10.3, Step 4: Updated for ICC connection Section 2.11.4: Replaced "Resonance Module" with "Magnet Assembly" Table 2-15: Updated values Table 2-17: Added ICC and pICC, removed HEC; updated weight for PGR and PEN Table 2-18: Added ICC and pICC, removed HEC; updated weight for PGR and PEN Table 2-21: Added ICC components, removed HEC pump and HEC blower
Chapter 3	<ul style="list-style-type: none"> Section 3.5.4: Removed steps related to service hatch Figure 3-9: Removed Blower Box Section 3.5.5: Removed figure for service hatch location Figure 3-17: Updated dimensions for distance between wheels

Revision 14, February 2024, English Document review and approval per DOC0841068	
Location	Description
Chapter 4	<ul style="list-style-type: none"> • Section 4.1: Added ICC, removed HEC • Figure 4-1 and 4-2: Updated for ICC • Section 4.2.2: Updated for new MDP model M7100ZA and M7100ZB • Section 4.3: Updated weight • Figure 4-4: Updated CoG dimensions • Section 4.4: Updated weight • Figure 4-5: Removed blower box, updated CoG dimensions • Section 4.5: Added entire section for ICC • Section 4.7.2, step 1: Updated for ICC
Chapter 7	<ul style="list-style-type: none"> • Section 7.1.1: Added ICC, removed HEC and CRY • Table 7-3: Updated L1, L3 and L11 (Changed HEC connections to ICC connections). Added L17, L18 and L19 • Table 7-10: Added ICC, removed HEC • Table 7-11: Added ICC, removed HEC • Figure 7-10: Updated for ICC
Revision 13, July 2023, English Document review and approval per DOC0841068	
Location	Description
	Throughout manual: <ul style="list-style-type: none"> • Added information for AR magnet
Revision 12.0, April 2023, English Document review and approval per DOC0841068	
Location	Description
Chapter 1	<ul style="list-style-type: none"> • Table 1-2: Updated manual numbers
Chapter 2	<ul style="list-style-type: none"> • Section 2.7.2, Step 1: Added additional products • Section 2.9.2: Added step 3 • Table 2-10: Removed "Standby Consumption" • Section 2.10.2, Table 2-12 and 2-13: Updated values for row "F" • Section 2.10.3: Added note to explain exempt countries • Figure 2-23 : Updated height and width • Table 2-16: Removed "Minimum Circuit Ampacity" • Table 2-18: Updated weight for Operator Workspace Display
Chapter 3	<ul style="list-style-type: none"> • Section 3.4, Step 7: Updated sentence for clarity • Section 3.5.4, Step 1.f.i: Updated sentence for clarity
Chapter 4	<ul style="list-style-type: none"> • Section 4.2.1: Added requirement for an electrical outlet • Section 4.2.2: Added sentence that only exempt countries can supply their own MDP • Section 4.3, Step 1: Updated weight (SPR HCSDM00714302) • Figure 4-3: Updated CoG dimensions (SPR HCSDM00714302) • Section 4.5, Step 1: Updated weight (SPR HCSDM00714305) • Figure 4-5: Updated CoG dimensions (SPR HCSDM00714305) • Section 4.7, Step 1: Updated weight (SPR HCSDM00714309)

Revision 12.0, April 2023, English Document review and approval per DOC0841068	
Location	Description
Chapter 5	<ul style="list-style-type: none"> Section 5.1.3: Changed anchor size from 9.5mm to M10 Section 5.1.4: Changed weight from 11.8 kg to 6.7 kg
Chapter 7	<ul style="list-style-type: none"> Table 7-4: Updated cable lengths for "PEN, top edge to RF common ground stud" and "SPW, top edge to RF common ground stud" Table 7-8: Updated "From" "To" entries and cable lengths Section 7.2.5, Step 1: Updated sentence for clarity
Chapter 8	<ul style="list-style-type: none"> Section 8.1: Removed definitions "Exclusion Zone", "Gradient" and "Security Zone"
Revision 11.0, June 2022, English Document review and approval per DOC0841068	
Location	Description
Chapter 2	<ul style="list-style-type: none"> Section 2.2.2: Added MRE Magnetic Resonance Elastography Section 2.2.3: Added Physiological Acquisition Transceiver (PAT) Charging Station Section 2.3: Expanded section with additional information Table 2-1: Updated values for Equipment Room Figure 2-9: Added Note to identify magnet isocenter Table 2-2: Corrected units (nanotesla instead of microtesla) Section 2.6.3: Corrected units (μT instead of mT) Section 2.9.2: Expanded section with additional information Section 2.10.1: Added Step 1 b i, Step 2 and Step 7 Table 2-9: Added MRU, MON, Optional MRE, OXY and Physiological Acquisition Transceiver (PAT) Charging Station Table 2-21: Updated specifications and added new column for Component Location Section 2.11.4: New section for Temperature and Humidity Storage Requirements
Chapter 3	<ul style="list-style-type: none"> Section 3.5.5: New section for Penetration Panel Closet Specifications Section 3.5.9: New section for Storage Cabinets Figure 3-15: Added asterisk for table length with handle in the extended position
Chapter 4	<ul style="list-style-type: none"> Section 4.1: Expanded section with additional information and Notes Figure 4-1: Updated dimensions Section 4.2.1: New section for MDP Requirements Figure 4-4: Corrected figure Figure 4-6: Updated dimensions Section 4-10: Updated weight
Chapter 5	<ul style="list-style-type: none"> Section 5.1.3: Updated center of gravity and weight Section 5.3: New section for Physiological Acquisition Transceiver (PAT) Charging Station
Chapter 7	<ul style="list-style-type: none"> Table 7-2: Added Physiological Acquisition Transceiver (PAT) Charging Station Table 7-7: Removed cable part numbers Section 7.1.4: Added section Brainwave Option Section 7.2.1: Added steps 6 and 7 Section 7.2.2: Expanded step 5a Section 7.2.4: Updated minimum height of cable supports to 2578 mm Table 7-12: Added row for OXY and Physiological Acquisition Transceiver (PAT) Charging Station

Revision 10.0, English Document review and approval per DOC0841068 R13, August 2020	
Location	Description
All chapters	<ul style="list-style-type: none"> Major content re-organization throughout the manual, and some content was moved to separate documents. See the Related Documents section in this manual for document numbers. New personnel matrix and icons added throughout the manual to help contractors and customers locate the content they need for their respective project work.
Chapter 2	<ul style="list-style-type: none"> MDP updates IEC requirements moved to a separate document Acoustic requirements moved to a separate document Grounding requirements related to RF shielding moved to a separate document
Chapter 3	<ul style="list-style-type: none"> Magnet Room venting requirements moved to a separate document RF shielded room requirements moved to a separate document
Chapter 4	MDP updates
Chapter 7	<ul style="list-style-type: none"> New specifications for available lengths Removed list of system cable specifications from the PIM
Appendix	<ul style="list-style-type: none"> Updated sample customer-supplied MDP schematic RF shielded room testing requirements moved to a separate document

Revision 9.0, Document review and approval per DOC0841068 R12, May 2019			
Location	Description		
		750w	Architect
Chapter 2, section 2.4	<ul style="list-style-type: none"> Updated the following tables for IEC EMC Edition 3.2: <ul style="list-style-type: none"> Guidance And Manufacturer's Declaration – Electromagnetic Emissions Guidance And Manufacturer's Declaration – Electromagnetic Immunity Guidance And Manufacturer's Declaration – Electromagnetic Proximity field Immunity 	x	x
Chapter 6	New chapter created for Digital Service and Connectivity Requirements to add InSite RSvP requirements.	x	x
Chapter 7, section 7.4	Added reference to InSite RSvP requirements for the GOC.	x	x

Revision 8.0, Document review and approval per DOC0841068 R10, April 2016				
Location	Description			
		450w	750	750w
All pages	Added Gemini product name to publication title.	x		x
Chapter 3, section 5.4.1	Correction of RF Shield Requirements: Removed "Cryo lines", which should not be installed when measuring SE.	x	x	x
Chapter 3, section 5.4.2	Moved items related to RF shield testing from 5.4.1 to 5.4.2, a new section called "RF Shield Test Requirements".	x	x	x
Chapter 3, section 6.5	Corrected dimensions on Magnet Mounting Detail illustration(s).	x	x	x
Chapter 3, section 6.5	Added a more accurate magnet image to the "Magnet Room Floor Levelness Area" illustration.	x	x	x
Chapter 3 section 8.5	Clarified information in "Cryogenic Venting" overview.	x	x	x
Chapter 3 section 8.6.5	Clarified vent construction requirements.	x	x	x

Revision 8.0, Document review and approval per DOC0841068 R10, April 2016				
Location	Description	450w	750	750w
		Chapter 4, section 7, illustrations	Corrected dimensions on PEN and SPW Wall illustrations; added "Pen Panel and SPW Wall Openings" illustration.	x
Chapter 6, section 3	Updated System Cable List table for DV26.	x	x	x
Revision 7.0, Document review and approval per DOC0841068 R9, March 2016				
Location	Description	450w	750	750w
		Chapter 2	Removed the compliance declaration information for Edition 2.1. Added "MR" before "system."	x
Chapter 2, Section 1	Changed the word "Upgrade" to "System Level Requirements for Installing into Existing MR Suite"	x	x	x
Chapter 2, Section 3, under Caution	Within the caution note, added "The ceiling service area should be kept clear of overhead items, including soffits, HVAC, plumbing components, and brackets."	x	x	x
Chapter 2, Section 3	Added system pictures to graphics.	x	x	x
Chapter 2, Table 2-1	Added W x D dimensions to Equipment Room size.	x	x	x
Chapter 2, Illustration 2-4 Energy Transmitted Through Vibro	Added the new vibration plot graphic.	x	x	x
Chapter 2, Illustration	For "MR System Water Cooling Block Diagram": added the shutoff valves for both supply and return.	x	x	x
Chapter 2, Section 11	Added M7000WL MDP option. Added emergency backup power requirements.	x	x	x
Chapter 2, Illustration	For "Gradient Coil Cart": Updated the gradient coil Diameter XRMb to 892mm.		x	
Chapter 2, Illustration	For "Gradient Coil Cart": Updated the gradient coil Diameter XRMw to 907mm.	x		x
Chapter 2, Table	For "MR System Component Replacement Shipping Specifications": Changed the PGR international shipping weight to 3654 lb., 1657 kg.	x	x	x
Chapter 3, Illustration 3-4	Modified the illustration to show the height of the stud from the floor to the top of the stud.	x	x	x
Chapter 3, Section 6.2	Changed text in step 6.2.3.g. from "SPW" to "PEN".			x
Chapter 3, Section 6.2	Added "maintain the minimum service area outside 200G" to step 6.2.1.b.	x	x	x
Chapter 3, Illustration(s)	Updated Magnet Mounting Detail illustration(s) and specs for clarification, and added Ergopad illustration and specs.	x	x	x
Chapter 3, Illustration 3-7	Changed the metric measurement from 2487mm to 2362mm.		x	
Chapter 3, Illustration 3-10	Removed words in the illustration, and replaced them with callouts.	x	x	x
Chapter 3, Section 8	Added new vent graphics for clarity. Added combined vent illustration and details.	x	x	x
Chapter 4	Updated Equipment Room examples.	x	x	x

Revision 7.0, Document review and approval per DOC0841068 R9, March 2016					
Location	Description	Applies To			
		450w	750	750w	
Chapter 4, illustrations 4-8, 4-10, 4-11	Changed the distance between the Pen and SPW from 5.4 to 7.8 inches.	x	x	x	
Chapter 6, Section 1.2, Table 6-4	Changed the usable length value for "PEN cabinet, top panel" to "GOC, rear panel" from "141.7 (29000)" to "1141.7 (29000)".	x			x
Chapter 6, section 2	Added clarity: Removed redundant information in the requirements, illustrations and tables; re-organized requirements into smaller sections; cross-referenced from dimensions in the requirements to dimensions shown in the illustrations.	x	x		x
Chapter 6, Illustration 6-3, 6-4	Added dimension detail from the edge of the cable tray to the magnet isocenter to make sure the gradient cable aligns with the magnet isocenter.	x	x		x
Chapter 6, Illustration 6-1	Added detail to specify the dimensioned offset from magnet isocenter and cable tray edge (if stacked) or between trays (if side by side).	x	x		x
Chapter 6, Illustration 6-7	Transposed dimensions 13 and 7.25 to their correct positions on the graphic.	x	x		x
Chapter 6, Section 4	Added backup power specs for cryocooler.	x	x		x
Chapter 7, Section 4.3.8	Added a figure to illustrate the setup to measure ground isolation resistance. Added step 4 at the end of the procedure to tell the installer to reconnect the lines to ground after the measurement is complete.	x	x		x

Revision 6.0, Document review and approval per DOC0841068 R8, March 2014					
Location	Description	Applies To			
		450	450w	750	750w
Chapter 8, 3	Clarified requirement of position of the emergency air intake. The emergency air intake must be at the highest point of the finished ceiling. Additionally, there must be an emergency air intake between the finished ceiling and the RF ceiling.	x	x	x	x
Chapter 8, 3, Illustrations	Updated Magnet Room Exhaust Fan Schematic and Magnet Room Exhaust Fan Schematic with Optional Oxygen Monitor illustrations to reflect correct position of the emergency air intake.	x	x	x	x

Revision 5.0, Document review and approval per DOC0841068 R6, July 2013					
Location	Description	Applies To			
		450	450w	750	750w
Chapter 2, 2.2.6	Corrected profiling issue to remove a blank requirement.	x	x	x	x
Chapter 2, 3	Updated values in the Minimum Room Size table to reflect decreased cabinet installation area requirements. Corrected inches to mm conversions.	x	x	x	x

Revision 5.0, Document review and approval per DOC0841068 R6, July 2013					
Location	Description	Applies To			
		450	450w	750	750w
Chapter 2, 3 Note	Replaced text in Note after the Minimum Room Size table to read "Contact PMI to ensure correct ramp leads have been ordered."	x		x	
Chapter 2, 6-1	Changed the table title to "Acoustic Specifications (Under Ambient Conditions)"	x	x	x	x
Chapter 2, 6.2 Illustration 2-4, Note at end of 3.6.2	Changed titles on illustration to "Vibration Transmitted through VibroAcoustic Mat." In Note, corrected references of "energy" to "vibration."	x	x	x	x
Chapter 2, 7.4	Clarified the information about non-MR System equipment in the vicinity of the MR System and the potential sensitivity to magnetic fields. Emphasized that the table listing the types of equipment is for reference only and the listed gauge limit may not be accurate for the customer specific equipment.	x	x	x	x
Chapter 2, 10.1	Updated the Facility Liquid Coolant Requirements table with the inlet temperature for the wide bore systems.			x	x
Chapter 2, 10-1, Illustration 2-14	Create new Allowable Facility Water Temperature and Flow illustration with temperature range for wide bore systems.			x	x
Chapter 2, 11.2	Added new requirement for Main Disconnect Panel (MDP). For the emergency-off circuit, the wire size should be 12-22 AWG and is supplied by the customer.	x	x	x	x
Chapter 2, 12.3	Removed Note from MR System Component Replacement Shipping Specifications table. Added it to weight tables in Magnet chapter.	x	x	x	x
Chapter 2, 12	Checked conversions from inches to millimeters and pounds to kilograms. Corrected where needed on tables Delivery Route Requirements, MR System Component Shipping Specifications, MR System Optional Component Shipping Specifications, and MR System Component Replacement Shipping Specifications.	x	x	x	x
Chapter 3, 2.3	Added table to clarify measurements for Magnet Steady State Vibration Specifications diagramed in Illustration 3-2.	x	x	x	x
Chapter 3, 5.3.1	Added +/- 0.5 MHz tolerance to frequencies listed in first requirement. Added Cyro Lines to the list of items that must be installed before measuring shielding effectiveness.	x	x	x	x
Chapter 3, 5.3.2	Added new requirement for the RF Shield Test Report. Added the need to perform a dock anchor pull test. Verify pull that is greater than 600 lb.. Clarified Ground isolation test results by adding the sentence "this test must be performed before the site is turned over to the general contractor for electrical installation." Add requirement "Blank pen wall is installed."	x	x	x	x

Revision 5.0, Document review and approval per DOC0841068 R6, July 2013					
Location	Description	Applies To			
		450	450w	750	750w
Chapter 3, 7.1	Expanded magnet weight totals in table listing major components of magnet assembly. Added Note from MR System Component Replacement Shipping Specifications table.	x	x	x	x
Chapter 3, 7.1	Added 450w GEM magnet cover dimensions.		x		
Chapter 3, 7.3	In MRU Facility Power Requirements table, added qualifier to connection type requirement to read "Hard wired/ permanently wired directly to facility power, no plugs or connectors allowed."	x	x	x	x
Chapter 3, 8.3.2	Changed "All exhaust vent system components must be non-magnetic" to "All items within the RF enclosure must be non-magnetic."	x	x	x	x
Chapter 3, 8.4.1	Added qualifier to requirement to read "A pressure equalizing vent is required in the Magnet Room ceiling or in the wall, at the highest point possible."	x	x	x	x
Chapter 3, 8.6.5	Added "ONLY" to the beginning of the requirement "ONLY one dielectric break in the vent system (that is, Ventglas) is required in the Magnet Room"	x	x	x	x
Chapter 3, 8.7.2	Added requirement "GE Engineering recommends that the cryogen vent be constructed to the same specification as required inside the Magnet Room."	x	x	x	x
Chapter 3, 8.7.3 Illustration 3-15	Updated the Cryogenic Venting (Interior) illustration to show the required distance for the vent stack between the waveguides.	x	x	x	x
Chapter 4, 3 Illustration 4-3	In the Power, Gradient, RF (PGR) Cabinet illustration, corrected the meter measurement of the width of the PGR Cabinet Side View.	x	x	x	x
Chapter 6, 2-1 Illustration 6-2	On the Cable Bends and Obstructions illustration, moved the cables closer to the leading edge of the cable tray.	x	x	x	x
Chapter 6, 2-2, Illustration 6-3	Added an inset to the Cable Tray Requirements (Side by Side) illustration, showing a side view of cable trays and the cable and lines as they are routed around the non-ferrous cable support.	x	x	x	x
Chapter 6, 2-2, Illustration 6-4	Added an inset to the Cable Tray Requirements (Stacked) illustration, showing a side view of cable trays and the cable and lines as they are routed around the non-ferrous cable support.	x	x	x	x
Chapter 6, 2-3 Illustration 6-6	In the Gradient Cable Drop to SPW (Side View) illustration, included the dimension from the bottom of the cable tray to the floor.	x	x	x	x
Chapter 6, 2.4	In the table for Minimum Cable Tray Width, changed the component type for the PGR from electrical to air/water.	x	x	x	x
Chapter 6, 4	In the note at the end of Facility Supplied System Interconnects Specifications, add a paragraph providing more details about the RF door switch.	x	x	x	x

Revision 5.0, Document review and approval per DOC0841068 R6, July 2013					
Location	Description	Applies To			
		450	450w	750	750w
Chapter 7 4.3.8	Change the requirement "To help prevent personal hazard, it is necessary for the enclosure to be properly grounded." to a warning."	x	x	x	x
All	In the continuing effort to prepare these documents for translation, some illustrations were updated with call-outs and descriptive tables.	x	x	x	x

Revision 4.0, Document review and approval per DOC0841068 R5, November 2012					
Location	Description	Applies To			
		750	450	450w	750w
Ch 2, Sec 3, Table 2-1	Updated the Minimum Room Size values to reflect decreased cabinet installation area requirements.	X	X	X	X
Ch 2, Sec 11.1	Reorganized and reworded the requirements to support MDP requirements	X	X	X	X
Ch 2, Sec 11.1, Table 17	<ul style="list-style-type: none"> Removed Power Demand information Changed Requirement for Power Availability to read "Facility power is required at all times for operation of the Cryocooler (CRY) to minimize cryogen consumption." Removed MDP information from Power Availability row. Removed Emergency Off Button information. 	X	X	X	X
Ch 2, Sec 11.1, Table 18	Added table to provide specifications for System Power Demand.	X	X	X	X
Ch 2, Sec 11.2	Added new section to add requirements for MDP electrical requirements.	X	X	X	X
Ch 3, Sec 5.3.3.5.I	Corrected Newton conversion for torque value specified for the dock anchor.	X	X	X	X
Ch 3, Sec 8.7.3, Illustration 3-15 and Table 3-11	Added Item 12 to illustration to show that height from floor to bottom of contractor-supplied wave guide is a maximum of 116.9 inches.	X	X	X	X
Ch 3 Sec 8.7.3, Illustration 3-16	Revised cryovent drawing for 450w.			X	
Ch 4, Sec 1, Illustration 4-1	Updated Equipment Room Layout Example Illustration to reflect decreased cabinet installation area requirements.	X	X	X	X
Ch 4, Sec 1, Illustration 4-1 and Table 4-1	<ul style="list-style-type: none"> Updated the Equipment Room layout notes to support new illustration (. Added requirement "If you use the minimum service area, then riggers are required to install the equipment." 	X	X	X	X
Ch 4, Sec 2	Added "optionally" to the introductory sentence.	X	X	X	X
Ch 4, Sec 7, Table 4-2	Changed Item 1 to show minimum ceiling height.	X	X	X	X
Ch 4, Sec 7, Illustration 4-10, 4-11, 4-12	Updated Closet Service Hatch Space illustrations to correct standard to metric conversions.	X	X	X	X

Revision 4.0, Document review and approval per DOC0841068 R5, November 2012					
Location	Description	Applies To			
		750	450	450w	750w
Ch 6, Sec 1.2, Table 6-3	Added a paragraph and table to explain the order configurations for cables.	X	X	X	X
Ch 6, Sec 2.2, Table 6-11, Item 9	Added a reference to Illustration 6-3, Item 8.	X	X	X	X

Revision 3.0, Document review and approval per DOC0841068 R4, May 2012					
Location	Description	Applies To			
		750	450	450w	750w
Ch 2, Sec 3, Illustration 2-1	Rear magnet dimension shown as 84". Changed to 73".			X	
Ch 2, Sec 9, Note 3	Removed "(contact install base support in Florence SC for altitudes from 8100 to 11808 feet (2468.9 to 3600 m) above sea level"	X	X	X	X
Ch 2, Sec 11	Updated Voltage/Frequency power requirements to include 415/400/380 VAC 60 ±3 Hz	X	X	X	X
Ch 2, Sec 11	Updated daily voltage variation to +10% / -10%	X	X	X	
Ch 2, Sec 11, Table 2-16	Clarified text to: "Optional customer supplied" backup power must be provided the facility breaker ...	X	X	X	X
Ch 2, Sec 12.3	Updated magnet shipping dimensions and weights; added lifting bar details; clarified additional shipping details	X	X	X	X
Ch 3, Sec 2.2	Clarified environmental steel limits	X	X	X	X
Ch 3, Sec 2.3, Illustration 2-3	Changed vibration specification to start at 0.5 Hz	X			X
Ch 3, Sec 5.3.3.1.5.l	Clamping force is 600 lb. (2669 N)	X	X	X	X
Ch 3, Sec 5.3.3 Figure 3-4	Updated drawing to show removable portion of dock anchor	X	X	X	X
Ch 3, Sec 6.2	Clarified text: "minimum service area" changed to "200 Gauss"	X	X	X	X
Ch 3, Sec 8.6.5.1.c	Removed note for second Ventglas joint requirement	X	X	X	X
Ch 4, Sec 7, Table 4-2	Clarified text: Magnet Room service areas for PEN and SPW must not overlap "200 G" in magnet minimum service area	X	X	X	X
Ch 7, Sec 8.2.f	Changed R/D = 0.5 for standard sweep elbows to R/D = 1.0 for standard sweep elbows	X			X

Revision	Date	Description		
2.0	Mar2012	Document review and approval per DOC0841068 R3		
		Reason for Change	Location	Description
			Various Illustrations	Replace text callouts with numbers and in-line notes to facilitate translations

Revision	Date	Description		
		Specification Update	Ch 2, Sec 1	Clarified HSS assessment for upgrades. HSS may now be used instead of vibration testing. Eliminated redundant requirements and focused on essential requirements
		Specification Update	Ch 2, Sec 3	Simplified room size table
		Specification Update	Ch 2, Sec 3	Updated Minimum Magnet Room Size for consistency across platforms (also revised illustration to show new design tradeoff area)
		Specification Update	Ch 2, Sec 6.1	Added note to clarify customer need to contain noise within each of the three rooms
		Specification Update	Ch 2, Sec 6.2	Added chart to show energy transmitted through VibroAcoustic Mat
		Specification Update	Ch 2, Sec 7.2	Updated moving metal sensitivity distances for cars and trucks
		Specification Update	Ch 2, Sec 9.2	Added CAD Stream heat load to System Heat Output for Air Cooling; Added Watts for MRE
		Specification Update	Ch 2, Sec 10.1	Added requirement for HEC, CRY, and PGR to be located on the same floor
		Specification Update	Ch 2, Sec 11	Updated Voltage Variation, Voltage Transients, and Regulation. Removed requirement for input power neutral wire. Added restriction to disallow use of corner grounded DELTA source. Added note to MDP illustration for backup power requirement
		Specification Update	Ch 2, Sec 12.3	Added delivery dimensions for the GEM Table, 750w service tool, and large FRUs
		Specification Update	Ch 3, Sec 4	Added detail to allow customer acoustic engineer to design an MR Suite to contain the air and structure borne noise levels
		Specification Update	Ch 3, Sec 5 & Ch 7, Sec 4	New RF Shield Effectivity test requirements with a focus on the test methods, equipment, test procedures and report format
		Specification Update	Ch 3, Sec 6	Added subsection for PEN closet requirements
		Specification Update	Ch 3, Sec 6.2 & Ch 4 Sec 7	Added requirements for Closet service hatch (modified illustration 4-8 and added 4-10, 4-11, 4-12; Updated 15 in. to 381 mm. Added maximum distance between PEN and SPW
		Specification Update	Ch 3, Sec 6.3	Added reference to cable concealment kit
		Specification Update	Ch 3, Sec 6.4	Changed floor levelness requirement from 0.3125 in. over 120 in. to 0.125 over entire floor levelness area
		Specification Update	Ch 3, Sec 7.2	Added GEM Table per 450w/750w GEM programs
		Specification Update	Ch 3, Sec 7.1	Updated magnet specifications
		Specification Update	Ch 3, Sec 8	Added note to prohibit use of magnetic parts in Magnet Room
		Specification Update	Ch 3, Sec 8.3, 8.5, 8.6, 8.7	Simplified requirements and updated Cryogen Venting illustrations (3-15, 16, 17) eliminating nonessential details, focusing on requirements

Revision	Date	Description		
		Specification Update	Ch 3, Sec 9.2	Added note to prohibit use of magnetic parts in Magnet Room
		Specification Update	Ch 4, Sec 1	Added note that optional equipment is not shown
		Specification Update	Ch 4, Sec 2, 3, 4	Added location of seismic anchors and updated service install access
		Specification Update	Ch 4, Sec 4	Removed lifting hole detail from HEC illustration and updated service clearances
		Specification Update	Ch 4, Sec 2, 3, 4	Added location of seismic anchors and updated service install access
		Specification Update	Ch 4, Sec 5	Dimensioned blower box hose service clearance
		Specification Update	Ch 4, Sec 6	Added center of gravity dimensions to SPW illustration
		Specification Update	Ch 4, Sec 7	Updated/added PEN and SPW opening illustrations to show required service area
		Specification Update	Ch 5, Sec 5.3.3, Ch 3, Sec 6.4.5	Removed magnet floor mounting requirements. All magnets must use the VibroAcoustic dampening option
		Specification Update	Ch 6, Sec 2	Updated allowable obstruction distance to 7 in. above the bottom of the tray
		Specification Update	Ch 6, Sec 2.2	Updated cable tray height (104 in to 101.5 in) and distance from back of magnet (44.5 to 48.9)
		Specification Update	Ch 6, Sec 2.4	Added minimum cable tray widths between Equipment Room components
		Specification Update	Ch 7, Sec 8	Updated the Magnet Cryogenic Venting Pressure Drop Table. 3T LC pressure drop matrix now matches the 3T G3/LCC pressure drop matrix
1.0	Nov2010	Initial Release. Document review and approval per DOC0841068 R1 V5		



GE HealthCare

www.gehealthcare.com